

AD-A216 847

**CHEMICAL
RESEARCH,
DEVELOPMENT &
ENGINEERING
CENTER**

**CRDEC-CR-051
(GC-TR-89-1728-001)**

**TRACE ORGANIC ANALYSIS
OF MICROENCAPSULATED MATERIALS**

**Shirley A. Liebman
GEO-CENTERS, INC.
Newton Centre, MA 02159**

**A. Peter Snyder
Dennis J. Reutter, Ph.D.
Michael B. Wasserman
Leon J. Schiff**

RESEARCH DIRECTORATE

November 1989

**DTIC
ELECTE
JAN 13 1990
S E D**

**U.S. ARMY
ARMAMENT
MUNITIONS
CHEMICAL COMMAND**



Aberdeen Proving Ground, Maryland 21010-5423

DISTRIBUTION STATEMENT A

**Approved for public release;
Distribution Unlimited**

90 01 16 133

Disclaimer

The findings in this report are not to be construed as an official Department of the Army position unless so designated by other authorizing documents.

Distribution Statement

Approved for public release; distribution is unlimited.

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE

REPORT DOCUMENTATION PAGE				Form Approved OMB No. 0704-0188	
1a. REPORT SECURITY CLASSIFICATION UNCLASSIFIED			1b. RESTRICTIVE MARKINGS		
2a. SECURITY CLASSIFICATION AUTHORITY			3. DISTRIBUTION/AVAILABILITY OF REPORT Approved for public release; distribution is unlimited.		
2b. DECLASSIFICATION/DOWNGRADING SCHEDULE			5. MONITORING ORGANIZATION REPORT NUMBER(S)		
4. PERFORMING ORGANIZATION REPORT NUMBER(S) CRDEC-CR-051 (GC-TR-89-1728-001)			7a. NAME OF MONITORING ORGANIZATION		
6a. NAME OF PERFORMING ORGANIZATION GEO-Centers, Inc.		6b. OFFICE SYMBOL (If applicable) SMCCR-RSC-C		7b. ADDRESS (City, State, and ZIP Code)	
6c. ADDRESS (City, State, and ZIP Code) 7 Wells Avenue Newton Centre, MA 02159			9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER DAAA15-87-D-0007-0001		
8a. NAME OF FUNDING/SPONSORING ORGANIZATION CRDEC		8b. OFFICE SYMBOL (If applicable) SMCCR-RSC-C		10. SOURCE OF FUNDING NUMBERS	
8c. ADDRESS (City, State, and ZIP Code) Aberdeen Proving Ground, MD 21010-5423			PROGRAM ELEMENT NO.	PROJECT NO.	TASK NO.
			WORK UNIT ACCESSION NO.		
11. TITLE (Include Security Classification) Trace Organic Analysis of Microencapsulated Materials					
12. PERSONAL AUTHOR(S) Liebman, Shirley A. (GEO-CENTERS, INC.); Snyder, A. Peter; Reutter, Dennis J., Ph.D.; and Wasserman, Michael B.; and Schiff, Leon J. (CRDEC)					
13a. TYPE OF REPORT Contractor		13b. TIME COVERED FROM 87 Jan to 89 Jan		14. DATE OF REPORT (Year, Month, Day) 1989 November	
				15. PAGE COUNT 106	
16. SUPPLEMENTARY NOTATION COR: J. Pistritto, SMCCR-RSC-C, (301) 671-2443					
17. COSATI CODES			18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number)		
FIELD	GROUP	SUB-GROUP	Trace organic analysis Dynamic headspace		
15	06	03	Automated sample concentrators Automated thermal processing		
			Analytical pyrolysis On-line analysis		
19. ABSTRACT (Continue on reverse if necessary and identify by block number) Detection and identification of trace target species in complex polymeric mixtures requires application of advanced analytical instrumentation and methodology. The major goal of this task is the development of analytical systems and/or methodologies to detect/identify/quantitate trace organic chemical species in complex microencapsulated mixtures. To accomplish this goal, the approach is based on sample processing, concentration, separation, detection/identification, and data manipulation/processing. This task emphasized the application and development of commercially available analytical systems. Specific objectives were: (1) Develop current awareness of the extensive microencapsulation research and technology base; (2) Provide detection and characterization of microcapsules with their chemical/physical microstructure; (3) Establish detection and identification schemes for signature analysis and classification of core/shell formulated products, (continued on reverse)					
20. DISTRIBUTION/AVAILABILITY OF ABSTRACT <input checked="" type="checkbox"/> UNCLASSIFIED/UNLIMITED <input type="checkbox"/> SAME AS RPT <input type="checkbox"/> DTIC USERS			21. ABSTRACT SECURITY CLASSIFICATION UNCLASSIFIED		
22a. NAME OF RESPONSIBLE INDIVIDUAL SANDRA J. JOHNSON			22b. TELEPHONE (Include Area Code) (301) 671-2914		22c. OFFICE SYMBOL SMCCR-SPS-T

18. Subject Terms (continued)

Supercritical fluid technology
SF extraction/desorption
SF reaction/chromatography
Box-Behnken experimental design
Microencapsulated pesticides
Surfactants
Emulsifiers
Polymer shell/walls
Microcapsule cores
Fiber optic monitors (FOM)
Applied artificial intelligence (AI)
Neural network
Pattern recognition
EXMAT
MicroEXMAT

19. Abstract (continued)

including thickeners, additives, and carriers; (4) Determine the quantitative amounts of microcapsule samples (milligrams) at trace parts per million, parts per billion levels; (5) Encode information in an applied artificial intelligence format.

Four analytical instrumental systems were assembled to examine microencapsulated materials: System 1 - Pyrolysis atmospheric pressure chemical ionization (APCI) tandem mass spectrometer; System 2 - Pyrolysis/concentrator interfaced to GC-MS and GC-MS/MS units; System 3 - GC interfaced to Fourier transform infrared (FTIR) spectrometer; System 4 - Supercritical fluid chromatograph for extraction, desorption, reaction and chromatography. These analytical systems integrate automated and spectral detection/identification units.

Exploratory research was conducted for prototyping trace detection/identification instrumentation, which included fiber optic monitors, AI expert systems, and neural networks for in-field threat analysis in diverse environments.

PREFACE

The work described in this report was authorized under Project No. DAAA15-87-D-0007-0001. This work was started in January 1987 and completed in January 1989.

The use of trade names or manufacturers' names in this report does not constitute an official endorsement of any commercial products. This report may not be cited for purposes of advertisement.

Reproduction of this document in whole or in part is prohibited except with permission of the Commander, U.S. Army Chemical Research, Development and Engineering Center (CRDEC), ATTN: SMCCR-SPS-T, Aberdeen Proving Ground, Maryland 21010-5423. However, the Defense Technical Information Center and the National Technical Information Service are authorized to reproduce the document for U.S. Government purposes.

This report has been approved for release to the public.

Acknowledgments

We appreciate the background information microencapsulation supplied by Dorothea Paterno, the initial activities in supercritical fluid technology provided by Emory W. Sarver and Leon J. Schiff of CRDEC, and the cooperative research by R. Fifer of the U.S. Army Ballistic Research Laboratories in the fiber optic monitor development. Loan of the Guided Wave unit by Ted Novak and technical assistance by Vikki Henderson of CRDEC are also appreciated in addition to the cooperation provided by Raymond Herd and Phyllis Gray, Analytical Division, CRDEC, and Barbara Ball and Sol Heller, Audio-Visual, CRDEC. The outstanding continued assistance by Linda Jarvis of GEO-Centers, Incorporated, is gratefully acknowledged. In-depth technical review and general project management by Louis Isaacson, Geo-Centers, Incorporated, are sincerely appreciated.



Acc.	
Index	
Dist.	
Uncl.	
J.	
Ex.	
Info.	
Att.	
Dist.	
A-1	

Blank

CONTENTS

	Page
1. INTRODUCTION	13
1.1 Definition of the Problem and Program Goals	13
1.2 Experimental Design	14
1.3 Data Analysis	15
2. MICROENCAPSULATION RESEARCH AND TECHNOLOGY BACKGROUND.....	15
2.1 General Considerations	15
2.2 CRDEC Microencapsulation Workshop, December 3-4, 1986	16
3. INSTRUMENTATION BACKGROUND	17
3.1 General Trace Organic Analysis	17
3.2 Analytical Pyrolysis Technology	17
3.3 Supercritical Fluid Technology	19
4. EXPERIMENTAL DESIGN AND APPLIED ARTIFICIAL INTELLIGENCE.....	20
4.1 Experimental Design and Method Development	20
4.2 Applied Artificial Intelligence and Expert System Networks	20
4.3 Technical Resources: Materials, Instrumentation, and Software	21
4.3.1 Materials	21
4.3.2 Instrumentation	23
4.3.3 Software: Artificial Intelligence and Experimental Design	24
5. RESULTS OF TRACE ORGANIC ANALYSIS METHOD DEVELOPMENT	25
5.1 Experimental Design	25
5.2 Experimental Data and Discussion	25
5.2.1 System 1. Py/APCI/MS/MS	25
5.2.2 System 2. Pyrolysis-Concentrator GC-MS, and -MS/MS	27
5.2.3 System 3. Pyrolysis/Concentrator GC-FTIR	29
5.2.4 System 4. Supercritical Fluid Instrumentation	29
5.2.4.1 SFE. Supercritical Fluid Extraction	29
5.2.4.2 SFR. Supercritical Fluid Reaction	31
5.2.4.3 SFD. Supercritical Fluid Desorption	31
5.3 Supercritical Fluid Extraction-Chromatography SFE/SFC	32

	Page
6. EXPLORATORY RESEARCH AND TECHNOLOGY ASSESSMENT	32
6.1 Fiber Optics	32
6.2 AI and Neural Networks	32
6.3 Inverse GC and SFC	33
7. SUMMARY AND CONCLUSIONS	34
8. RECOMMENDATIONS	34
9. REFERENCES	36
10. BIBLIOGRAPHY	38
10.1 Bibliography on Microencapsulation Research and Technology	38
10.2 Bibliography on Supercritical Fluid Technology	40
10.3 Bibliography on Sensors/Artificial Intelligence	41
10.4 Bibliography on Inverse GC and SFC	44
10.5 Bibliography of Geo-Centers/CRDEC/BRL Presentations/ Publications Resulting from Task and Allied Program Related to Polymer/Microencapsulation/ Trace Organic Analysis	44
10.5.1 Presentations	44
10.5.2 Publications	50

LIST OF FIGURES

	Page
1. Integrated Intelligent Instruments (I ³).....	53
2. Trace Organic Analysis of Analytes in Gas, Liquid, Solid Matrices.....	53
3. Microencapsulated Systems and the I ³ Approach.....	54
4. CRDEC Instrumentation Schematic.....	54
5. Typical Analysis of Variance Output from Box-Behnken.	54
6. a. Polymer Characterization and Degradation (1971)..	55
b. Problem-Solving in Material Science (1986).....	55
7. a. Instrumentation for Sample Processing and On-Line Analysis.....	56
b. Advanced Sample Processing and Analysis.....	56
8. Polymer Microstructure and Composition.....	57
9. Pulse and Programmed Pyrolysis of Polyethylene	
a. 800°C/20 sec.....	57
b. 1°C/min to 800°C.....	57
10. a. CDS 320 Sample Concentrator Volatiles in Pharmaceuticals (Lot #1).....	58
b. CDS 320 Sample Concentrator Volatiles in Pharmaceuticals (Lot #2).....	58
11. Pyrolysis and Reproducibility of Glass-Lined System CDS Pyroprobe Model 124	
a. Human Hair.....	59
b. Raw Wool.....	59
c. Microorganisms.....	60
d. Silk	60
e. Cellulose.....	60
f. Chitin.....	60
12. Box-Behnken Experimental Design Used in Analytical Pyrolysis Method Development.....	61

	Page
13. Box-Behnken Experimental Design - SFE Method Development.....	61
14. General Processes Involved in Organic Analysis.....	62
15. a. MultiLevel Expert Systems Outline.....	62
b. EXMAT, A Linked Network of Expert Systems.....	63
c. Development of EXMAT, A Linked Network of Expert Systems for Material Analysis.....	63
d. Pyrolysis/Concentrator with Capillary GC-MS/DS Formatted for Entry into Database Management System, EXDBM.....	64
e. Factor Analysis - Test Data for Chemometrics.....	64
f. Factor Rotation - Test Data for Chemometrics.....	64
16. EXMAT Network, MICROCAP Expert System.....	65
17. EXMAT - HELP "CHOICE SYS 8/CT".....	65
18. Box-Behnken Design for SFE - Block 1.....	66
19. Box-Behnken Design for SFE - Block 2.....	66
20. a. Time-Resolved Analytical Pyrolysis APCI MS/MS Schematic.....	67
b. Pyrogram Pulse 700°C/20 sec Polyethylene (Air, N ₂), Oxidized PE (Air, N ₂) and PEO (Air, N ₂).....	68
21. a. Pyrograms - Microcap #1 (180°C/min to 400°C) No Heat.....	68
b. Pyrograms - Microcap #1 (180°C/min to 400°C) With Heat.....	69
c. Pyrograms - Microcap #1 (Pulse to 700°C/20 sec)..	69
d. Pyrograms - Microcap #1 (240°C/min to 500°C).....	70
e. Pyrograms - Microcap #1 (300°C/min to 200°C), 4 min interval.....	70
22. Pyrogram Polymer A 180°C/min to 400°C.....	71
23. a. Pyrogram Microcap #1 pulse to 400°C/20 sec.....	71
b. Pyrogram Microcap #2 pulse to 400°C/20 sec.....	72

	Page
c. Pyrogram Microcap #3 pulse to 400°C/20 sec.....	72
d. Pyrograms of Microencapsulated Protein, Gelatin Powder and Microencapsulated Gelatin.....	73
24. a. Instrumentation for Trace Organic Analysis - Analytical Pyrolysis/Concentrator GC-MS/DS.....	73
b. Pyrograms - Concentrator/GC-MS (HP 5985B) Data Generation - Groups I, II, III Organic Amine Salts.....	74
c. MS Patterns of TIC Peaks from Pyrograms and Library Match.....	75
25. Amine Salt Structures for Test Method (EXMAT Series).....	76
26. a. High Density PE, 800°C/20 sec Pyrogram - RIC.....	77
b. Scansets 650-950, Carbon Triplets ~C ₁₃ -C ₁₅ of 26(a).....	77
c. Oxidized PE, 800°C/20 sec Pyrogram - RIC.....	78
d. Scansets 650-950, ~ C ₁₃ -C ₁₅ of 26(c).....	78
e. Polymer Bulk Ref. 60°/min to 400°C.....	79
f. Scansets 650-950, ~ C ₁₃ -C ₁₅ of 26(e).....	79
27. GC-MS Experimental Conditions (Finnigan TSQ GC/MS)...	80
28. a. Diazinon-Pyrogram, 240°/min to 700°C, 4 min (Scanset #672).....	81
b. EI-MS of Scanset #672 - Diazinon.....	81
29. a. Microcap #1, 240°/min to 700°C, 4 min (Scanset #657).....	82
b. EI-MS of Scanset #656 - Diazinon core.....	82
30. a. Microcap #1 RIC 180°C/min to 600°C (Scanset #662).....	83
b. EI-MS of Scanset #662 - Diazinon core.....	83
c. Microcap #1 RIC 240°/min to 700°C (Scanset #668).....	84

	Page
31. a. Microcap #1 RIC Pulse 500 ⁰ /20 sec (Scanset #753).....	84
b. EI-MS of Scanset #753 - Diazinon Core.....	85
c. MS-MS Daughter Ions from Parent 304 Ion of Microcap #1 Sample.....	85
32. a. Microcap #4 - Cypermethrin Core Pulse 500 ⁰ /20 sec (Scanset #608).....	86
b. EI-MS of Scanset #608 - Cypermethrin Core.....	86
33. a. Microcap #4 (RIC) Cypermethrin Core (Scanset #395).....	87
b. EI-MS of Scanset #395 from Microcap #4 RIC.....	87
34. a. Microcap #5 (RIC) Dursban Core (Scanset #831)....	88
b. EI-MS of Scanset #831 Dursban Core.....	88
c. MS-MS Daughter Ions from Parent 199 Ion of Microcap #5 Sample.....	89
35. a. Pyrogram - RIC of Dodge Gelatin in PS-PES Shell..	90
b. Pyrogram - RIC of Emulsifier Gelatin (Bend).....	90
c. Pyrogram - RIC of Gelatin (P. Snyder).....	91
d. Pyrogram - RIC of Gelatin (~ 500 μ Particles - Dodge).....	91
e. Pyrogram - RIC of Gelatin (~ 500 μ Particles - Dodge) Desorb 260 ⁰ C/5 min (Scanset #206).....	92
f. EI-MS of Scanset #206 - Gelatin Particles.....	92
36. Pyrolysis GC-FTIR	
a. Total Response Chromatogram (TRC) Microcap #1 (Diazinon Core).....	93
b. FTIR of ~ 12.5 min Peak with Library Match to Phosphorothioic Derivatives.....	93
37. SFE-SFC of Propellant Formulations.....	94

	Page
38. SFE-SFC of Microcap #1, Two Samples from Same Batch.....	94
39. SFE Method Development	
a. Box-Behnken Experimental Design (Runs 2,4,6).....	94
b. SFE Profiles for Microencapsulated Pesticides....	95
c. SFE Profiles for Microcap #1 - Box-Behnken Experimental Design.....	95
d. SFE Profiles Microcap #2 - Box-Behnken Experimental Design.....	96
e. SFE Profiles for Microcap #3 - Box-Behnken Experimental Design.....	96
40. SFE-SFC of Emulsifiers A and B.....	97
41. SFR-SFC of RDX and RDX and Catalyst.....	97
42. a. SFD (SFE)-SFC of Simulant in Soil.....	98
b. SFD-SFC and Capillary GC of Charcoals from Air Sampling.....	99
c. Conclusions of SF Technology Study for Microencapsulated Materials.....	100
43. Remote Spectroscopy - Fiber Optic Monitor (FOM) Prototypes for CCS Model 5000 SFC/GC (Supercritical Fluid Extractor/Accumulator Module Configurations 1, 2 and 3).....	101
44. a. Trace Aromatics - Configuration 2; Trace Toluene, Methyl Naphthalene, Added Fluorene, (Quartz/Cross Union).....	102
b. Trace Aromatics - Configuration 2; Trace Toluene, Fluorene, Anthracene Added.....	102
45. Mid-Infrared FOM Prototype in Configuration 2 (Modified) for Fabric Softener on ZnSe Rod and Fiber Optic Cable.....	102
46. Summary and Conclusions - Remote Spectroscopy and FOM Applications.....	103

47. Applied AI - MicroEXMAT

- a. Decision Structure for Expert System -
Configuration "EXCONFIG"..... 104
- b. Partial Decision Structure for Expert
System - Method Development "EXMETH"..... 105
- c. Embedded Intelligence and Applied AI Networks.... 106

TRACE ORGANIC ANALYSIS OF MICROENCAPSULATED MATERIALS

1. INTRODUCTION

1.1 Definition of the Problem and Program Goals.

Detection and identification of trace target species in complex polymeric mixtures requires application of advanced analytical instrumentation and methodology. In particular, microencapsulated toxic substances present a demanding challenge for their detection and identification in the environment. The accomplishment of this task requires knowledge of microencapsulation technology including the formulation and processing of complex materials, so that analytical characterization can be designed in a comprehensive and meaningful manner. The initial activity on this task utilizes background information in microencapsulation research and technology that involves the food, agricultural, pharmaceutical, consumer product, and biomedical industries.

The major goal of this task is the development of analytical systems and/or methodologies to detect/identify/quantitate trace organic chemical species in complex microencapsulated mixtures. The analytical chemistry required to characterize formulated materials has been developed over many years in context to polymeric substances, both natural and synthetic origins. The instrumentation and methods developed and applied have focused on analytical pyrolysis interfaced to chromatographic and spectral separation/detection/identification systems (1). Trace organic analysis involves the use of enrichment or concentration units which are manufactured as highly automated systems. This Integrated Intelligent Instrument (I³) approach, outlined in Figure 1, is based on sample processing, concentration, separation, detection/identification, and applied artificial intelligence.

The approach selected for this task emphasizes selected analytical tools that are commercially available and used for detection and identification of trace amounts of toxic species that are in formulated polymeric forms, such as immobilized carriers or microencapsulated materials. Development of the I³ approach focused on two major technologies; namely, analytical pyrolysis and supercritical fluids, with advanced data analysis and applied artificial intelligence (AI) methods. Figure 2 outlines how formulated (e.g. microencapsulated) materials are studied within the I³ approach.

Specific program objectives are outlined in the following five statements:

1. Develop current awareness of the extensive microencapsulation research and technology base.

2. Provide detection and characterization of microcapsules with their chemical/physical microstructure.

3. Establish detection and identification schemes for signature analysis and classification of core/shell formulated products, including thickeners, additives, and carriers.

4. Determine the quantitative amounts of microcapsule samples (milligrams) at trace ppm, ppb levels.

5. Encode information in an applied artificial intelligence format.

The means to accomplish these objectives include (a) assessment of microencapsulated R&D as background for developing analytical schemes; (b) assessment and assembly of selected analytical instrumentation to qualitatively and quantitatively determine component levels in complex mixtures; and (c) document and outline the information with applied AI. Our approach is shown schematically in Figures 3 and 4.

The long-term goal includes extension of laboratory results from these studies to the in-field, battlefield environment. This longer-term goal was considered in the initial design and execution of the task. The role of highly automated, expert-system oriented analytical instrumentation is significant in this regard. Specialty hardware/software designed to accomplish the detection and identification of trace toxic materials in complex matrices can be realistically proposed based on results from this and allied studies. Following the current emphasis on technology assessment and transfer methods for R&D links to commercially significant ends, such specialty tools based on the I^3 approach have significance throughout our industrial processing and environmental fields. This contract research provides the technical basis for the design and manufacture of a prototype highly automated portable, modular threat analyzer to meet Army needs.

1.2 Experimental Design.

Although there are numerous methods available to statisticians in the conduct of experimental design for chemical studies (2-5), the choice of a Box-Behnken design was made for specific reasons important to the interpretation stage. The Box-Behnken is a central composite, blocked design to randomize experimental error. The design provides an excellent multivariate tool for evaluation of experimental parameters. The pure experimental error is determined, in addition to measuring how well the variables account for the system variance (or lack of fit to the model). A typical output is shown in Figure 5. The response surface is determined over the instrumental range under study, so as to detect main variables responsible for the measured response, and importantly, any second order or quadratic interactions. Use

of experimental designs are emphasized in current industrial R&D to provide the most information in the shortest possible time, thereby resulting in productive research. Background experience in this area led to the specific choice of the Box-Behnken and also led to extensions of its use for multivariate response data needed in our studies.

1.3 Data Analysis.

The role of chemometrics has increased significantly over the past decade (6-8) to increase the information content from data generated by highly automated instruments. Only such reproducible, accurate data can be handled by chemometric techniques, including factor analysis and pattern recognition software. The mathematical and statistical treatments common in classification or multivariate analysis were further developed in this project to include applied artificial intelligence. Use was made of AI shell software packages, as well as specialty software, to develop a prototype expert system network, EXMAT. The need for combined numeric and symbolic computing was recognized early in our program and ideas were presented at a Workshop held in Seattle, WA in August 1985, sponsored by the American Association for Artificial Intelligence and Boeing Computer Services. Extension of this approach is being made to current developments in neural networks addressing both hardware and software aspects in cooperation with General Research Corporation, McLean, Virginia.

2. MICROENCAPSULATION RESEARCH AND TECHNOLOGY BACKGROUND

2.1 General Considerations.

Microencapsulation technology has developed over the past decade in a rapid and diverse manner throughout the food, chemical (pesticides, consumer products, etc.) agricultural, and pharmaceutical industries (9). Aerosol carriers, agent thickeners, and additives are composed of varied organic substances, many of which are polymeric in nature. Complex mixtures result when such materials (surfactants, surface and viscosity modifiers, emulsifiers, plasticizers, humectants, etc.) are formulated to achieve specific performance and compatibility requirements. Wall materials or specialty membranes and polymers may surround the encapsulated core or serve as transport agents in novel ways (10). These multicomponent products are formulated to provide "tailored" physicochemical behavior over a range of environmental stabilities, controlled release mechanisms, and transfer/distribution means. The advantages of providing reactive or labile species in timed-release/encapsulated formulations are given extensively in the agrochemical, biomedical, and pharmaceutical industries.

Fabrication techniques needed to produce such materials are likewise varied, including microencapsulation processes; e.g.

Wuster coating, centrifugal extrusion, vapor deposition (Parylene), liquid walls, interfacial polymerization, phase separation, complex coacervation (thermal, nonsolvent, interfacial), in-situ and matrix polymerization. Many processes are commercially developed and available for specialty applications (9).

Characterization of microcapsules thus requires detailed, comprehensive analytical methods for microstructure, composition, and physicochemical information. Concerted organic trace analysis conducted in the industrial R&D sector ensures that needed production quality and in-field applications are achieved over the planned range of formulation, processing, and end-use conditions. Furthermore, prediction of physical/mechanical properties for the intended lifetime of the material requires knowledge of degradation mechanisms and associated kinetics. Such controlled release performance critically depends on polymeric microstructure and compositional details of the thickeners, additives, and carriers used in the process. Thus, instrumentation designed for analysis of complex materials are required in our study.

2.2 CRDEC Microencapsulation Workshop, December 3-4, 1986.

Of singular importance to the Army was the Microencapsulation Workshop that was organized and held at CRDEC in December, 1986. Collection of the Workshop presentations has been published as "Macromolecules in Microencapsulation Research and Technology," CRDEC-SP-87022, August 1987.

The purpose of the CRDEC Microencapsulation Workshop was to encourage technology transfer among active industrial R&D, academic and U.S. Government investigators. Emphasis was placed on the physicochemical characterization of natural and synthetic macromolecules used in the microencapsulation field, with the associated synthesis, formulation, and processing aspects. Specialty design of polymeric microstructure and formulation/processing variables required for wide-ranging performance was emphasized with regard to controlled release, transport mechanisms, microporosity, and environment stability.

Over 50 companies, 20 academic/research institutes, and 8 government agencies participated. Nearly 160 persons attended the Workshop. Topics of research presentations included preparation, chemical/physical modification, and characterization of specialty polymers, surfaces, and membranes; design and monitoring of selected diffusion, permeability, adsorption/absorption or stabilization phenomena; and macromolecules and associated materials used as carriers, stabilizers, or thickeners in microemulsions and colloidal systems.

Directly relevant for Army applications is research that involves the production, delivery, and detection of substances encapsulated/immobilized in complex polymeric/inorganic systems.

The analytical capabilities needed for these materials also impact industrial manufacturers in food, chemical, agricultural, pharmaceutical and consumer product fields, as well as forensic sciences. Environmental and legal requirements, as well as quality control/assurance operations require that modern analytical systems be integrated into most corporate R&D programs.

Relevant literature in the microencapsulation field is given in Bibliography 10.1, (page 38), including references in microemulsions, microcapsules, polymeric/liposome delivery systems, physicochemical research, and specialty polymer information.

3. INSTRUMENTATION BACKGROUND

3.1 General Trace Organic Analysis.

The concerted organic analysis approach has been discussed (11) and developed over nearly two decades for applied research and problem-solving in industrial R&D (1). A schematic (Figure 6a) outlines the approach as it was presented in 1971 (12) and updated in Figure 6(b) in 1986. Figure 7(a, b) show the background information for sample processing and on-line analysis. The specific advantage of our study over previous and on-going efforts in polymer characterization is that of developing and/or applying new and automated instrumentation that is now available. Analysis of synthetic and natural polymeric mixtures have benefitted from the advantages of combined dynamic headspace analysis (DHS) with pulse and/or programmed pyrolysis. Careful thermal processing under controlled conditions is accomplished on small (milligrams or micrograms) amounts of sample. For many cases, this thermal method replaces solvent extraction or time-consuming isolation methods. For samples that are thermally labile, nonthermal treatment is used, and supercritical fluid (SF) instrumentation (e.g., extraction) achieves this goal. Our general problem-solving in materials science uses these two major instrumental fields as outlined in Figure 6(b).

3.2 Analytical Pyrolysis Technology.

The pyrolysis instrumentation available for performing this task is the Chemical Data Systems (CDS) Model 123 Pyroprobe, configured with a programmable pyrolyzer interfaced to the Model 320 Sample Concentrator. For samples that may be thermally labile and for ultra-trace analyses, the Model 330 Sample Concentrator is used with cryogenic enrichment and direct transfer of the trapped volatiles (splitless mode) to the separation or identification system (GC, FTIR, MS, etc.). The latter system, the Pyroprobe Model 124, provides one of the most sensitive means for trace organic analysis of complex materials. The acquisition and method development performed as part of this task established the significance of the pyrolysis technology to study the selected microcapsules (see Section 5).

In a Dynamic Headspace (DHS) experiment, a small (mg) sample is placed in the thermal desorber unit of the sample concentrator and heated to the desired temperature which is below the degradation level for the sample. Volatiles (residual monomers, water, additives, solvents, etc.) are released from the matrix and trapped or concentrated by means of the on-line internal trapping system. This enrichment stage uses either solid adsorbents and/or cryogenic means such as liquid N₂. When desired, the trapped volatiles are pulse-heated off/out of the trap and delivered to the analytical system. In some cases, DHS may not be needed and the pyrolyzate (from pulse or programmed treatment) is delivered directly to the GC or analytical unit. Hence, pure materials for reference patterns, as well as those that must be manipulated from complex matrices are examined in an absolutely repeatable, reliable manner in the microprocessor-controlled system under the conditions to obtain the most useful diagnostic information.

Furthermore, in an example of DHS combined with pyrolysis, Figure 8 illustrates results from thermal processing of a high molecular weight polyethylene-vinyl acetate (PE-PVA) copolymer with He and air atmospheres, respectively. The DHS analysis revealed low-level volatiles (ppm amounts) that were removed from the polymer matrix in the low-temperature (120°C for 5 min) mode and the pulse pyrolysis step (750°C) showed the characteristic triplets from PE, as well as the PVA fragments.

Importantly, programmed pyrolysis or "time-resolved" pyrolysis exhibits advantages over pulse methods for detection of microstructural details and mechanistic information (1). In Figure 9, programmed pyrolysis results (1°C/min to 850°C) are compared to pulse data (800°C/20 sec) on an NBS reference polymer, high density PE. A major difference is seen in the fragmentation pattern which has been used for diagnostic purposes. Development of linear, low-density PE (LLDPE) systems provides improved processing and performance behavior by means of incorporating small amounts of comonomers such as butene, pentene, and octene with ethylene. Compounding co- or terpolymers results in "tailor-made" materials for desired ranges of performance which are illustrative of the microencapsulation field.

Common polymers which are employed as wall materials, carriers, thickeners, or stabilizers for microencapsulation applications have been studied with analytical pyrolysis over the past two decades and are discussed in reference 1, (e.g., see Chapter 5 (Synthetics) and Chapter 8 (Biopolymers)). For natural polymers used in the food and pharmaceutical industries, pyrolysis methodology has proved to be important. Figure 10 (a, b) illustrates typical results from a concentrator/capillary GC study obtained from a DHS analysis of volatiles from two lots of a formulated product. Lot #1 (Figure 10a) shows three major peaks that were reduced significantly in Lot #2 (Figure

10b). Further, quantitative assays may be needed, such as pharmaceutical products that are processed with various solvents that must be removed to very low levels. A DHS experiment provides the needed information on residual solvents in the part-per-million or part-per-trillion range in this automated system (13).

Other naturally occurring materials, such as protein in hair samples or microorganisms are studied in crime laboratories using analytical pyrolysis (1). The need for glass-lined systems is demonstrated in the pyrograms of human hair, raw wool, microorganisms, silk, cellulose, and chitin shown in Figures 11 (a-f). Naturally-occurring materials are often more thermally/catalytically sensitive and require non-reactive surfaces to ensure high reproducibility. It is thus concluded, analytical pyrolysis technology has emerged as a powerful means to study diverse materials over wide concentration and thermal ranges for industrial R&D, and the materials, environmental and forensic sciences.

3.3 Supercritical Fluid Technology.

In addition to using thermal methods (e.g., DHS, pyrolysis) to assist in characterizing complex polymeric materials, nonthermal methods may be necessary for difficult samples, such as microcapsules with thermally sensitive core targets. The supercritical fluid (SF) treatments that can be used involve extraction (SFE), desorption (SFD), reaction (SFR), and on-line SF chromatography (SFC) or capillary GC. Background information in this field is given in Bibliography 10.2 (page 40).

Since supercritical fluid (SF) densities resemble liquid densities, high solvating power can be achieved with SF, yet the diffusivity/viscosity property of the fluids are gas-like for fast mass-transfer behavior. Temperature influences solubilities of solutes in ways that depend on the involved region of the solubility diagram. The general statement may be made that "the dissolving power of the SF is essentially proportional to its density". Therefore, significant changes in solvent properties are brought about by variation in temperature and pressure. Instrumentation with highly controlled temperature and pressure (density) parameters was used in these analytical applications. A beta-site arrangement with Computer Chemical Systems, Avondale, PA, made possible the development of this technology for microcapsule analysis at CRDEC.

For all analyses, the supercritical mobile fluid employed was CO₂ (critical temp = 31°C, critical pressure = 1070 psi), since it has properties that permit safe use at moderate temperatures and pressures (1100 psi to 6000 psi). Therefore, SFE or SFD analyses were conducted with complex materials, including consumer products, microcapsules, charcoals, soils, etc. Samples were placed into the extraction/desorber unit (up to ca. 400 mg of sample) and treated under SF conditions, a typical set of conditions being 3000 psi, 100°C, 15 min with mobile fluid CO₂. The effluent containing the

solvated species is transferred automatically to the on-line chromatographic unit for separation in packed or microbore columns. Simultaneously, the SFE effluent is monitored by a flame ionization detector (FID) to establish extraction profiles of the SF process. Detectors other than those used with GC instruments may also be on-line with the SF unit such as HPLC-type detectors, including ultra-violet (UV), visible (VIS), or spectral systems (MS, FTIR, or MS/MS). Exploratory research described in Section 6 (page 32) describes activities in fiber optic monitors (FOM) for the SFE unit throughout various spectral regions. This permits real-time, on-line detection and identification of SFE effluents (14 b, c).

4. EXPERIMENTAL DESIGN AND APPLIED ARTIFICIAL INTELLIGENCE

4.1 Experimental Design and Method Development.

Development of specific analytical methods proceeded on the pyrolysis and SF systems, which in two cases used an experimental design for assistance in method parameter screening, or defining a "response surface". The chosen design for this purpose was a five-factor, three-level Box-Behnken program (15). Extension was made of this well-known approach (which uses a single response) to a multivariate response because of the unique needs of our analytical data. The 46-run design used for analytical pyrolysis and SFE method development is shown in Figures 12 and 13. Patterns from GC, MS, or FTIR data were the needed "responses" upon which our interpretations were based and, therefore, the object of our attention. Method development is designed to provide the most information in an efficient manner, so generation of the "most informative" pattern is the goal in this situation. To quantitate these patterns, an extended Box-Behnken design was developed using a multivariate response derived from factor analysis and principle component analysis (16).

4.2 Applied Artificial Intelligence and Expert System Networks.

The concerted analysis of complex organic materials with modern instrumentation involves many judgemental decisions to be made in analytical strategy, instrumental configurations, data analysis and interpretation (Figure 14 (12)). Computer-assisted experimentation with microprocessor-based tools often includes embedded intelligence in the diagnostics, data acquisition, advanced data analysis and failsafe operations. However, need exists for the analyst to approach problem-solving for difficult materials such as microcapsules using the instrumental capabilities combined with some problem-solving expertise in the strategy and interpretive areas. An AI approach was developed that aids in the experimental determination of chemically significant information for these materials. A commercial shell, TIMM, (General Research Corp., McLean, VA) was used to develop the prototype network,

EXMAT, a linked network of expert systems for materials analysis (17). The prototype system also included the database management expert system, EXDBM, and the expert system for chemometrics EXMATH. The outline is given in Figure 15 (a-f) (17). AI-based work at CRDEC involved a test series of organic amine salts analyzed with analytical pyrolysis/concentrator/GC-MS instrumentation. Completion of the prototype development was shown by the use of the AI software and by chemometric analysis software to provide classification of the three groups of structurally similar organic amine salts (Figure 15 (e, f)) (18).

The EXMAT prototype network was extended for application to the microencapsulation field. This development included a decision structure outline for microencapsulated materials (MICROCAP) using an analytical strategy based on nine available and functional instrumental capabilities. In addition to the four major systems noted above, other analytical systems include FTIR, high pressure liquid chromatography (HPLC), thermal analysis (DSC, TGA, TMA), microscopy (SEM, TEM), X-ray fluorescence, and microchemical C,H,N,O elemental determinations.

Background in sensors and artificial intelligence topics is given in Bibliography 10.3 (page 41).

4.3 Technical Resources: Materials, Instrumentation, and Software.

4.3.1 Materials.

Necessary reference and control samples were obtained for method development within each instrumental field. Bend Corporation (Bend, OR) made available three different core pesticides with polymer shells/walls, solvent, emulsifier and surfactant. R.T. Dodge Corporation provided additional test samples with gelatin and isocyanurate cores inside variously sized polymer walls. This series is different from the Bend series not only in the type of shell and core material, but in the relative amounts; the former contained large amounts of target core species (50-80% by wt.), whereas the Dodge series had smaller amounts of different targets (2-5% by wt.). Selected test materials were used to establish analytical ranges of detection and sample loadings needed for the different instrumental techniques. The samples are representative of the materials to be analyzed and were available alone and as fully formulated products. Experimental designs require representative materials of varied process, treatment, or compositional parameters and the two sources provided the initial test materials.

Bend Research: Specific Sample Inventory: Microcapsules and Ingredients

Batch 1 - polycarbonate/Diazinon; Batch 2 - polysulfone/Diazinon; Batch 3 - polyetherimide (Ultem)/Diazinon;


Batch 4 - polysulfone/Cypermethrin; Batch 5 - polycarbonate/Dursban. Thus, the requested three polymer shells and three different pesticide cores were supplied along with the formulation ingredients: emulsifier A = gelatin; emulsifier B = (Pluronic)-polyol surfactant (31R2-BASF); solvent = methylene chloride; and water.

Bend core pesticides are: Diazinon, mol. wt. 304, decomposes above 120°C, Cypermethrin, mol. wt. 416, semisolid mp. 60-80°C, and Dursban mol. wt. 350, granular white solid, mp. 41-42°C. Their structures are shown below:

Pesticide Structures:


Diazinon

(Merck Index)

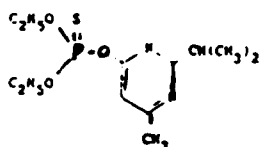

Liquid. Fast ester-like odor. d_4^{20} 1.116-1.118. n_D^{20} 1.4478-1.4481. Vapor pressure at 20° = 1.4×10^{-4} , at 40° = 1.1×10^{-3} (about 5 times vapor pressure of malathion). Volatility at 20° = 2.4 mg/cubic meter; at 40° = 17.6 mg/cubic meter. Decomposes above 120°. Sol. in water at 20° = 0.004%. Miscible with alcohol, ether, petrol, cyclohexane, benzene, and similar hydrocarbons.


Cypermethrin

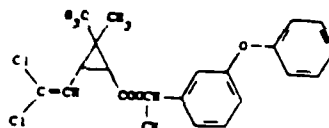

Commercial product is a mixture of eight isomers. Viscous semi-solid, mp 60-80°. Insol in water, sol in methanol, acetone, xylene, methylene dichloride.

USE Insecticide

Dursban

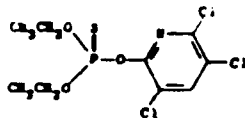

White granular crystals, mp 41-42°. Vapor press at 25° = 7×10^{-6} mm Hg. Sol. in water 2 ppm, insol in octane.

R.T. Dodge, Co.: Specific Sample Inventory:
Microcapsules and Ingredients: (20 gm)

18 samples: 2% gelatin/33.3% polyethersulfone (PES)/64.7% polysulfone; 5% gelatin/16.7% PES/78.3% (PS) polysulfone; 3% caffeine HCl/97% PS; 6% caffeine HCl/94% PS; 6% ACI 59 (Potassium Dichloroisocyanurate)/94% PS; 3% ACI 59/97% PS; 3% gelatin/97% PS; 6% gelatin/94% polyurethane (PU) (label); (3% caffeine HCl/97% PU); (3% ACI 59/97% PU); 6% ACI 59/94% PU; Gelatin; Caffeine HCl; ACI 59; PS; PES; PU.

The microcapsules contain core materials (bleach (ACI 59) or caffeine HCl); Wall materials - PS, PES or PU; Core sizes ranged by a factor of 3 to 10 (e.g., 20 microns/200 microns); wall thickness expressed in volume percent, 10 and 20%; isolation method was conducted with and without fumed silica.

In addition to these reference microencapsulated materials, numerous consumer products were available, such as the pesticide in the polymeric "No-Pest Strip", and over-the-counter medications as formulated in final products (Contac, Bufferin, etc.).

4.3.2 Instrumentation.

Four analytical instrumental systems were assembled to examine microencapsulated formulations:

System 1. Pyrolysis Atmospheric Pressure Chemical Ionization (APCI) Tandem Mass Spectrometer (Sciex TAGA 6000 unit).

System 2. Pyrolysis/Concentrator System (Pyroprobe Model 124) interfaced to (a) Hewlett-Packard 5985B GC/MS and to (b) Finnigan TSQ in MS and MS/MS modes. These modules allowed the following variations of operation: dynamic headspace, programmed and pulse pyrolysis, cryogenic and sorbent bed concentration/enrichment.

System 3. Pyrolysis GC with Fourier Transform Infrared Spectrometer (FTIR) (Hewlett-Packard IRD) at BRL.

System 4. Supercritical Fluid Chromatograph (CCS Model 5000) for extraction, desorption, reaction and chromatography.

Detection and identification were conducted using the instrumentation and protocols for development of analytical capability at CRDEC. These analytical systems combine sample processing (thermal and nonthermal) with concentrator and chromatographic units interfaced to spectral detection/identification systems. Also, data needed in the concerted analysis approach, included an elemental analysis capability. Conversion of an existing microchemical reaction unit (Chemical Data Systems Model 820) was completed for on-line C,H,N,O elemental

determinations on gas, liquid, or solids at microgram quantities. However, the recent announcement in 1989 by Hewlett-Packard of an atomic emission detector (AED) for GC effluents makes this unit obsolete. It is evident that elemental data are important for overall integration into the concerted analytical approach, whether achieved by a microchemical system or an AED unit.

The general schematic shown in Figure 4 outlines the experimental setup.

4.3.3 Software: Artificial Intelligence and Experimental Design.

Analytical strategy, data analysis, and interpretation were topics included in development of applied artificial intelligence (AI) capabilities at CRDEC (Figure 15a). The expert system network, EXMAT, is outlined in Figure 15 (b). One of the seven outlined knowledge bases was developed (ES#1, Figure 15c) for analytical strategy to include over two hundred rules. The other six were only roughly outlined and remain to be properly developed and/or modified into a functional expert systems network. The prototype development, however, provided invaluable experience in opening up this field to Army applications. These activities are discussed in detail in references 17 and 18. The testing of the database management expert system EXDBM (A.M. Harper, (18)) in the EXMAT network is illustrated (Figure 15 d) with data from the analytical pyrolysis study on a group of twenty samples of three structurally similar groups (discussed in Section 5.2.2). Classification was achieved and extended analysis (factor rotation and plotting of covariance data) are shown in Figure 15 (e, f).

Within the EXMAT network (Figure 16) is a specific expert system, MICROCAP, which includes a prototype decision structure and a few test rules in the knowledge base. The instrumentation needed to determine such complex materials and the accompanying method details are part of the main EXMAT network. The software includes descriptive commentary under "HELP" or "Verbose Version" options for each of the analytical systems developed for the CRDEC facilities, as well as those at Ballistics Research Laboratory (BRL) (pyrolysis/concentrator GC-MS and FTIR units). General information is grouped under "HELP" for analytical pyrolysis/concentrator and chromatographic/spectral units with varying levels of detail. For example, the TAGA 6000 triple quadrupole MS/MS unit includes specific settings for a "typical" run using the Pyroprobe Model 122 sample treatment interfaced to it (Figure 17). The documentation of this applied AI activity for development of CRDEC short- and long-term capability is given in Bibliography 10.3 (page 41), and 10.5.1 (page 44), Presentations.

5. RESULTS OF TRACE ORGANIC ANALYSIS METHOD DEVELOPMENT

5.1 Experimental Design.

Work reported from BRL (15) permitted adaptation of the Box-Behnken experimental design to the microencapsulation task. The Box-Behnken computer program was adapted for the Compaq computer following minor changes in the Pascal code. Use of the design was made in a screening mode as in the pyrolysis runs, rather than a strict statistical manner, because of experimental limitations at that time. Similarly, SFE method development was approached for the first time with this same Box-Behnken experimental design for a complete 46-run series (see Figures 18 and 19).

The 46-runs in this central composite design include six replicates of the centerpoint to evaluate experimental reproducibility. The use of this experimental design for both analytical pyrolysis and SFE method development brought out an important limitation in the approach as noted from the BRL work. The statistical data can be evaluated using only a single response, rather than the resultant GC, MS, or SFC PATTERN response that is informative for the method development. Extended study has pointed out the benefits of extended Box-Behnken designs, if use is made of a multivariate response that represents the pattern generated from the multiparameter study (16). This extended software, now operational on the CRDEC Compaq, is ready for application studies to determine its utility in general method development involving multivariate or pattern response, rather than the single response format.

5.2 Experimental Data and Discussion.

5.2.1 System 1. Py/APCI/MS/MS

The combination of the Pyroprobe Model 122/TAGA 6000 provided analytical pyrolysis (pulse and programmed to 1400°C) and detection/identification with atmospheric pressure chemical ionization MS/MS. Applications include detection and identification of volatiles in solids or semi-solids, volatile liquids and occluded solvents or degradation products from thermal treatments.

The APCI-MS/MS TAGA unit received volatiles from a pyrolyzer unit, Model 122, which provides only thermal processing, but no enrichment/concentration capability. The inert N₂ or air atmospheres were provided by modification of the Pyroprobe coil insert when interfaced to the APCI MS/MS unit, as discussed elsewhere (19), and shown in Figure 20(a).

In context to a general method for determination of polymer microstructure, decomposition products or low-level thermal

treatments applied to various materials, the following polymers were studied with pyrolysis APCI MS/MS and examined in the single quadrupole mode: poly(methylmethacrylate) (PMMA), poly(ethyleneimine) copolymer, poly(vinylchloride) (PVC), formulated PVC (Shell No-Pest Strip), polyethylene - high density (HDPE), hydroxy-terminated polybutadiene (R-45M) (HTPB), formulated R-45M (HTPB-1, -2, -3), carbon fibers (pitch, polyacrylonitrile-carbonized PAN), nylon 6/6 poly(hexamethylenedipamide), nylon 6/9 poly(hexamethylenenonanediamide), nylon 6/10 poly(hexamethylenedipamide sebacamide), oxidized PE, polypropylene isotactic (iso-PP), glycogen, polyethylene terephthalate (PET), sodium alginate, PE-oxidized 2, Pellethane (Upjohn) polyurethane, Estane 5703 (B.F. Goodrich) polyurethane, ethyl cellulose, polyacrylamide (high carboxyl-modified and low carboxyl-modified), polyethylene oxide (PEO), poly(butyl methacrylate) (PBMA), ethylene-vinyl acetate copolymer (18% VA) (PE-VA), butyl methacrylate-isobutyl methacrylate copolymer (PBMA-IBMA), polyacrylic acid, polyacrylamide.

Pyrograms shown in Figure 20(b) demonstrate the value of comparing inert (N_2) vs. oxidative (air) pyrolyses for a PE; an oxidized PE, and an "in-chain" oxygen polymer PEO (CH_2CH_2-O)_n. The APCI MS/MS has excellent sensitivity to oxidative products produced in these studies and is a preferred method for such analyses.

The microcap #1 sample and the three polymer shells A, B, and C from Bend Research were studied with this system under a range of pyrolysis conditions, pulse and programmed, and appropriate reference (blank) conditions (240°C/min to 700°C) under air flow of ca. 150 ml/min with the interface configured as in Figure 20a. Figure 21 (a, b) gives the observed patterns for microcap #1 under (1) 180°C/min to 400°C control run with no heat (shows no molecular ion for diazinon at m/z 304) (2) 180°C/min to 400°C (with heat) for a comparison of volatiles detected under identical conditions of flow, etc. to those with a specific low-level DHS-like treatment (although no enrichment is possible in this configuration). The 18 m/z dalton loss in Figure 21b is likely the H₂O loss from the microencapsulated formulation and could not be obtained in other systems with FID detection (either PGC-FID or SFE-SFC). Other experiments include: (3) pulsing the same sample from (2)-treatment to 400°C/20 sec to give an excellent m/z 304 molecular ion peak for the core Diazinon and an m/z 153 minor peak; a second pulse at 400°C/20 sec on the same sample from (3)-treatment did not give a significant peak at m/z 304, but peaks at m/z 153, 127, and 99. A pulse to 375°C/40 sec on microcap #1 gave an easily distinguished m/z 304 peak, a minor m/z 121, as did pulsing to 400°C/20 sec. A pulse to 700°C/20 sec (Figure 21c) gave no significant m/z 304 ion, but a major fragment was observed at m/z 153 and a minor peak at m/z 137. Programming at 240°C/min to 500°C (Figure 21d) gave a minor m/z 304, strong m/z 153, and minor m/z 127, not unlike the 700°C/20 sec pulse. Programming 300°C/min to 200°C (Figure 21e) (4 min interval) gave a moderate m/z 304, and

m/z 288, and major m/z 153, 139 and 111 fragments. This last run compares the fast rate (300°C/min) to get a relatively low (200°C) final temperature to the pulse treatment of 400°C in other runs not shown) which yield a predominant m/z 304 core peak. Several in the series are given in Figures 21 (a-e).

This series of runs would conclude that the pulse at 375°C/20 sec or 400°C/20 sec is most informative as to the detection of the m/z 304 molecular ion for Diazinon and the shell polymer fragments. The polymer shell fragments are seen separately in Figure 22 when the polycarbonate shell was programmed to 400°C at 180°C/min. Hence, both shell and core were detected by their MS fragment patterns when thermally treated. A further comparison is shown in Figure 23 (a, b, c) for microcaps #1, 2, and 3 all with the Diazinon core, pulsed at 400°C/20 sec to give a strong m/z 304 m/z and less intense fragments at m/z 288 and 152/153. Figure 23 (d) presents a microencapsulated protein and its reference run to demonstrate MS pattern comparisons for protein detection in microencapsulation (19b).

5.2.2 System 2. Pyrolysis-Concentrator GC-MS, and -MS/MS

Pyrolysis/Concentrator Pyroprobe 124 with Hewlett-Packard 5985B GC-MS and with Finnigan TSQ in MS and MS/MS modes were used, including DHS, pulse and programmed pyrolysis, cryogenic and/or sorbent bed concentration/enrichment. Sample processing can be conducted under inert or reactive atmospheres, with either a coil or ribbon pyrolyzer sampling system or thermal desorption module for aerosol samples. Not included in this study was the demonstrated use of the purge and trap module for analysis of dissolved targets in water/liquids or from bulky solids placed in glass desorption vessels up to 50 ml size.

Figure 24 a is the experimental schematic. Figure 24b and c illustrate the results from conducting a DHS/pyrolysis GC-MS experiment with the Pyroprobe Model 124 interfaced to the Hewlett-Packard GC-MS on a series of twenty organic amine salts from three structurally similar groups and relevant referenced compounds (Figure 25).

The same Pyroprobe 124 system was then interfaced to the Finnigan TSQ GC/MS unit. Microcapsule components were examined for their behavior in the new configuration. Control runs were made on a reference polymer HDPE, as well as oxidized PE and short-branched PE to compare "standard" PGC pyrograms on the capillary DB-5 30 meter x .33mm fused silica column. Figures 26 (a-f) show the excellent triplet pattern and short-branch and oxidation information that can be obtained with the TSQ unit on the reference polymer (20). More detailed pyrograms are illustrated in Figures 26 (b, d, and f) for the region between 650 and 950. Figure 27 records the experimental conditions under which the GC-MS spectra were obtained.

Selected microencapsulated materials from the Bend and Dodge Co. series were examined over a range of pyrolysis conditions. The core ingredient, Diazinon (Bend Corp.), was examined under a $240^{\circ}\text{C}/\text{min}$ to 700°C final temperature with the thermal desorber at 150°C to give the RIC seen in Figure 28(a). EI/MS examination of the reconstructed ion chromatogram (RIC) for scanset 672 is shown in Figure 28(b) and served to identify the pesticide as Diazinon. Its definitive EI-MS pattern consists of the molecular ion at m/z 304 and fragments at 276, 179, 152, 137, etc. Figure 29 (a, b) shows the detection of Diazinon in a shell of polycarbonate (microcap #1) examined under the same conditions as Figure 28. Several combinations of pulsed/programmed pyrolysis were examined for behavior of the core/shell degradation patterns, e.g., compare microcap #1 run at $180^{\circ}\text{C}/\text{min}$ to 600°C in an 8-min interval (Figure 30 a, b) to that obtained at $240^{\circ}\text{C}/\text{min}$ to 700°C in a 4 min interval (Figure 30c). Both runs gave easy definition of the pesticide core from the polymer/ingredients fragments in scansets 409, 437, 455, 479, 481, 515, and 581 regions. The intensity of the core (scanset #668) in Figure 30c in the faster $240^{\circ}\text{C}/\text{min}$ to 700°C indicated it to be the preferred programmed pyrolysis condition for maximum Diazinon release with minimum thermal degradation. Separate runs under pulse conditions ($500^{\circ}\text{C}/20\text{sec}$) gave the RIC and MS patterns for Diazinon and the MS/MS run of m/z 304 ion seen in Figure 31 (a, b, c).

Under the pulse mode of $500^{\circ}\text{C}/20$ sec, the Bend Microcap #4 (Cypermethrin) gave the patterns in Figure 32 (a, b). The EI/MS pattern from scanset #608 yields fragments at m/z 208, 173, 163, 122, 109, 91, etc. Programmed pyrolysis ($240^{\circ}\text{C}/\text{min}$ to 400°C) of the Microcap #4 gave scanset #395 which shows an EI/MS pattern of m/z 208, 184, 173, 163, 127, 109, 91, etc. (Figure 33 a,b). Microcap #5 containing Dursban in a polycarbonate shell gave the RIC shown in Figure 34a. The EI/MS pattern of scanset #831 is shown in Figure 34b. The MS/MS run of the m/z 199 parent MS fragment (Figure 34c) yields a m/z 171 and minor m/z 163, 100 fragments. The Dursban core differs from Diazinon by the polychlorinated pyridine ring linked to the phosphorothioate portion, whereas Diazinon has the alkylated pyrimidyl thiophosphate structure (see page 22).

The gelatin emulsifier (5%) in a polysulfone-polyethersulfone shell (Dodge) gave an RIC with scansets at 168, 202, 320 and others (Figure 35a). Gelatin from different sources gave RIC results seen in Figure 35 (b), (c), and (d) with comparison of the Bend gelatin, P.S. and Dodge (500 micron size range) gelatins. The P.S. gelatin RIC (Figure 35c) gave scansets 196, 206, 258, and 301 when treated at $240^{\circ}/\text{min}$ to 700°C , 8 min. interval, while the Dodge gelatin (Figure 35d) gave scansets 191, 276, 345, 387. The Dodge gelatin desorbed 5 min at 260°C gave the RIC of Figure 35e. The EI-MS of scanset 206 yields fragments at 207, 191, 133, etc. (Figure 35f).

5.2.3 System 3. Pyrolysis/Concentrator GC-FTIR

System 3 is a Pyroprobe Model 122 coupled to a GC with a Fourier transform infrared spectrometer (FTIR). A representative total response chromatogram (TRC) is shown (Figure 36a) for the Bend microcapsule #1 and the associated FTIR pattern (Figure 36b) (also shown are the library search results). Since Diazinon was not present in the library used in this search (EPA/Sadtler), a direct match was not obtained, nor was a separate reference run conducted to add it to the library. Importantly, however, the major functionality (phosphorothioic acid derivative) was identified. The significance of this complementary data to the MS spectral pattern cannot be overstated. Identification of "true" unknowns requires the combined information from GC-MS/GC-FTIR instrumentation.

5.2.4 System 4. Supercritical Fluid Instrumentation

5.2.4.1 SFE. Supercritical Fluid Extraction

Supercritical Fluid Chromatograph Model 5000 SFC/GC is configured for extraction, desorption, reaction, and/or chromatography. The initial studies with SFE and SFR involved a polymer matrix containing many propellant formulation ingredients and a sensitive component, RDX, as the target species, albeit not in a microencapsulated form. Figure 37 illustrates the resultant SFE profiles of formulated binders 1, 2, 3 with species extracted under non-optimized conditions. Automatic transfer of the extractables into the on-line SFC unit provided the separation of components in a standard microbore chromatographic column. Many feasibility runs on components in polymeric materials (polyurethanes, polyethylenes, etc.) demonstrate the potential for full SFE-SFC analyses for temperature sensitive substances in microcapsule or immobilized form. Figure 38 presents the analogous SFE-SFC analyses for two samples from the same batch of formulated microcap #1.

The Bend Research microcap series of known compositional variations provided the background for parameterization of the SFE method. Within the experimental limits of the instrument, the task was to determine the "optimum" combinations of pressure, temperature, density, flow, time, etc. to yield the most comprehensive information on such complex samples. Hence, the design was used in a non-statistical "response surface" screening mode.

Results from this series are non-statistical because the response of concern is an extraction profile, not a standard chromatographic or spectral pattern. Figure 39 (a-c) gives a range of typical profiles obtained in the Box-Behnken series. The information content resides in the core target extraction

"efficiency" and the pesticide core separation, as well as associated formulation ingredients. This screening or experimental response surface showed that the SFE method permits wide ranges of pressure, temperature, flows, and residence times for satisfactory detection/analysis of these representative core targets in the selected polymeric shells.

Since the payload or core of the test microcaps was high (~ 50-80% by wt) and only a few milligrams of sample were used in each run, it is safe to say that a typical SFE experiment for these types of studies would be conducted with microgram amounts of sample under nominal conditions: i.e., 1200 psi initial CO₂ mobile fluid, 8 min. extraction time at 100°C, and pressure programming at 200 psi/min to 4000 psi (Figure 39b). Presence of emulsifiers/surfactants are not necessarily detected under these conditions, and therefore, varied conditions are needed to ensure a more complete analysis of all the formulation ingredients. Perhaps a two-step pretreat/extraction or reactive extraction may be used for this purpose, not unlike the route taken in analytical pyrolysis method development: i.e., dynamic headspace thermal preprocessing of volatiles, followed by pulse or programmed pyrolysis.

Since SFE is only one step in the detection of a core target in a polymeric or complex matrix, the identification step must involve more than a proper extraction time. The automatic SFE-SFC system provides chromatographic separation and detection systems for more complete characterization.

Figure 39a shows the SFE profiles of microcapsule #1 from three runs in the experimental design. Five samples of microcapsules examined under a single set of "reasonable" (as seen in the Box-Behnken) experimental settings are given in Figure 39b. Microcaps 1, 2, and 3 contain Diazinon within three different shells...polycarbonate, polysulfone, and polyetherimide (Ultem), respectively. Microcap 4 contains Cypermethrin inside a polysulfone shell and Microcap 5, Dursban inside a polycarbonate shell. Each of the three pesticides are readily extracted under these conditions, but the profiles differ according to the shell/formulation content (Figures 39c-e). Variability between microcapsule formulation or relative ease of extraction through the different shells cannot be determined from this series. Repeat runs (see Figures 39(d) and (e)) with a different sample from the same batch illustrate this point...in-group variability can be monitored and compared to profiles from different batches or to different shell/core formulations. The data only serve to demonstrate that the target cores and some formulation ingredients are readily analyzed with SFE. Also, the three selected target cores are easily differentiated by their extraction behavior in this trial series.

Furthermore, Figure 40 shows SFE-SFC data for gelatin (emulsifier A) (used in the Bend microcapsule formulation) and a typical polyol emulsifier/surfactant (emulsifier B). The complex profile with several individual peaks seem to indicate solubilization of at least some of the gelatin material. The proteinaceous gelatin serves as a reference for developing the analytical SFE-SFC method.

These microcapsules of known compositional variations thus provide the background for the parameterization of the SFE method. Within the experimental limits of the instrument, the task remains to determine what are the combinations of pressure, temperature, density, flow, time, etc. to provide the most comprehensive and quantitative information on such complex materials that may vary over wide ranges of composition.

5.2.4.2 SFR. Supercritical Fluid Reaction

Extension to SFR was also shown with the RDX formulated system, since research on RDX decomposition mechanisms had resulted in significant new information in an important area of research (19). The exploratory runs with RDX decomposition under CO₂ SF conditions and with superimposed thermal input (up to ca. 250°C) in the extractor gave SFR results shown in Figure 41. These runs demonstrated that monitoring/analyzing decomposition products over a range of extended reaction conditions were possible and enabled comparisons to be made for catalytic vs. noncatalytic decompositions. These RDX mechanism studies demonstrated extended capability for thermally sensitive and reactive targets to be analyzed over a range of conditions, including high pressure regimes not previously available with commercial analytical systems with such automated capabilities. Extension to "reactive extraction" of different biologicals, toxins, etc. is an important direction to follow.

5.2.4.3 SFD. Supercritical Fluid Desorption

A simulant, bis-(ethylhexyl)phosphonate, was examined by (SFE) SFD/SFC in trial runs. Results are shown in Figure 42a indicating successful analysis is possible when the simulant is in a soil matrix at ppm levels. Extrapolation of the data to expected detection levels in a method development mode indicates detection in the nanogram range with many adjustable parameters to extend it even further. Hence, targets such as the simulant phosphonate appears readily extracted/desorbed at trace detection levels. Identification again would require more than the retention time in either SFC or GC on-line analysis. An important illustration of SFD is given in Figure 42(b) showing different charcoal samples analyzed in the system. The on-line capillary GC analyses are preferred for these studies, since the complex hydrocarbon adsorbents are readily separated using typical environmental analysis protocols; i.e., spectral identification by SF interfaced

to an FTIR or MS unit. The fiber optic feasibility studies would be particularly useful in this regard...use of a fiber optic monitor (FOM) in the mid-IR region for identification using the IR pattern and/or information on chemical functionality (see Section 6.1)

5.3 Supercritical Fluid Extraction-Chromatography SFE/SFC.

SFE-SFC configurations, as well as SFE-, SFD-, or SFR-capillary GC are important integrated analytical tools for detection and identification. SF extraction with minimal thermal input is needed for the target species that are thermally labile followed by on-line automated separation (SFC or capillary GC). These analytical steps have been demonstrated for microcapsules (Figure 38) in this series. Integration with spectral detection (MS and/or FTIR) and enrichment or concentration stages provide the means for identification of trace targets or classification of unknowns with the aid of chemometric approaches. The newly emerging fields of neural network/applied AI and fiber optic sensors are directly relevant to a capability for in-field, portable analytical monitors of target species in complex environmental scenarios. Results from the microEXMAT development (the expert system network for supercritical fluid technologies) are reported in context to the importance of applied AI to analytical chemical instrumentation with "embedded" intelligence. General conclusions of the SF technology series with microencapsulated pesticides are given in Figure 42(c).

6. EXPLORATORY RESEARCH AND TECHNOLOGY ASSESSMENT

6.1 Fiber Optics.

Exploratory research was conducted with fiber optic monitoring of supercritical fluid extractables in the UV, VIS, NIR, and mid-IR regions using probe/cables with suitable spectrometers. These latter units included (1) Guided Wave's Model 100 spectrometer (CRDEC); (2) Mattson's FTIR spectrometer (R. Fifer, BRL); and (3) Infrared Fiber Systems' potential use of a solid state spectrometer based on an acousto-optic tunable filter (AOTF). The feasibility studies were documented (14b, c) and included special configurations that were tested for the SFE-FOM interface (Figure 43). Feasibility runs on the system for detection of trace aromatics are shown in Figure (44 a, b). Extension to the mid-IR is shown in Figure 45. The summary/conclusions of this fiber optic development are given in Figure 46 (a and b).

6.2 AI and Neural Networks.

Extension of EXMAT to microEXMAT (an expert system for supercritical fluid technologies) was documented (21) and is

illustrated in Figure 47. Figure 47 (a) shows the "EXCONFIG" outline and 47 (b) the outline of "EXMETH". These are two of the initial expert systems in the network.

Technology assessment of the AI field indicated that neural networks (noted in the 1987 literature) have realistic potential for applications of particular significance in pattern recognition. The hardware and software were tied into sensor systems and adaptive training algorithms for classification and prediction. Decisions are directly analogous to those developed using the AI software of TIMM. Neural networks, as is TIMM, are used on microcomputers, not just mainframes such as the VAX 8600 or MicroVAX II computers that run the existing TIMM software. This fast-paced development within the AI/computer industry must be followed closely to make use of relevant technology for the CRDEC mission; particularly in the areas of pattern recognition and information fusion from selected chemical sensors (see Bibliography 10.3 on Sensors/AI, (page 41), (see ref. 42).

6.3 Inverse GC and SFC.

Inverse GC has been shown to be a useful tool for the study of surfaces and subsurfaces of polymers (1). The method involves use of GC/SFC instrumentation for close control of time, temperature, flows, and sensitive detection units just as in normal GC, but the nature of the column function is changed. Rather than inject unknown liquids/gases into the separation column so that they elute according to the chromatographic process, injection of known "molecular probes" is made into a chromatographic column containing the substrate film, coated particulates, or other solid materials (fibers, etc.) which are to be studied. Behavior of the selected molecular probes in the mobile phase as they interact with the stationary phase in a classic chromatographic separation experiment provides detailed information on thermodynamics of mixing and interaction; e.g., free energy, enthalpy and entropy of adsorption, partial molar free energy, enthalpy and entropy of mixing, solubility parameters, heats of solution, diffusion coefficients, surface concentration at monolayer coverage, chemical potentials, and other data not easily obtained by other analytical methods for polymeric materials.

Extension to the SF field provides a new opportunity to examine fibers, surfaces, and subsurfaces of materials with the added dimension of high pressure and unique properties of SF mobile phases. One potential application involves SF technology applied to preparation and analysis of specialty treated charcoals. Solubilizing oligomeric/polymer coatings and depositing them upon desired charcoal or solid surfaces presents a new way of preparing coated high surface-area materials. Evaluation of these coated materials then is achieved by challenging the treated samples in-situ with the appropriate vapors and measurement of

classic breakthrough/quantitation information. Use of the SFD mode as discussed above for analysis of charcoal traps from air pollution studies involved desorption from charcoal placed in the extractor/desorber accessory module. The inverse GC or SFC mode may or may not use the accessory, but uses the chromatographic portion of the system most appropriate for the nature of the information needed.

7. SUMMARY AND CONCLUSIONS

The intended task to establish at CRDEC the capability to perform trace organic analysis in the field of microencapsulated materials and complex polymeric materials of concern in threat analysis has been outlined and in part completed. The selected instrumental tools and methods emphasizing analytical pyrolysis and supercritical fluid technologies were demonstrated and made operational at CRDEC for a variety of test materials. These concerted analyses on microencapsulated samples successfully showed that trace detection and identification can be carried out for target cores in a range of polymer matrices and for targets that are immobilized in soils, adsorbed on charcoal surfaces or other particulate carriers/solid matrices.

During the initial work, only qualitative and semi-quantitative analyses were conducted. It is necessary to have a full quantitative study on a selected system that would demonstrate the necessary detailed methodology in terms of the reproducibility at the sampling stage, data generation (pyrolysis or extraction), and in data analysis/interpretation. Extension of the applied AI decision structure and initial input to the knowledge base in this specific domain was made in context to the EXMAT network. Extension to initial assessment of neural networks was undertaken in relation to pattern recognition and multi-sensor fusion which have the potential for real-time detection and identification.

8. RECOMMENDATIONS

Existing instrumentation and method development should be used and extended to establish full qualitative and quantitative determination of trace amounts of microencapsulated materials important in threat analysis. This must be combined with appropriate sampling common to multimedia environmental protocols. Pattern recognition algorithms with some AI treatment should be developed from the background in AI and multivariate data analysis that exists from the CRDEC/BRL studies. The prototype EXMAT requires an improved, current decision structure, as well as development of the associated knowledge bases on the microVAXII system. Of the existing seven linked expert systems, only one

("strategy") has a sizeable knowledge base (250 rules). Persons knowledgeable (i.e., experts) in environmental and trace organic analysis could develop the needed decision structure/knowledge bases with reasonable effort (approximately 8-10 months).

Development of microEXMAT illustrated the capability of a PC-level shell for a focused field (supercritical fluid technology). The associated chemometric software developed at CRDEC can thus be used in conjunction with the GRC-AI shell (microTIMM) if purchased for this purpose. In this case, development of the chemometrics/AI expert system would require not only expertise in the analytical chemistry applications, but specific capability for efficient microTIMM (Fortran 77)/chemometric software integration. General Research Corporation has demonstrated such embedded intelligence software for a range of multi-sensor applications (Figure 47c). The down-sizing of computer hardware/software now used for multi-sensor data treatments and pattern recognition/classification studies fits directly into the threat analysis area in a very realistic manner.

Importantly, extension of these laboratory units/methods to microtechnology in manufacturing would make possible a modular "universal organic threat analyzer" in a highly automated, miniaturized form. Sampling modules for air, liquid, or solids can be integrated into the unit and thermal/nonthermal sample processing would provide "extractables" automatically delivered to the multidetection/identification unit, with or without the immediate concentration or separation. The resultant set of descriptor patterns from the key GC, IR, and MS detection units for the selected target species would be combined with elemental data and provide input to a knowledge base trained for the unit. Expertise in handling a range of experimental and interpretive decision-making stages would be included in the applied AI format, using advantages of neural network advances. This potential system would be designed as an "on-site", in-field portable threat analyzer for screening a selected range of targeted materials, but also providing information on unknowns not in the library "trained" for detection of specific species. Knowledge base extensions would be made accordingly. The key datasets would provide the most efficient, reliable, and versatile screening approach to handle difficult unknowns, including threat agents formulated within microencapsulated materials in multimedia.

Thus, results from microcapsule method development combined with exploratory research efforts in fiber optics, applied AI and neural networks provide information needed to design and fabricate a modular portable threat analyzer for military, forensic, and environmental uses. This combined user-manufacturer development is one of the most productive ways for technology transfer to take place (22).

9. REFERENCES

1. a. Pyrolysis and GC in Polymer Analysis, Liebman and Levy (Eds.), M. Dekker, NY 1985.
b. Munson, T.O., Vick, J., J. Anal. Appl. Pyrolysis, 8, 493-501 (1985).
2. Montgomery, D.C., Design and Analysis of Experiments, J. Wiley and Sons, NY (1976, 1984)
3. Dodge, Y., Analysis of Experiments with Missing Data, J. Wiley & Sons, NY (1985) in Wiley Series in Probability and Mathematical Statistics.
4. Dowdy, S. and Wearden, S., Statistics for Research, J. Wiley and Sons, NY (1983) in Wiley Series of Probability and Mathematical Statistics.
5. Malinoski, E.R., Howerly, D.G., Factal Analysis in Chemistry, J. Wiley-Interscience, NY (1980).
6. Sharaf, M.A., Illman, D.L., Kolwaski, B., Chemometrics, Vol. 82, Chemical Analysis, John Wiley and Sons, NY (1986).
7. Box, G.E.P. and Draper, N.R., Empirical Model-Building and Response Surfaces, J. Wiley & Sons, NY (1987).
8. Strouf, O., Chemical Pattern Recognition, Research Studies Press, J. Wiley & Sons, England, 1986.
9. Microencapsulation Course, Center for Professional Advancement, NJ, Robert E. Sparks, Course Director, March 10-12, 1986, Chicago, IL.
10. Dialog Search 399, Microencapsulation, D. Paterno, E.W. Sarver, D.J. Reutter, CRDEC Library.
11. a. Liebman, S.A., et al., J. Polym. Sci., Part A-1, 9, p. 1823, 1921 (1971).
b. Liebman, S.A., Levy, E.J., NBS Spec. Pub. No. 519 (1979), Hertz and Chester (Eds.), U.S. Govt. Printing Office, Washington, D.C., p. 753.
12. Liebman, S.A., Amer. Lab., pp 17-19, Dec. 1971.
13. a. Wampler, T.P., Bowe, W.A., Higgins, J. and Levy, E.J., Amer. Lab., August 1985 and J. Chromat. Sci. 23, 64-67 (1985).

b. Gargus, A.G., Liebman, S.A., Watterson, C.,
Amer. Lab., pp 114-120, May 1981.

14. a. Levy, E.J., et al. Amer. Lab, 19, 66 (1987).

b. Liebman, S.A., Levy, E.J., Lurcott, S.,
O'Neill, S., Guthrie, J., Ryan, T., and Yocklovich, S.,
J. Chromat. Sci, 27, 118-126 (1989).

c. Liebman, S.A., Fifer, R.A., Griffiths, P.,
Lurcott, S., Bergman, B. and Levy, E.J., Pittsburgh Conference,
Atlanta, GA, 1989, paper no. 1545.

15. Liebman, S.A., Duff, P.J., Fickie, K.D., Schroeder,
M.A., and Fifer, R.A. J. Hazardous Materials, 13 51-56 (1986).

16. Snyder, A.P., Windig, W. and Liebman, S.A., work to
be published... Also, earlier report, Liebman, S.A., Ahlstrom, D.H.,
Hoke, A.T., Chromatographia, 11, 427 (1978).

17. Liebman, S.A., in Applications of AI in Analytical
Chemistry, ACS Sympos. Series 306, ACS, Washington, DC, 1986, p
365-384.

18. Wasserman, M.B., Reutter, D.J., Snyder, A.P.,
Leonard, J.W., Harper, A.M., Liebman, S.A. and Windig, W., 1987
CRDEC Conference on Chemical Defense Research, November 1987, and
Pittsburgh Conference, New Orleans, LA, 1988.

19. a. Liebman, S.A., Snyder, A.P., Kremer, J.H.,
Reutter, D.J., Schroeder, M.A., and Fifer, R.A., J. Anal. Appl.
Pyrolysis, 12, 83-95 (1987).

b. Snyder, A.P. and Windig, W., CRDEC-TR-011,
February 1989.

20. Liebman, S.A., Ahlstrom, D.H., Starnes, W.H., Jr.
and Schilling, F.C., J. Macromol. Sci.-Chem, A17(6) pp 935-950
(1982).

21. Liebman, S.A., Ryan, T., Yocklovich, S.,
Lurcott, S., Watkins, J., and Levy, E.J. The Pittsburgh Conference
and Exposition, March 1989, paper no. 109.

22. Specific analytical instrumentation developed in
this manner include the manufactured products from (a) Chemical
Data Systems, Oxford, PA, for the analytical pyrolysis sample
concentrators, micropilot plant reactor systems, and (b) Computer
Chemical Systems, Avondale, PA, for supercritical fluid
technologies. All commercially available products are
microprocessor-based units developed over the past fifteen years.

10. BIBLIOGRAPHY

10.1 Bibliography on Microencapsulation Research and Technology.

1. (a) Microencapsulation Course, Center for Professional Advancement, NJ, Robert E. Sparks, Course Director, March 10-12, 1986, Chicago, IL.
(b) Dialog Search 399, Microencapsulation, D. Paterno, E.W. Sarver, D.J. Reutter.
2. Reutter, D.J., Shuely, W.J., and Liebman, S.A., Workshop on Macromolecules in Microencapsulation Research and Technology, December 3-4, 1986, CRDEC-SP-87022, U.S. Army Chemical Research, Development and Engineering Center, Aberdeen Proving Ground, MD.
3. Short Course, "How to Make Microcapsules", C. Thies, Chem. Engin. Dept., Washington Univ., St. Louis, MO.
4. P. Guiot, P. Couvreur, Polymeric Nanoparticles and Microspheres, CRC Press, Inc., Boca Raton, FL, 216 pp, 1986.
5. L. Krowczynski, Extended-Release Dosage Forms, CRC Press, Inc., Boca Raton, FL, 1987.
6. K. Shinoda, S. Friberg, Emulsions and Solubilization, Wiley & Sons, Inc., NJ, 186 pp, 1986.
7. D.S.T. Hsieh, Controlled Release Systems: Fabrication Technology, CRC Press, Inc., Boca Raton, FL, 200 pp, 1988.
8. P. Molyneux, Water-Soluble Synthetic Polymers: Properties and Behavior, Vol. 1, CRC Press, Inc., Boca Raton, FL, 1986.
9. S.E. Friberg, P. Bothorel, Microemulsions: Structure and Dynamics, CRC Press, Inc., Boca Raton, FL, 232 pp, 1987.
10. Kondo, A., Van Valkenburg, J.W., Microcapsule Processing and Technology, Marcel Dekker, Inc., 192 pp, 1979.
11. Gutcho, M.H., "Microcapsules and other Capsules 1979 (Advances Since 1975)", Noyes Data Corp., Park Ridge, NJ, 340 pp, 1979, Chem. Tech. Review No. 135.
12. Deasy, P.B., Microencapsulation and Related Drug Processes, M. Dekker, NY, 1984.

13. Microencapsulation -- International Symposium, J.R. Nixon, Ed., M. Dekker, NY, 1976.
14. Kricka, L.J., Thorpe, G.H.G., "Immobilized Enzymes in Analysis", Trends in Biotechnology, Vol. 4, No. 10, Oct 1986, Elsevier Sci. Publ., NY.
15. Polymer Surfaces and Interfaces, W.J. Feast, H.S. Munro (Eds.), J. Wiley & Sons, NY, 1987.
16. "Ordered Media in Chemical Separations," W.L. Hinze, W.L., D.W. Armstrong, Eds., ACS Sympos. Series No. 342, Amer. Chem. Soc., Washington, DC, 1987.
17. Wasan, D.T., Ginn, M.E., Shah, D.O., Surfactants in Chemical/Process Engineering, M. Dekker, NY, 544 pp, 1988.
18. Bailey, J.E., Ollis, D.F., Biochemical Engineering Fundamentals, McGraw-Hill Book Co., NY, 1986.
19. Kesting, R.E., Synthetic Polymeric Membranes, J. Wiley & Sons, NY, 1985.
20. Macromolecules as Drugs and as Carriers for Biologically Active Materials, D.A. Tirrell, L.G. Donaruma, A.B. Turek, Eds., Annals NY Academy Sci., NY Academy Sci., NY, 1985.
21. Starzak, M.E., Physical Chemistry of Membranes, Academic Press, NY, 1984.
22. R.S. Langer, and Wise, D.L., Medical Applications of Controlled Release, Vol. 1, CRC Press, Inc., Boca Raton, FL, 1984.
23. Effect of Polymers on Dispersion Properties, Th.F. Tadros, Ed., Academic Press, London, 1984.
24. Baker, R., Controlled Release of Biologically Active Agents, J. Wiley & Sons, NY, 1987.
25. Biological Applications to the Controlled Delivery of Drugs, R.L. Juliano, Ed., Annals of NY Academy Sci., NY Academy Sci., NY, 1987.
26. Physics of Complex and Supramolecular Fluids, S.A. Safran, N.A. Clark, Eds., Exxon Monograph, J. Wiley & Sons, NY 1987.
27. "Controlled Release Technology", Chem. Abst. Selects, SVC 05C, Chemical Abstracts Service, Columbus, OH.
28. Microproducts/3M, Imaging Systems Div., St. Paul, MN, Product Bulletins.

29. Rohm and Haas Co., Philadelphia, PA, Specialty Chemicals - Surfactants Technical Bulletins.

30. Moleculon BioTech, Inc., Cambridge, MA, Controlled Release Technology Information.

31. Detergents.....Special Report, C&E News, 21, Jan. 25, 1988.

32. Micro Membranes, Inc., Technology and Product Bulletins.

SPECIALTY POLYMER REFERENCES

33. Controlled Release of Drugs-Polymers and Aggregate Systems, M. Rosoff (Ed.), VCH Publications, NY (1988).

34. Carbon Fibers, Technology, Uses and Properties: Plastics and Rubber Institute, London (Ed.), Noyes Publications, NJ, 1986.

35. High Temperature Polymer Matrix Composites, T.T. Serafini, Noyes Publ., NY. (1988)

36. "Polymers for Fibers and Elastomers", J.C. Arthur, Jr. Ed., ACS Sympos. Series 260, Amer. Chem. Soc., Washington, DC. 1983.

37. Journals: Macromolecules; Ind. & Eng. Chem. Research; Modern Plastics; Polymer Processing; Preprints Polymeric Materials Science and Engineering; Polymer Preprints; SAMPE (J. Soc. for Advancement of Material and Process Engineering); J. Composites Technology and Research.

10.2 Bibliography on Supercritical Fluid Technology.

1. "Supercritical Fluid Extraction", McHugh and Krukonis, Johns Hopkins Univ., 1987. Course Contents and book, Principles and Practice, Butterworths, Boston, MA, 1986.

2. "Supercritical Fluid Extraction and Chromatography", ACS Sympos. Series 366, B.A. Charpentier, M.R. Sevenants, Eds., Amer. Chem. Soc., Washington, DC, 1988.

3. Chemical Engineering at Supercritical Fluid Conditions, M.E. Paulaitis, J.M.L. Penninger, R.D. Gray, Jr., P. Davidson Eds., Ann Arbor Science, Ann Arbor, MI, 1983.

4. "Supercritical Fluids; Chemical and Engineering Principles and Applications", T.G. Squires, M.E. Paulaitis, Eds., ACS Sympos. Series 329, Amer. Chem. Soc., Washington, DC, 1987.

5. "Supercritical Fluid Chromatography", R.M. Smith, Ed., RSC Chromatography Monographs, Royal Society Chem., London, 1988.

10.3 Bibliography on Sensors/Artificial Intelligence.

1. Solid State Chemical Sensors, J. Janata, R.J. Huber, Eds., Academic Press, NY, 1985.

2. "Fundamentals and Applications of Chemical Sensors", ACS Sympos. Series 309, Amer. Chem. Soc., Washington, DC.

3. "Sensors, The Journal of Machine Perception", Helmers Publ., Peterborough, NH.

4. "AI Magazine", Amer. Assoc. Artif. Intelligence, Stanford Univ., Palo Alto, CA.

5. "Expert Systems", Intern. J. Knowledge Engineering, Learned Information, Inc. Publ., Oxford, NJ.

6. Internat. J. Pattern Recognition and Artificial Intelligence, World Scientific Publ. Co. Pte Ltd. Singapore.

7. IEEE Expert, Institute of Electrical and Electronic Engineers, Inc., Computer Soc. IEEE Publ., NY.

8. Strouf, O., Chemical Pattern Recognition, Research Studies Press, J. Wiley & Sons, England, 1986.

9. James, M., Classification Algorithms, J. Wiley-Intersci. NY, 1985.

10. Jurs, P.C., Computer Software Applications in Chemistry, J. Wiley & Sons, NY, 1986.

11. Computer-Supported Spectroscopic Databases, J. Zupan, Ed., Halsted Press, J. Wiley, NY, 1986.

12. Principles of Database Design, S.B. Yao, Ed., Vol. 1, Logical Organizations, Prentice-Hall, NJ, 1985.

13. Pal, S.K., and Majumder, D.K.D., Fuzzy Mathematical Approach to Pattern Recognition, Halsted Press, J. Wiley & Sons, NY, 1986.

14. G.L. Simons, Expert Systems and Micros., Nat. Computing Centre Ltd., England and Halsted Press, NCC, National Centre for Information Technology, 1985.

15. Sharaf, M.A., Illman, D.L., and Kowalski, B.R., Chemometrics, Wiley-Interscience, Chemical Analysis Monographs, Vol. 82, 1986.

16. Wolff, D.D., and Parsons, M.L., Pattern Recognition Approach to Data Interpretation, Plenum Press, NY, 1983.
17. Handbook of Artificial Intelligence, Barr, A., Feigenbaum, E.A., Eds., Vols. 1-3, W. Kaufmann, Inc. Los Altos, CA, 1982.
18. Booth, G.M., Design of Complex Information Systems: Common Sense Methods for Success, McGraw-Hill Co., NY, 1983.
19. Harmon, P., and King, D., Expert Systems, AI in Business, J. Wiley & Sons, NY, 1985.
20. Negoita, C.V., Expert Systems and Fuzzy Systems, Benj. Cummings Publ. Co., Menlo Park, CA, 1985.
21. Kandel, A., Fuzzy Techniques in Pattern Recognition, Wiley-Intersci., NY, 1982.
22. The AI Business; Commercial Uses of AI, P.H. Winston, K.A. Prendergast, Eds., MIT Press, Cambridge, MA, 1984.
23. Pearl, J., Heuristics; Intelligent Search Strategies for Computer Problem-Solving, Addison-Wesley Publ. Co., Reading, MA, 1984.
24. Abelson, H., and Sussman, G.J., Structure and Interpretation of Computer Programs, MIT Press, Cambridge, MA, 1985.
25. Inose, H., and Pierce, J.R., Information Technology and Civilization, W.H. Freeman and Co., NY, 1984.
26. Hillis, W.D., The Connection Machine, MIT Press, Cambridge, MA, 1985.
27. Schutzer, D. Artificial Intelligence; An Applications-Oriented Approach, Van Nostrand Reinhold Co., NY, 1987.
28. Teorey, T.J., and Fry, J.P., Design of Database Structures, Prentice-Hall, Inc., NJ, 1982.
29. Bramer, M., and Bramer, D., The Fifth Generation, Addison-Wesley Publ. Ltd., England, 1984.
30. Hayes-Roth, F., Waterman, D.A., and Lenat, D.B., Building Expert Systems, Vol. 1, Addison-Wesley Publ. Co., Inc., Reading, MA, 1983.
31. Buchanan, B.G., and Shortliffe, E.H., Rule-Based Expert Systems; MYCIN Experiments of the Stanford Heuristic Programming Project, Addison-Wesley Publ. Co., Reading, MA, 1984.

32. Martin, J., "Fourth-Generation Languages", Vol. II, Representative 4GLs, Prentice-Hall, NJ, 1986.

33. Ioffe, I.I., Application of Pattern Recognition to Catalytic Research, Research Studies Press Ltd., Wiley & Sons, NY, 1988. (GEO-CENTERS, INC.)

34. Gray, N.A.B., Computer-Assisted Structure Elucidation, J. Wiley & Sons, Inc., NY, 1986.

NEURAL NETWORKS REFERENCES.

35. AI Week and Neural Network News, P.O. Box 2513, Birmingham, AL 35201.

36. "NEUROCOMPUTERS", Neuroengineering News, Sixth Generation Reviews, Gallifrey Publishing, P.O. Box 155, Vicksburg, MI 49097.

37. Wasserman, P.D., and Schwartz, T., Parts 1 and 2: "Neural Networks...What Are They and Why is Everybody Interested in Them Now?", IEEE Expert, Spring 1988.

38. Wihl, M.G., "Neural Networks Make Cybernetic Multisensor Systems Possible", Sensors, 5, No. 9, 11, Sept. 1988.

39. Lerner, E.J., "Computers that Learn", Aerospace America, 33, June 1988.

40. Pennisi, E., "Neural Nets are Sparking Heated Debates....", The Scientist, 2, No. 21, Nov. 14, 1988, Neural net research funding by industry....Oct. 17, 1988.

41. Johnson, R.C. and Brown, C., Cognizers - Neural Networks and Machines That Think, J. Wiley & Sons, NY, 1988.

42. DARPA Neural Networks Study, October 1987 - February 1988, AFCEA International Press, Fairfax, VA (1988) 800-336-4583 or 6190.

43. Widman, L.E., Loparo, K.A. and Nielsen, N.R., Artificial Intelligence, Simulation and Modelling, J. Wiley & Sons, NY (1989).

44. PAO, Yoh-Han, Adaptive Pattern Recognition and Neural Networks, Addison-Wesley Publishing Co., Reading, Ma (1989).

45. Johns Hopkins University, Baltimore, MD course on Neural Networks (March 1989).

46. Souček, B. and Souček, M., "Neural and Massively Parallel Computers," J. Wiley and Sons, NY, 1989.

10.4 Bibliography on Inverse GC and SFC.

1. Aspler, J.S., "Theory and Applications of Inverse GC", Chapter 9, Pyrolysis and GC in Polymer Analysis, Liebman and Levy, Eds., M. Dekker, Inc., 1985.
2. Vukov, A.J., and Gray, D.G., "Inverse GC Study of Carbon Fiber Surfaces", Preprints, Polym. Materials: Sci. Engin. Div., Third Chem. Congress of N. Amer., Toronto, Canada, 58, pp 917-21, 1988.
3. Raymer, J.H., Cooper, S.D., and Pellizzari, E.D., "Studies into Analyte Competition on Polyimide Sorbents Using a Deuterated Tracer Pulse Technique", Ibid, pp 923-927.
4. Guillet, J.E., "Studies of Polymer Structure and Interactions by Automated Inverse GC", Ibid, 644-49.
5. Wesson, S.P., "Surface Energetics of Plasma Treated Carbon Fiber Reinforcements", Ibid, pp 650-654.
6. Inverse GC Topics, by various research groups, Ibid, pp 655-751.
7. Liebman, S.A., Ahlstrom, D.H., Foltz, C.R., "Inverse GC Studies of Poly(vinyl chloride)", J. Chromatogr., 67, 153 (1972).
8. Course on High Resolution GC and SFC, Center for Prof. Advancement, Prof. Peter C. Uden, P.O. Box H, E. Brunswick, NY 08816.

10.5 Bibliography of Geo-Centers/CRDEC/BRL Presentations/ Publications Resulting from Task and Allied Program Related to Polymer/Microencapsulation/Trace Organic Analysis.

10.5.1 Presentations.

"Artificial Intelligence Applications to Materials Characterization", S.A. Liebman, to Research Directorate, CRDEC, APG, MD, 9 May 86.

"An Expert System Network for Chemical Characterization of Materials", S.A. Liebman (GEO-CENTERS, INC.), E.W. Sarver, D.J. Reutter (CRDEC), R.A. Fifer (BRL), A.M. Harper (NWC), Amer. Assoc. Artificial Intelligence Workshop, Philadelphia, PA, May 1986, and presentation at Workshop on Coupling Symbolic and Numeric Computing in Expert Systems, Boeing Computer Services Co., sponsored by AAAI, August 1985.

"Analytical Pyrolysis Technology in an Expert System Network for Materials Characterization", S.A. Liebman, A.P. Snyder, R.A. Fifer, A.M. Harper, Federation of Anal. Chem. & Spec. Soc. (FACSS), St. Louis, MO, 28 Sep - 3 Oct 86.

"Time-Resolved Analytical Pyrolysis Analytical Pyrolysis Studies of Nitramine Decomposition with a Triple Quadrupole MS/MS System", Federation of Anal. Chem. & Spec. Soc. (FACSS), St. Louis, MO, 28 Sep - 3 Oct 86.

"Artificial Intelligence Applications to Materials Characterization", S.A. Liebman, Eastern Analysis Symposium, New York City, NY, 19-20 Oct 86.

"How Far Can You Go with Small Expert Systems?", A.M. Harper, S.A. Liebman, Scientific Computing & Automation Conference & Exposition, Atlantic City, New Jersey, 5-7 Nov 86.

"The EXMAT Project: An Example of a Networked Expert System for Chemical Research", A.M. Harper, S.A. Liebman, R.A. Fifer, M. Shen, Scientific Computing & Automation Conference and Exposition, Atlantic City, NJ, 5-7 Nov 86.

"Development of an Expert System Network for Materials Analysis", S.A. Liebman, J. Leonard, D.J. Reutter, E.W. Sarver, R.A. Fifer, A.M. Harper, CRDEC Conference on Defense Research, 18 November 1986.

"Analyte Detection in Polymer Matrices by Pyrolysis-Atmospheric Pressure Chemical Ionization MS/MS", A.P. Snyder, J.H. Kremer, R.A. Yost, S.A. Liebman, Workshop on Macromolecules in Microencapsulation Research and Technology, CRDEC, APG, MD, 3-4 December 1986.

"Hyphenated Instruments - A Problem Solver's Route", S. A. Liebman, Pittsburgh Conference and Exposition, Atlantic City, NJ, March 9-13, 1987.

"Analytical Instrumentation for Supercritical Fluid Extraction, Desorption, and Chromatography", S.A. Liebman, R. Smardzewski, D. Reutter, L. Schiff, and R.A. Fifer, ARO/ONR Workshop on Applications of Supercritical Fluid Technologies: Research Needs for DoD Applications, Seattle, Washington, May 6-8, 1987.

"Integrated Instruments and Expert System Networks", S.A. Liebman, Spring Symposium-Thermal Analysis Forum Delaware Valley, Wilmington, DE, May 21, 1987

"Analyte Detection in Solid Polymer Matrices by Pyrolysis-Atmospheric Pressure Chemical Ionization MS/MS", Snyder, Kremer (CRDEC, Yost (University of Florida), Liebman (GEO-CENTERS, INC.) 35th ASMS Conference on Mass Spectrometry and Allied Topics, Denver, CO, May 24-28, 1987, pg. 181.

"Development of an Expert System Network for Materials Analysis", S.A. Liebman (GEO-CENTERS, INC.), J.W. Leonard, D.J. Reutter, E.W. Sarver (CRDEC), R.A. Fifer (BRL), A.M. Harper (Naval Weapons Center, China Lake) Annual Mtg, Amer. Assoc. Artif. Intell., Philadelphia, PA, May 1987.

"Instruments for Advanced Sample Processing - Pyrolysis, DHS, Extraction/Reaction with On-Line Analysis", S.A. Liebman, presented at Chromatography Forum Meeting, Delaware Valley, 9-10 June 1987.

"Polymer Characterization by Pyrolysis-Atmospheric Pressure Chemical Ionization Mass Spectrometry", A.P. Snyder, J.H. Kremer (CRDEC), R.A. Yost (University of Florida), S.A. Liebman (GEO-CENTERS, INC.) Presentation at Gordon Research Conference, Plymouth, NH, July 1987.

Invited "Integrated Intelligent Instruments (I³) and Applied AI", S.A. Liebman, Dow Chemical Co., Analytical Division, Midland, MI, 16 Sep 87.

"Pyrolysis/Concentrator/GC-MS Studies of Amine Salts in an Artificial Intelligence Network", M.B. Wasserman, D.J. Reutter, A.P. Snyder, J.W. Leonard, S.A. Liebman, and W. Windig, 1987 CRDEC Conference on Defense Research, November 1987.

"Development of an Artificial Intelligence Prototype for Pyrolysis/Concentrator/GC-MS Analyses", A.M. Harper (NWC), D.J. Reutter, J.W. Leonard, M.B. Wasserman (CRDEC), S.A. Liebman (GEO-CENTERS, INC.), 1987 CRDEC Conference on Chemical Defense Research, November 1987.

"Advanced Pairs Expert System for Mixture Analysis within EXMAT Network", S.A. Tomellini, R.R. Smardzewski, D.J. Reutter, J.W. Leonard, R.A. Fifer, S.A. Liebman, 1987 CRDEC Conference on Chemical Defense Research, November 1987.

"Integrated Intelligent Instrumentation (I³) for Trace Organic Analysis", S.A. Liebman, Joint Chromatography Forum/Delaware Valley, 17 November 1987.

"Integrated Intelligent Instruments (I³): Supercritical Fluid Extraction, Desorption, Reaction and Chromatography". E.J. Levy, S. Lurcott, S. O'Neill, J. Gurthrie, S. Yocklovich, S.A. Liebman, The Pittsburgh Conference and Exposition, New Orleans, LA, 22-26 Feb 88.

"Trace Analyses with Supercritical Fluid Chromatography: Advanced Column Technologies for Packed Microbore and Capillary Systems", S. Yocklovich, S. O'Neill, J. Guthrie, E.J. Levy, and S.A. Liebman, The Pittsburgh Conference and Exposition, New Orleans, LA, 22-26 Feb 88.

"Analytical Pyrolysis/Concentrator GC-MS Studies using Applied Artificial Intelligence", M.B. Wasserman, D.J. Reutter, A.M. Harper, S.A. Liebman, The Pittsburgh Conference and Exposition, New Orleans, LA 22-26 Feb 88.

"A Fast Self-Modelling Curve-Resolution Method for Time-Resolved Mass Spectral Data", W. Windig and S.A. Liebman, A.P. Snyder, and M.B. Wasserman, The Pittsburgh Conference and Exposition, New Orleans, LA 22-26 Feb 88.

"Propellant Characterization with SFE-SFC in an Expert Systems Network", S.A. Liebman, E.W. Sarver, A.P. Snyder, R.A. Fifer, F. Shaw, A.M. Harper, Presented at The Joint International Symposium on Compatibility of Plastics and Other Materials with Explosives, Propellants, Pyrotechnics and Processing of Explosives, Propellants and Ingredients, (American Defense Preparedness Association) New Orleans, LA, 18-20 April 1988. BEST PAPER AWARD.

"Structures of RDX Pyrolysis Products Using Atmospheric Pressure Chemical Ionization Tandem Mass Spectrometry (APCI-MS/MS) With Isotopic Tagging", R.A. Fifer, S.A. Liebman, A.P. Snyder, J.H. Kremer, M.A. Schroeder, JANNAF Kinetics Panel Meeting, APL, 2-4 May 1988.

(Invited) "Integrated Instruments & Applied Artificial Intelligence", S.A. Liebman, Ohio Valley Thermal Analysis Symposium on Pyrolysis of Polymers and Dielectric Thermal Analysis", Dayton, OH, 11 May 1988.

"Integrated Intelligent Instruments in Materials Analysis", S.A. Liebman, NASA Lewis, Cleveland, OH, May 1988.

"Integrated Instruments and Expert System Networks", S.A. Liebman, Wright-Patterson AFB, OH, 12 May 1988.

"Integrated Intelligent Instruments (I^3) for Energetic Materials Research", S.A. Liebman, Sandia National Laboratory, Albuquerque, New Mexico, September 20, 1988.

"Integrated Intelligent Instruments (I^3) for Trace Organic Analysis", S.A. Liebman, Los Alamos National Laboratory, Los Alamos, New Mexico, September 21, 1988.

"Integrated Intelligent Instruments (I³) for Trace Organic Analysis", S.A. Liebman, NASA Ames, Moffett Field, CA, September 23, 1988.

"Integrated Intelligent Instruments (I³) for Energetic Materials Research", S.A. Liebman, Sandia-Livermore, CA, September 26, 1988.

"Analytical Chemistry of Microencapsulated Materials - Integrated Intelligent Instruments (I³) for CBW Defense", presented at CRDEC, Sep 88.

"An Expert System Network for Chemical Systems", A.M. Harper, S.A. Liebman and D.J. Reutter, Computers in Chemistry Division, American Chemical Society National Meeting, Los Angeles, CA, September 1988.

"Analytical Instrumentation and Applied Artificial Intelligence in Materials Science", S.A. Liebman, R.R. Smardzewski, E.W. Sarver, D.J. Reutter, A.P. Snyder, A.M. Harper, E.J. Levy, S. Lurcott, S. O'Neill, Polymeric Materials Science and Engineering Division, American Chemical Society National Meeting, Los Angeles, CA, September 1988.

"Supercritical Fluid Applications in Materials Science", S.A. Liebman, ASTM Committee E-19 on Chromatography, Baltimore, MD, 11 Oct 88.

"Artificial Intelligence" to French Representatives, CRDEC, APG, MD, 25 October 1988.

"Analytical Toxicology: An Overview of Trace Organic Analytical Instrumentation Methods", S.A. Liebman, Association of Government Toxicologists, "New Techniques and Concepts for Reducing Drug Toxicities", Rockville, MD, 26-27 October 1988.

"Integrated Intelligent Instruments (I³) for Trace Organic Analysis", Northeastern University, Boston MA, 28 Oct 1988.

"Integrated Intelligent Instruments (I³) and Applied AI in Analytical Chemistry", S.A. Liebman, W. Windig, S.L. Emery, D.J. Reutter, A.P. Snyder, E.J. Levy, S. Lurcott, Federation of Analytical Chemistry and Spectroscopy Studies Fifteenth Annual Meeting, Boston, MA, 30 Oct - 4 Nov 88.

"Explosives, Drug and Technical Applications of Pyrolysis/Atmospheric Pressure Ionization/Tandem Mass Spectrometry", S.A. Liebman, A.P. Snyder, R.A. Yost, Federation of Analytical Chemistry and Spectroscopy Studies Fifteenth Annual Meeting, Boston, MA, 30 Oct - 4 Nov 88.

"Analysis of Volatiles on Solid Sorbents by Thermal Desorption/Cryotrapping and GC/MS Detection", S.A. Liebman, M.B. Wasserman and M.E. Brooks, 1988 Chemical Conference on Chemical Defense Research, CRDEC, APG, MD, November 88 (Poster Session)

"Characterization of Microencapsulated Materials with Analytical Pyrolysis and Supercritical Fluid Technologies", S.A. Liebman, A.P. Snyder, M.B. Wasserman, S.L. Emery, M.E. Brooks, 1988 Conference on Chemical Defense Research, CRDEC, APG, MD, November 88.

"Chemical Agent Simulants for Testing Transparent Materials" by R.E. Lewis, S.A. Liebman, L. Isaacson, P. Grasso, and E.W. Sarver, 1988 Conference on Chemical Defense Research, CRDEC, APG, MD, November 1988.

"Ethical Elements of Chemical/Biological Defense Research", S.A. Liebman, presented at the Conference on "Ethical Issues Associated with Scientific Research for the Military", Arlington, VA, 26-28 January 1989.

"An Expert System Network for Supercritical Fluid Technologies", S.A. Liebman, T. Ryan, S. Yocklovich, S. Lurcott, J. Watkins and E.J. Levy, The Pittsburgh Conference and Exposition, March 1989, Paper No. 109.

"Supercritical Fluid Extraction-Chromatography: Profiles and Analysis for QC/QA and Process Development", S.A. Liebman, A.P. Snyder, S. Yocklovich, S. Lurcott, J. Watkins and E.J. Levy, The Pittsburgh Conference and Exposition, March 1989, Paper No. 920.

"Detection Systems for Supercritical Fluid/GC Instrumentation: Flame Ionization Detector (FID) and Fiber Optic Monitor (FOM) Units", S.A. Liebman, R.A. Fifer, P. Griffiths, S. Lurcott, B. Bergman, and E.J. Levy, The Pittsburgh Conference and Exposition, March 1989, Paper No. 1545.

"Detection and Identification of Microencapsulated Biological and Chemical Materials by Advanced Analytical Instrumentation", A.P. Snyder, S.A. Liebman, W. Windig, Chemical Symposium on "Agents of Biological Origin", Laurel, MD, March 1989.

"Integrated Intelligent Instruments (I^3) in Materials and Environmental Sciences", S.A. Liebman, A.P. Snyder, M.B. Wasserman, M.E. Brooks, J. Watkins, S. Lurcott, S. O'Neill, and E.J. Levy, "International Conference on Analytical Chemistry", 30 Jul - 5 Aug 89, University of Cambridge, UK.

"Microencapsulated Materials and Advances Analytical Instrumentation", Chemical Protection Symposium, Umea, Sweden, 11-16 June 1989.

10.5.2 Publications.

Liebman, S.A., Leonard, J., Reutter, D.J., Sarver, E.W., Fifer, R.A., Harper, A.M., "Development of an Expert System Network for Materials Analysis", Proceedings of the 1986 CRDEC Conference on Defense Research, November 18-21, 1986, CRDEC-SP-87008, compiled by M.D. Rausa.

Liebman, S.A., Duff, P.J., Pickie, K.D., Schroeder, M.A., Fifer, R.A., "Degradation Profile of Propellant Systems with Analytical/Pyrolysis/Concentrator/GC Technology", J. Hazardous Materials, 13, p 51-56 (1986).

Liebman, S.A., Duff, P.J., Schroeder, M.A., Fifer, R.A., Harper, A.M., "Concerted Organic Analysis of Materials and Expert Systems Development", ACS Sympos. Series No. 306, Artificial Intelligence Applications in Chemistry, Pierce and Hohne, Eds., Amer. Chem. Soc., Washington, DC, pp 365-384, 1986.

Liebman, S.A., Snyder, A.P., Kremer, J.H., Reutter, D.J., Schroeder, M.A., Fifer, R.A., "Time-Resolved Analytical Pyrolysis Studies of Nitramine Decomposition with the Triple Quadrupole Mass Spectrometer System", J. Analytical and Applied Pyrolysis, 12, pp 83-95 (1987).

Levy, E.J., Lurcott, S., O'Neill, S., Yocklovich, S., Choen, R., Pfeiffer, K., Wampler, T.P., and Liebman, S.A., "Advances in Supercritical Fluid Systems", Amer. Lab, pp 66-72, Aug 87.

Reutter, D.J., Shuely, W.J., and Liebman, S.A., "Workshop on Macromolecules in Microencapsulation Research and Technology, 3-4 December 1986", CRDEC-SP-87022, D.J. Reutter, W.J. Shuely, S.A. Liebman, August 1987.

Snyder, A.P., Kremer, J.H., Yost, R.A., and Liebman, S.A., "Analyte Detection in Polymer Matrices by Pyrolysis-Atmospheric Pressure Chemical Ionization MS/MS", Workshop on Macromolecules in Microencapsulation Research and Technology, 3-4 December 1986", CRDEC-SP-87022, August 1987.

Wasserman, M.B., Reutter, D.J., Snyder, A.P., Leonard, J.W., Liebman, S.A., and Windig, W., "Pyrolysis/Concentrator/GC-MS Studies of Amine Salts in an Artificial Intelligence Network", Proceedings of the 1987 CRDEC Conference on Defense Research, 17-20 November 1987, CRDEC-SP-88013, compiled by M.D. Rausa.

Harper, A.M., Reutter, D.J., Leonard, J.W., Wasserman, M.B., and Liebman, S.A., "Development of an Artificial Intelligence Prototype for Pyrolysis/Concentrator/GC-MS Analyses", Proceedings of the 1987 CRDEC Conference on Chemical Defense Research, 17-20 November 1987, CRDEC-SP-88013, compiled by M.D. Rausa.

Tomellini, S.A., Smardzewski, R.R., Reutter, D.J., Leonard, J.W., Fifer, R.A., and Liebman, S.A., "Advanced Pairs Expert System for Mixture Analysis within EXMAT Network", Proceedings of the 1987 CRDEC Conference on Chemical Defense Research, 17-20 November 1987, CRDEC-SP-88013, compiled by M.D. Rausa.

Snyder, A.P., Kremer, J.H., Liebman, S.A., Schroeder, M.A., and Fifer, R.A., "Characterization of Cyclotrimethylenetrinitramine (RDX) by N,H Isotope Analyses with Pyrolysis-Atmospheric Pressure Ionization Tandem Mass Spectrometry", Org. Mass Spec., 24, 15-21 (1989).

Liebman, S.A., Levy, E.J., Lurcott, S., O'Neill, S., Guthrie, J., Yocklovich, S., "Integrated Intelligent Instruments: Supercritical Fluid Extraction, Desorption, Reaction and Chromatography" J. Chromat. Sci., 27, 118-126 (1989).

Lewis, R.E., Liebman, S.A., Isaacson, L., Grasso, P.S., and Sarver, E.W., Chemical Agent Simulants for Testing Transparent Materials, CRDEC-CR-88069, May 1988.

Liebman, S.A., Smardzewski, R.R., Sarver, E.W., Reutter, D.J., Snyder, A.P., Harper, A.M., Levy, E.J., Lurcott, S., O'Neill, S., "Analytical Instrumentation and Applied Artificial Intelligence in Materials Science", Preprints, Polymeric Materials Science and Engineering Division, American Chemical Society National Meeting, Los Angeles, CA, September 1988.

Liebman, S.A., "An Overview; Analytical Toxicology and Integrated Intelligent Instruments (I³)", Association of Government Toxicologists, "New Techniques and Concepts for Reducing Drug Toxicities", Rockville, MD, 26-27 October 1988, submitted to Telford Press, 1988.

Liebman, S.A., Wasserman, M.B., and Brooks, M.E., "Analysis of Volatiles on Solid Sorbents by Thermal Desorption/Cyrotapping and GC/MS Detection", Proceedings of the 1988 Conference on Chemical Defense Research, 15-18 November 1988, in press, compiled by J. Williams.

Liebman, S.A., Snyder, A.P., Wasserman, M.B., Emery, S.L., Brooks, M.E., "Characterization of Microencapsulated Materials with Analytical Pyrolysis and Supercritical Fluid Technologies", Proceedings of the 1988 Conference on Chemical Defense Research, 15-18 November 1988, in press, compiled by J. Williams.

Lewis, R.E., Liebman, S.A., Isaacson, L., Grasso, P.S., and Sarver, E.W., "Chemical Agent Simulants for Testing Transparent Materials", Proceedings of the 1988 Conference on Chemical Defense Research, 15-18 November 1988, in press, compiled by J. Williams.

Windig, W., Liebman, S.A., Wasserman, M.B., and Snyder, A.P., "Fast Self-Modeling Curve-Resolution Method for Time-Resolved Mass Spectral Data", Anal Chem, 60, pp 1503-1510 (1988).

Liebman, S.A., Smardzewski, R.R., Snyder, A.P., Wasserman, M.B., Fifer, R.A., Shaw, F., Levy, E.J., Lurcott, S., Yocklovich, S., Ryan, T., "Materials Characterization with Integrated Analytical Instruments and Applied Artificial Intelligence", submitted to American Chemical Society Advances in Chemistry Series, December 1988.

Snyder, A.P., Cross, C.W., Liebman, S.A., Eiceman, G.A., Yost, R.A., "Thermal Desorption-Tandem Mass Spectrometry for the Detection and Identification of Formulated Products", J. Anal. Appl. Pyrolysis, 1989.

Liebman, S.A., "Ethical Elements of Chemical/Biological Defense Research", S.A. Liebman, submitted to The New York Academy of Sciences, New York, January 1989.

Liebman, S.A., Ryan, T., Yocklovich, S., Lurcott, S., Watkins, J., and Levy, E.J., "An Expert System Network for Supercritical Fluid Technologies", The Pittsburgh Conference and Exposition, Atlanta, GA, March 1989, Paper No. 109, in Proceedings.

Liebman, S.A., Snyder, A.P., Yocklovich, S., Lurcott, S., Watkins, J., and Levy, E.J., "Supercritical Fluid Extraction-Chromatography: Profiles and Analysis for QC/QA and Process Development", The Pittsburgh Conference and Exposition, Atlanta, GA, March 1989, Paper No. 920, in Proceedings.

Liebman, S.A., Fifer, R.A., Griffiths, P., Lurcott, S., Bergman, B., and Levy, E.J., "Detection Systems for Supercritical Fluid/GC Instrumentation: Flame Ionization Detector (FID) and Fiber Optic Monitor (FOM) Units", The Pittsburgh Conference and Exposition, Atlanta, GA, March 1989, Paper No. 1545, in Proceedings.

Snyder, A.P., Liebman, S.A., Windig, W., "Detection and Identification of Microencapsulated Biological and Chemical Materials by Advanced Analytical Instrumentation". Chemical Symposium on "Agents of Biological Origin", Laurel, MD, March 1989, in Proceedings.

Figure 1

INTEGRATED INTELLIGENT INSTRUMENTS (I³)

AUTOMATED/INTERFACED

- SAMPLE PROCESSING
- MONITORING/SAMPLING
- ACCUMULATOR/PREP SCALE
- SEPARATION
- DETECTION/IDENTIFICATION
- DATA ANALYSIS/CHEMOMETRICS
- APPLIED AI

Figure 2

TRACE ORGANIC ANALYSIS OF ANALYTES IN GAS, LIQUID, SOLID MATRICES

AGENTS, EXPLOSIVES, PROPELLANTS, ADDITIVES, SIMULANTS, TOXINS, PRIORITY POLLUTANTS, BIOPOLYMERS, SYNTHETIC POLYMERS, FORMULATED MATERIALS, DEGRADATION AND REACTION PRODUCTS, FORENSICS

INSTRUMENTATION AND METHODOLOGY **AUTOMATED SAMPLE PROCESSING/CONCENTRATOR TECHNOLOGY**

THERMAL TREATMENTS

DYNAMIC HEADSPACE (DHS)
 ANALYTICAL PYROLYSIS (PY)
 SIMULATED ENVIRONMENT
 DEGRADATION PROFILES

SUPERCRITICAL FLUID TREATMENT

SUB-, NEAR, SUPERCRITICAL FLUID
 EXTRACTION (SFE)
 DESORPTION (SFD)
 REACTION (SFR)

ON-LINE DETECTION AND IDENTIFICATION

CAPILLARY GC/FID
 GC-FTIR
 GC-MS/MS
 APCI MS/MS
 SF CHROMATOGRAPHY (SFC)

INTEGRATED SYSTEMS

DHS/PGC/FID
 DHS/PGC/MS/MS
 DHS/PGC/FTIR
 PY-MS/MS/MS
 PY-APCI MS/MS/MS
 SFE/SFD/SFR-SFC
 SFE/SFD/SFR-CAPILLARY GC/FID

ADVANCED DATA ANALYSIS

EXPERIMENTAL DESIGN
 CHEMOMETRICS
 APPLIED ARTIFICIAL INTELLIGENCE

Figure 3

MICROENCAPSULATED SYSTEMS AND THE 13 APPROACH

CHEMICAL SYSTEMS:

COMPONENTS

- o SHELL - NATURAL AND SYNTHETIC POLYMERS
- o CORE - PESTICIDES, TOXINS, SIMULANTS, BIOLOGICALS, VOLATILE/REACTIVE, ENERGETIC CHEMICALS
- o IMMOBILIZED SPECIES - ORGANICS AND INORGANICS IN/ON MATRIX
- o ADDITIVES - BINDERS, SURFACTANTS, EMULSIFIERS, THICKENERS, PROCESSING AIDS, STABILIZERS, SOLVENTS, INORGANICS, CATALYSTS, BUFFERS, SIZINGS/INTER-FACE - SURFACE TREATMENT CHEMICALS

CONCENTRATION RANGE

TARGET SPECIES: TRACE (MICROGRAM/MANOGRAM) TO 50-80% BY WT

SAMPLE (TARGET AND MATRIX) AMOUNT AVAILABLE

TRACE (MICROGRAM/MANOGRAM) TO UNLIMITED

SAMPLING PROTOCOLS - ENVIRONMENTAL


AIR/AEROSOLS, WATER, SOLIDS

INSTRUMENTATION AND METHODS

ENCODE ANALYTICAL STRATEGY AND SPECIFIC CONFIGURATIONS/METHODS INTO CRDEC AI EXPAT NETWORK

CRDEC Instrumentation Schematic

Figure 4

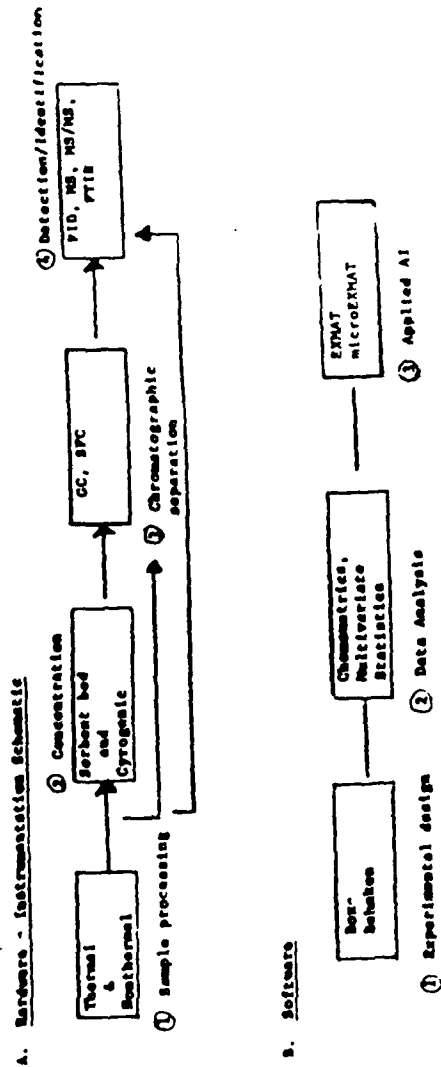

TYPICAL ANALYSIS OF VARIANCE OUTPUT FROM BOX-BEHNKEN


Figure 5

C:\PASCAL>TYPE OUTPUT.DAT

ANALYSIS OF VARIANCE

SOURCE OF VARIATION	SUM OF SQUARES	DEGREES OF FREEDOM	MEAN SQUARE
TOTAL SUM OF SQUARES	2561.0000	96	
DUE TO THE MEAN	2439.6739	1	
DUE TO FIRST ORDER TERMS	7.2500	5	
DUE TO SECOND ORDER TERMS			
INTERACTIONS	16.5000	10	
QUADRATICS	21.6592	5	
LACK OF FIT	67.0836	19	
ERROR SUM OF SQUARES	8.8333	6	
			3.5307
			1.4722

POOLED ESTIMATE OF VARIANCE = 3.0367

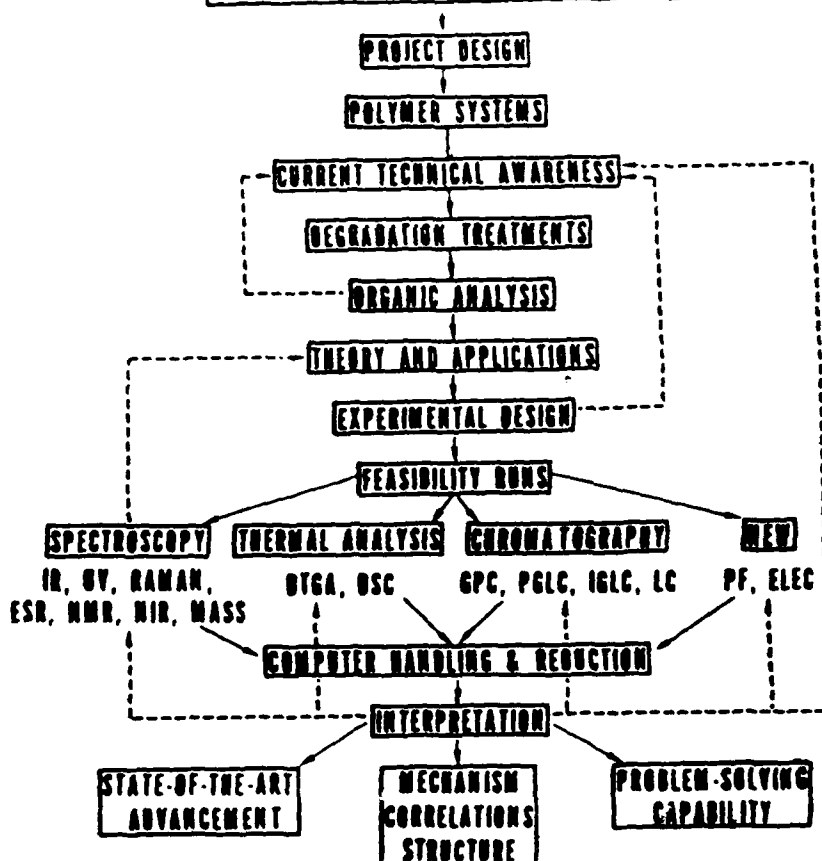

PROBLEM-SOLVING IN MATERIAL SCIENCE (1986)

Figure 6(b)

POLYMER CHARACTERIZATION AND DEGRADATION MECHANISMS

GENERAL PERSPECTIVE

- MICROPROCESSOR - BASED INSTRUMENTATION AND METHOD DEVELOPMENT
- SAMPLE PROCESSING WITH TRACE DETECTION/IDENTIFICATION
- CHEMOMETRICS AND APPLIED AI

SPECIFIC EXPERIMENTATION

- ANALYTICAL PYROLYSIS/CONCENTRATOR GC, GC-MS, MS/MS AUTOMATED SYSTEMS
- INVERSE GC FOR SURFACE, SUBSURFACE ANALYSIS THERMODYNAMIC CONSTANTS, SOLUBILITY PARAMETERS, DIFFUSION CONSTANTS
- SUPERCRITICAL FLUID TECHNOLOGY EXTRACTION, DESORPTION, REACTION, CHROMATOGRAPHY

CURRENT FOCUS

INTEGRATED INTELLIGENT INSTRUMENTS (I³)


EMBEDDED INTELLIGENCE AND DECISION-SUPPORT/INTERPRETATION

INSTRUMENTATION FOR SAMPLE PROCESSING AND ON-LINE ANALYSIS

1969 PYROPROBE ANALYTICAL PYROLYSIS CHEMICAL DATA SYSTEMS, INC.

EARLY AND MID '70s

TRACE ORGANIC ANALYSIS/ENVIRONMENTAL

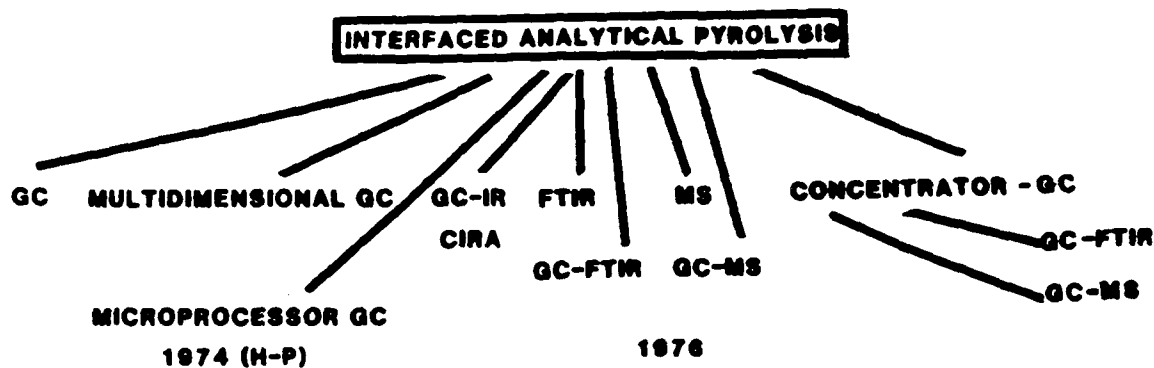

ADVANCED SAMPLE PROCESSING AND ANALYSIS

Figure 7(b)

FULLY AUTOMATED SYSTEMS

LATE '70s
EARLY '80s

**SAMPLE PROCESSING
AND CONCENTRATION**

● THERMAL TREATMENTS
PULSE, PROGRAMMED PYROLYSIS
CONTROLLED ATMOSPHERES

● CONCENTRATOR TECHNOLOGY
SORBENT AND CRYOGENIC ENRICHMENT
DYNAMIC HEADSPACE (DHS)

SEPARATION

● PACKED, FUSED SILICA CAPILLARY
MULTIDIMENSIONAL GC COLUMNS


**DETECTION AND
DATA ANALYSIS**

● GC DETECTORS
TCD, FID, PID, NPD, FPD, MED


● SPECTRAL DETECTORS
FTIR, MS, MS/MS

● CHEMOMETRICS AND AI

**PYROLYSIS /CONCENTRATOR GC/FID
PYROPROBE 123**

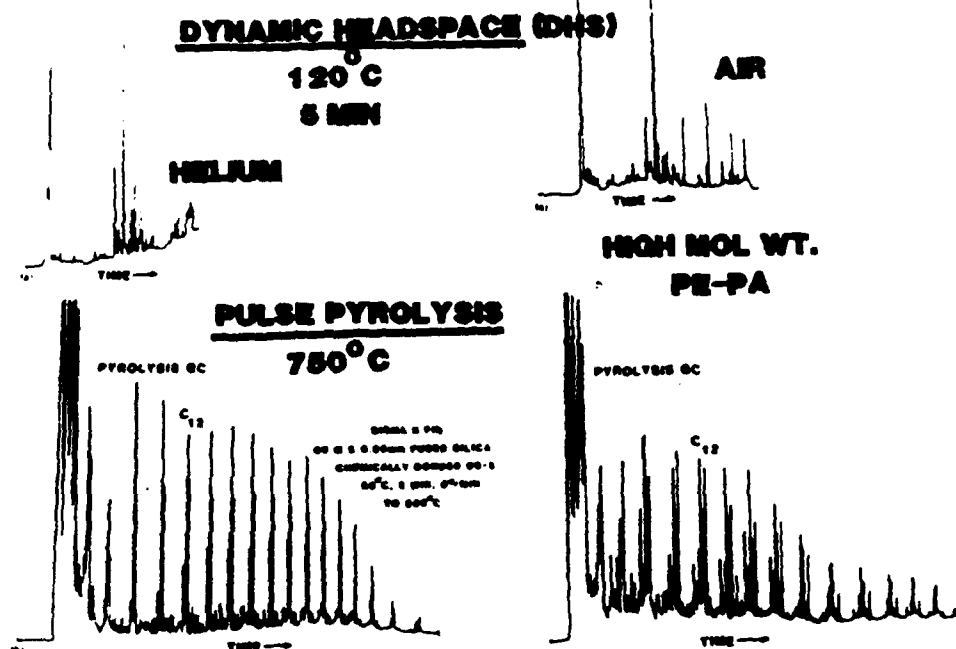


Pulse and Programmed Pyrolysis of Polyethylene

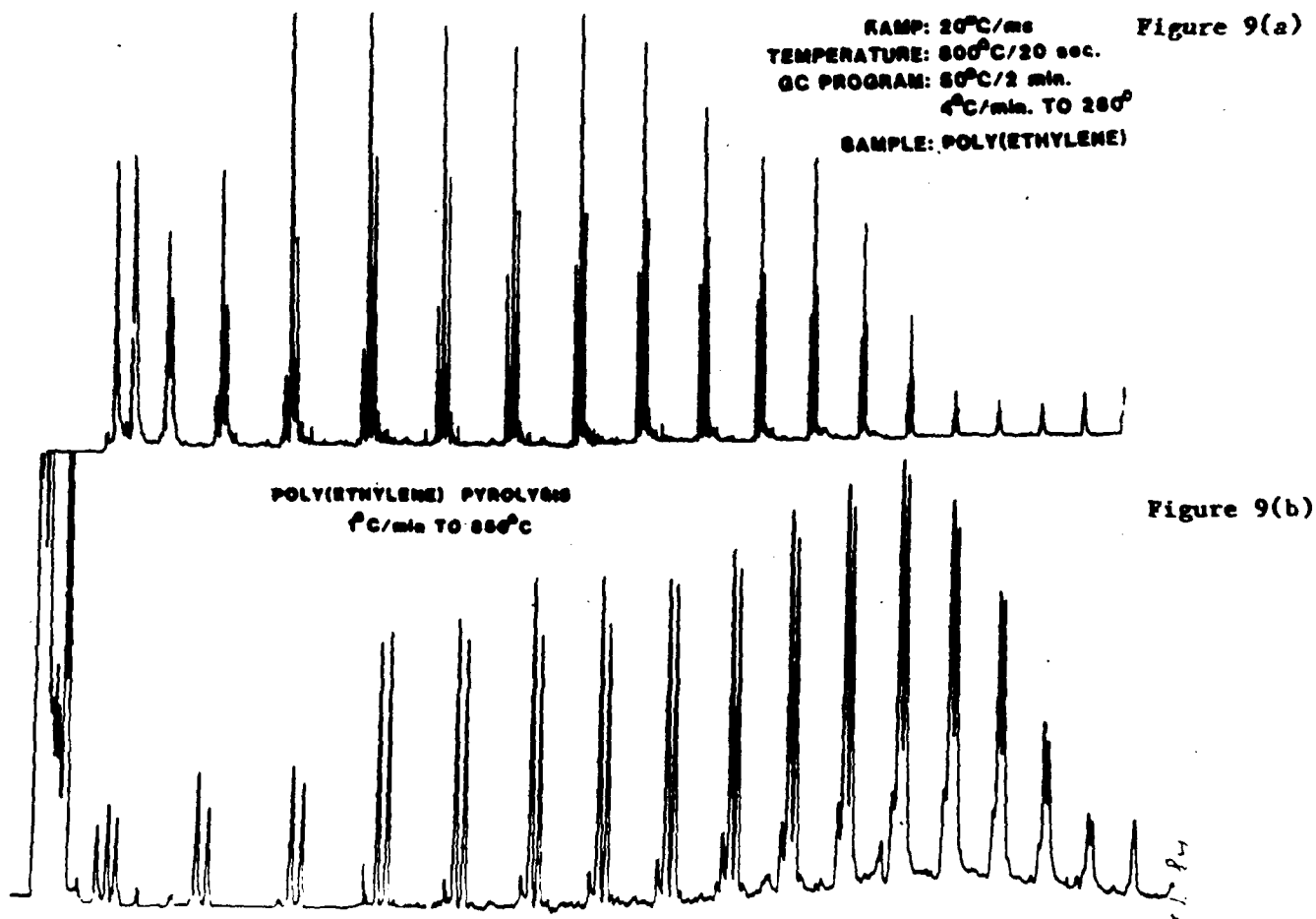


**CDS 320 SAMPLE CONCENTRATOR
VOLATILES IN PHARMACEUTICALS**

Figure 10(a)

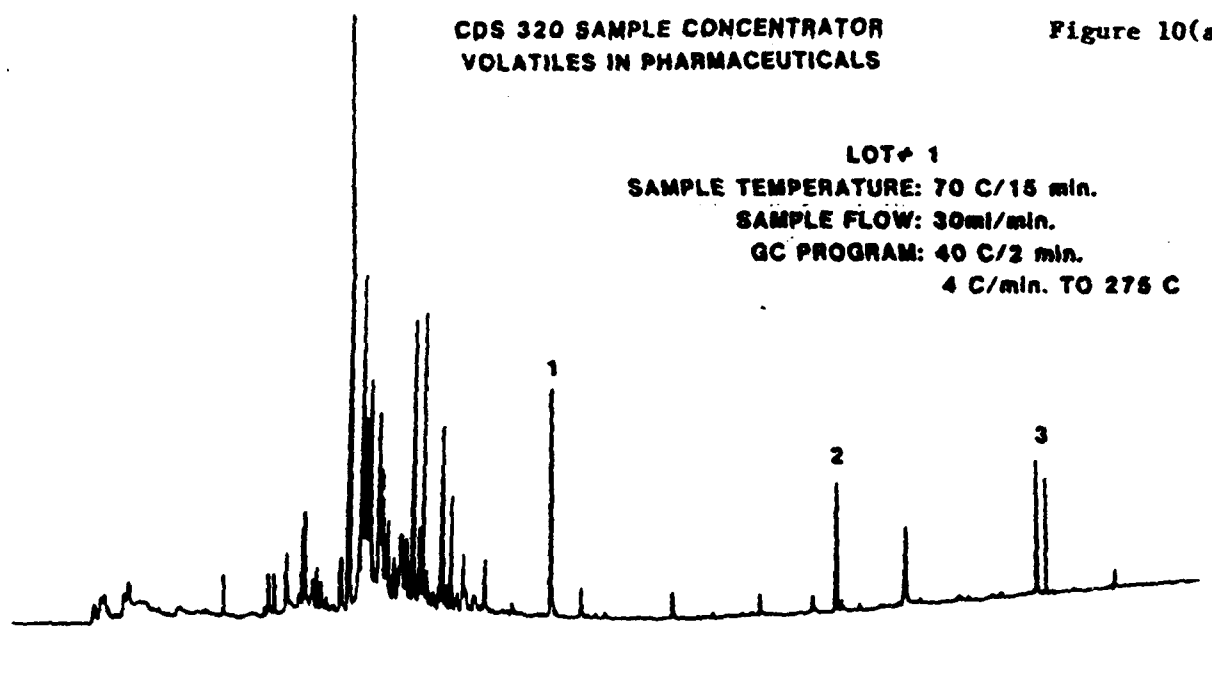


**CDS 320 SAMPLE CONCENTRATOR
VOLATILES IN PHARMACEUTICALS**

Figure 10(b)

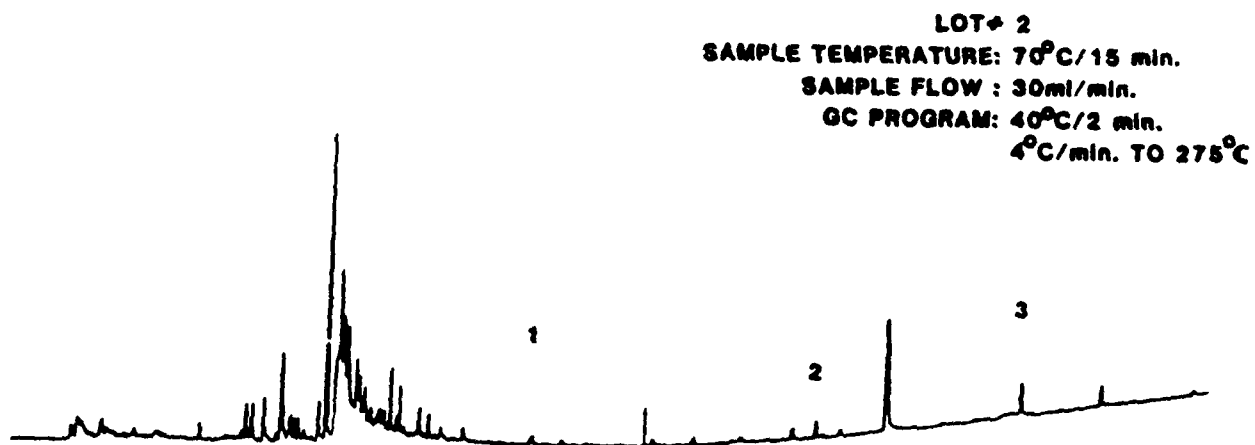

Pyrolysis and Reproducibility of Glass-Lined System CDS Pyroprobe Model 124
(Human Hair)

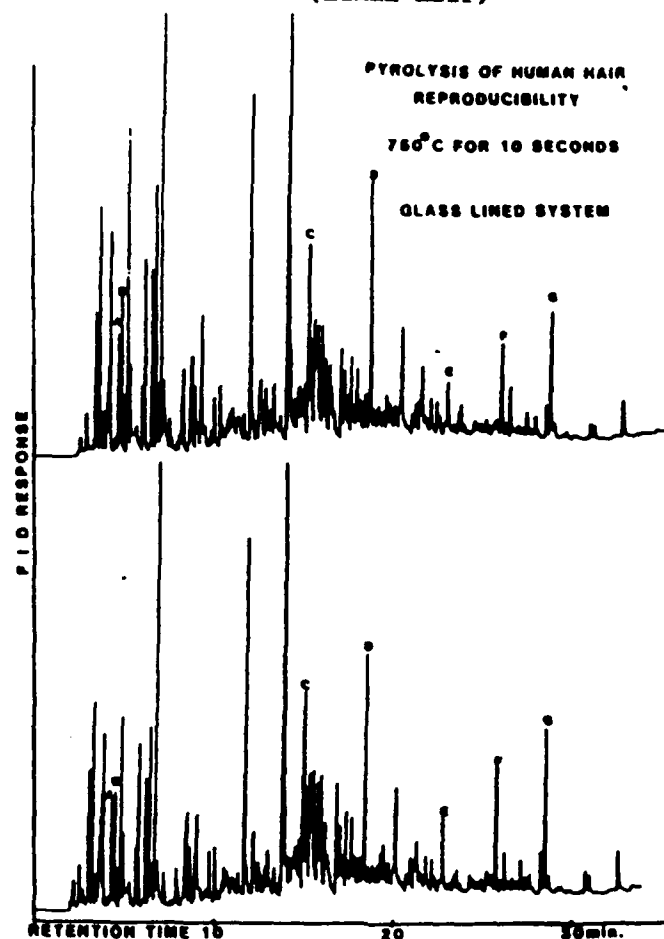

(Raw Wool)



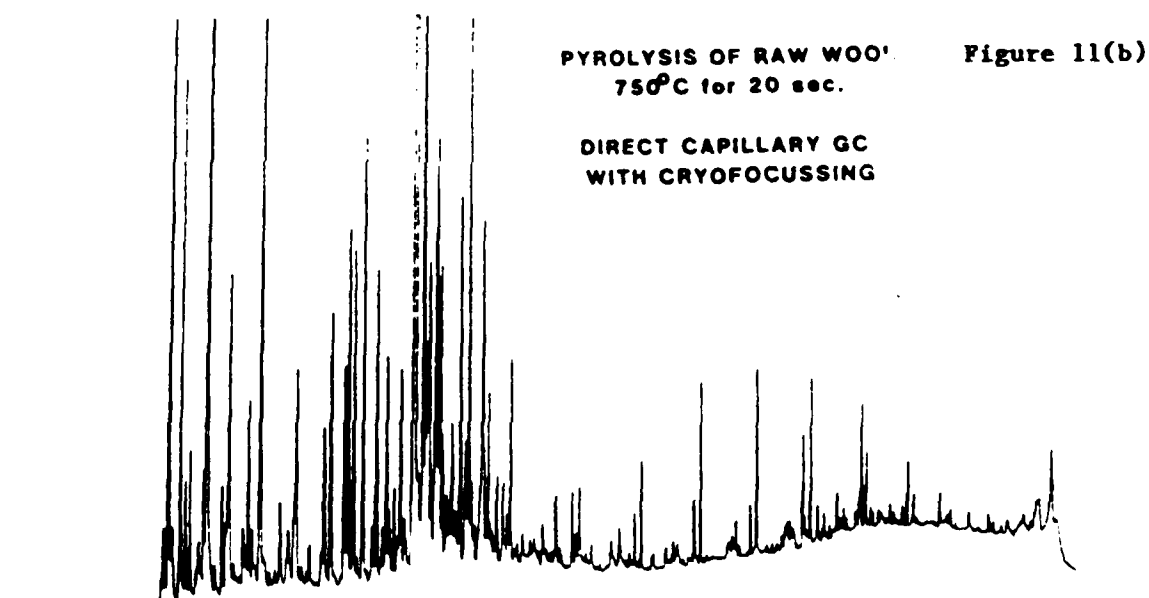




Figure 11(c)
(Microorganisms)

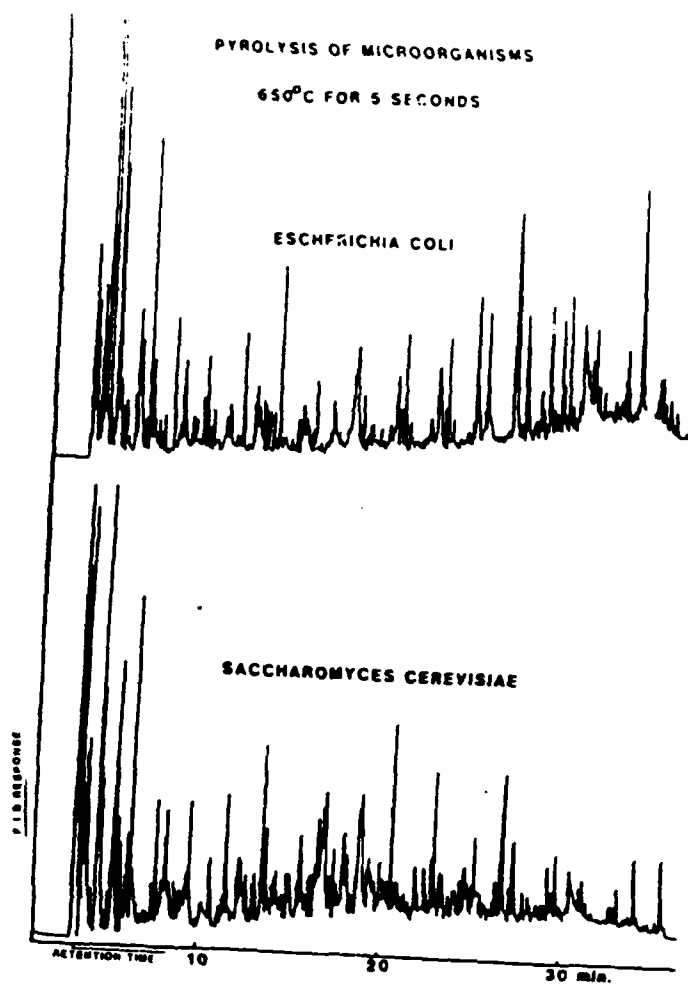

(Silk) Figure 11d)

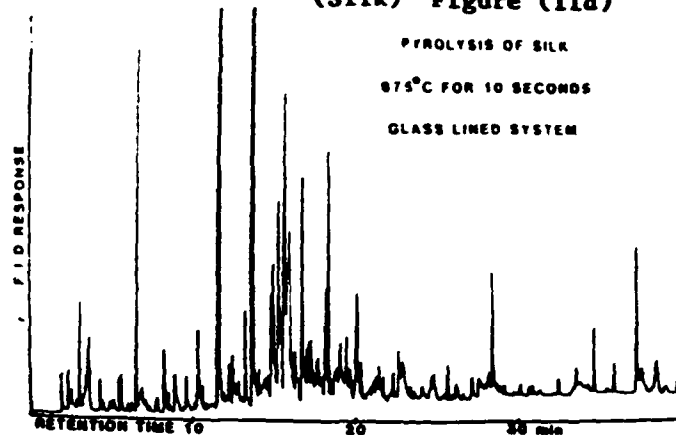

(Cellulose) (Figure 11e)

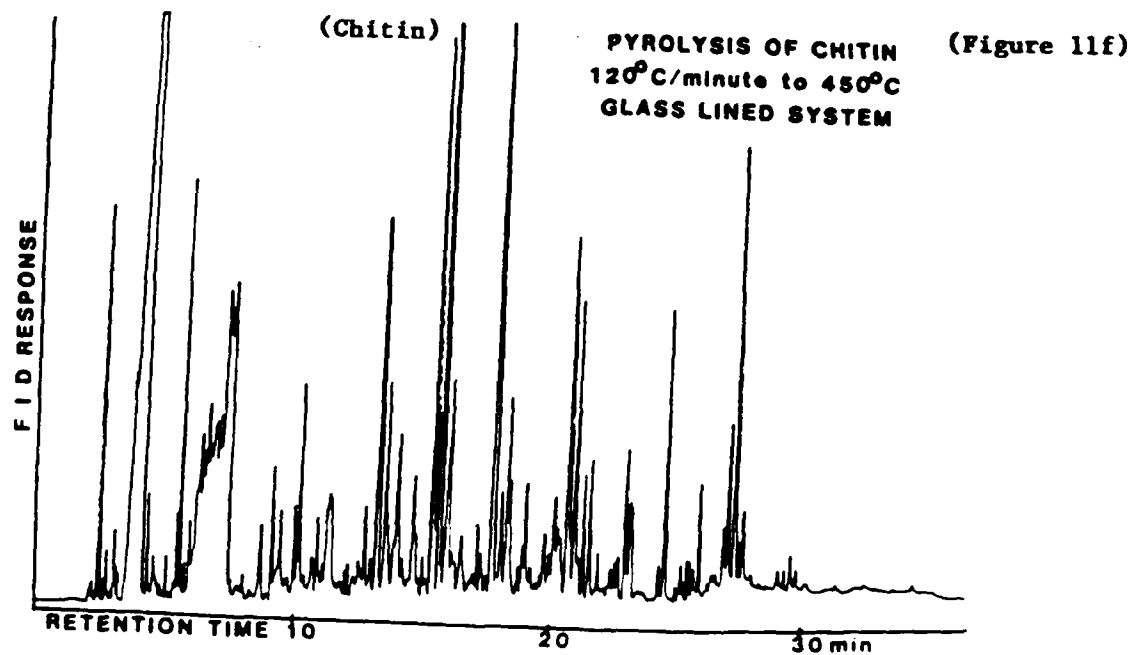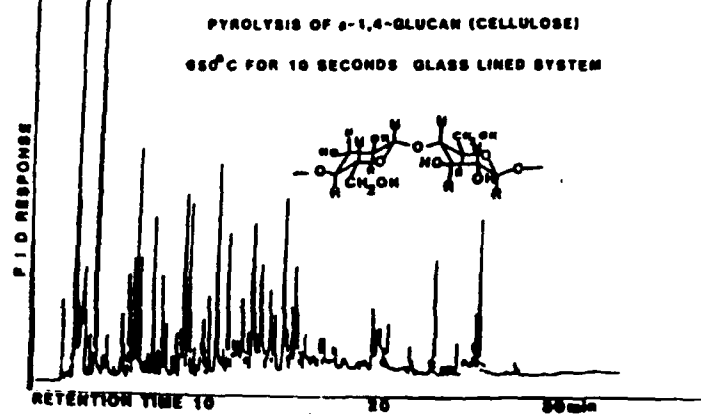

**BOX-BEHNKEN EXPERIMENTAL DESIGN
USED IN ANALYTICAL PYROLYSIS METHOD DEVELOPMENT**

Figure 12

- MODIFIED TURBOPASCAL (3⁵) PROGRAM
FOR MICROCOMPUTERS
- ESTABLISH EXPERIMENTAL ERROR IN
MULTIVARIATE PARAMETER STUDY
 - PYROLYSIS TEMP
 - PYROLYSIS MODE
 - PRETREAT TEMP
 - MATRIX
 - ANION TYPE

BOX-BEHNKEN EXPERIMENTAL DESIGN

Figure 13

SFE METHOD DEVELOPMENT

FIVE-FACTOR THREE-LEVEL (3⁵) OUTLINE

<u>INITIAL PRESSURE</u> (PSI)	<u>TIME</u> (MIN)	<u>TEMP</u> (C)	<u>RATE/FINAL</u> (PSI/MIN)	<u>MicroCAP</u>
1000	10	R.T.	100/3500	1
1200	8	100	200/4000	2
1500	6	120	300/4000	3
 <u>EXPERIMENTAL</u>				
#1 1200	8	100	200/4000	2
#2 1200	8	100	200/4000	2
⋮				
#16 1500	8	120	200/4000	2
⋮				
#22 1200	8	100	100/3500	3
⋮				
#36 1200	8	120	200/4000	1
⋮				
#46 1200	10	100	100/3500	2

General Processes Involved in Organic Analysis

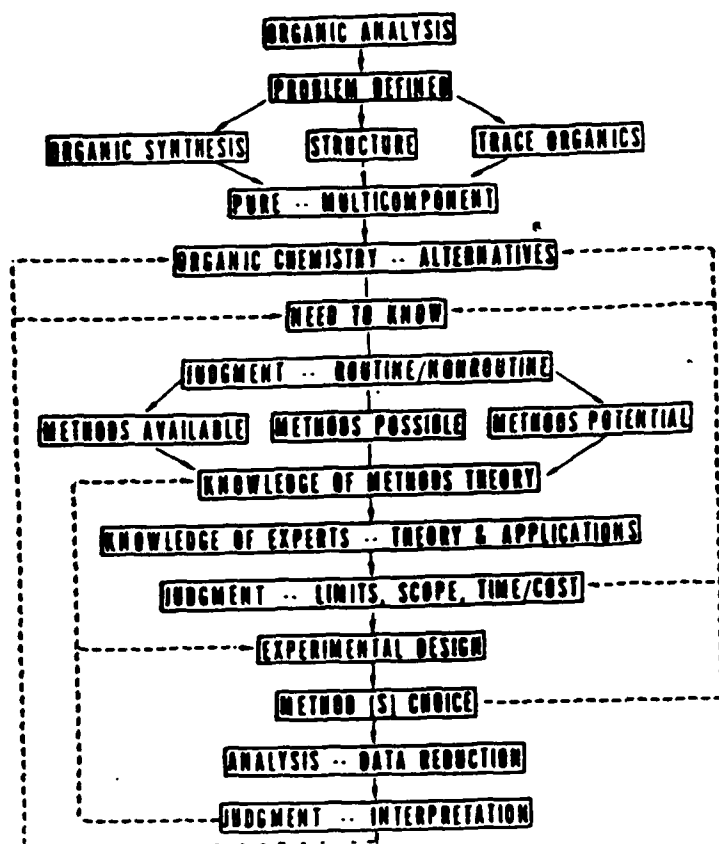


Figure 14


References

1. LIEBMAN, S. A., et al., *J. Poly. Sci. Part A-1*, 9, p. 1823, 1921 (1971).
2. LIEBMAN, S. A., et al., *Macromol.* 4, 134 (1971).
3. LEVY, E. J., et al., papers presented at Eastern Analytical Symposium, New York, Nov. 1971, and Pittsburgh Conference on Analytical Chemistry and Applied Spectroscopy, Cleveland, Ohio, March 1972.
4. HARRIS, W. E., *Anal. Chem.* 42, no. 13, 53A (1970).
5. BANKS, C. V. and BECK, R. F., *Anal. Chem.* 42, no. 14, 56A (1970).
6. LONG, F. A., *Science* 171, 962 (1970).
7. CRAM, S. P., *Research/Development*, 16 (July 1970).

MULTILEVEL EXPERT SYSTEMS OUTLINE

Figure 15(a)

TOTALLY AUTOMATED EXPERIMENTAL DESIGN AND DECISION-MAKING

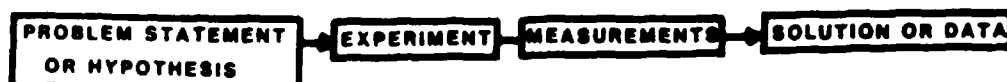

LABORATORY AUTOMATION USING EXPERT SYSTEMS DRIVERS

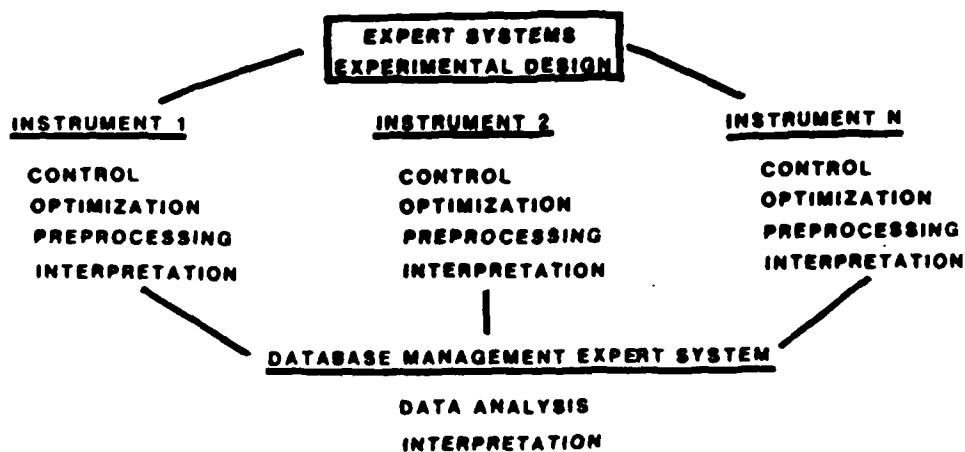


Figure 15(b)

EIMAT

A LINKED NETWORK OF EXPERT SYSTEMS
WITH PATTERN RECOGNITION AND SEARCH PROGRAMS
FOR MATERIALS CHARACTERIZATION

- | <u>COMPONENTS</u> | <u>ATTRIBUTES</u> |
|-------------------------------|---|
| 1. <u>DATABASE MANAGEMENT</u> | A. STORAGE OF PARAMETERS AND DATA ON SAMPLES FOR SELECTED INSTRUMENTAL TECHNIQUES |
| | B. RETRIEVAL OF SELECTED SAMPLES FORMING A DATA SET FORMATTED FOR MULTIVARIATE ANALYSIS |
| | C. CREATE, ADD, DELETE, HELP AND SHOW FUNCTIONS |

2. EXPERT SYSTEMS AND EMBEDDING SUBPROGRAMS-TIMM

- | |
|--|
| A. FORTRAN SOURCE CODE |
| B. EMBEDDING OF TIMM SYSTEM WITHIN USER PROGRAMS |
| C. CAPABLE OF HANDLING METRIC AND NON-METRIC INFORMATION |
| D. HEURISTIC DESIGN |

Figure 15(b)
(continued)

COMPONENT ATTRIBUTES

- | | |
|--|---|
| 3. <u>PATTERN RECOGNITION EXPERT SYSTEM - EXMATH</u> | A. HEURISTIC DESIGN |
| | B. SUPERVISED AND UNSUPERVISED PATTERN RECOGNITION, FACTOR ANALYSIS, PLOTTING |
| | C. EXPERTISE INCLUDES DATA PREPROCESSING AND EVALUATION OF RESULTS |
| | D. USER INTERVENTION FOR DATABASE MODIFICATIONS |
| | E. IMPLEMENTABLE AS JACKKNIFING PROCEDURE |

- | | |
|--|---------------------------|
| 4. <u>SPECTRAL SEARCH AND MATCH ALGORITHMS</u> | A. PAIRS INFRARED SPECTRA |
| <u>PARTIAL INTERPRETATION AIDS</u> | B. PSM-MASS SPECTRA |

DEVELOPMENT OF EXMATH Figure 15(c)

A LINKED NETWORK OF EXPERT SYSTEMS FOR MATERIAL ANALYSIS

ES #1 ANALYTICAL STRATEGY FOR DEFINED PROBLEM

ES #2 INSTRUMENTAL CONFIGURATION/CONDITIONS

ES #3 DATABASE GENERATION

ES #4 DATA TREATMENT

ES #5 DATA RESULTS

ES #6 DATA INTERPRETATION

ES #7 ANALYTICAL GOAL

Figure 15(d)

PYROLYSIS/CONCENTRATOR WITH
CAPILLARY GC-MS/DS

FORMATTED FOR ENTRY INTO
DATABASE MANAGEMENT SYSTEM,
EXDDM

SAMPLE ID = 0001

SAMPLE NAME = EA4004

SAMPLE REP. NO.=1

SAMPLE EXP. NO.=1

PYROPROBE 124-BP 5985B GC-MS

ANALYSIS TYPE = DBS/PGC-MS

TYPE = CORR. AREA, % OF TOTAL

SAMPLE FILE BASIS = S0087

NO. VARIABLES = 16

20
87/28/7
600C/20 SEC
175C/3 MIN
-50/4 MIN
-50/4 MIN
290C/3 MIN
290C/3 MIN
300/7 MIN
DB-5 (30M .25MM)
10 PSI
50C (3M) TO 300C
30C/MIN
EI (70EV)/3000A
200C
33 TO 500/SEC

PEAK1=6.285
15012
1.781
PEAK2=6.835
11391

Figure 15(e)

FACTOR ANALYSIS

Test Data for Chemometrics
FACTOR 1 VS FACTOR 2


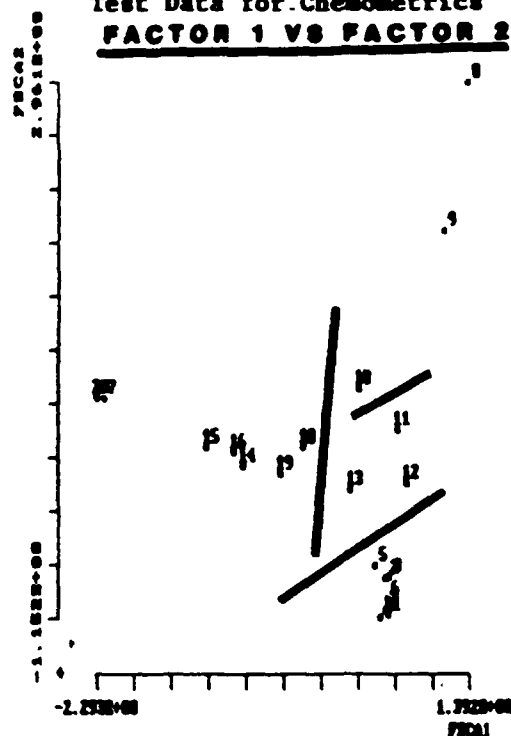


Figure 15(f)

FACTOR ROTATION

Test Data for Chemometrics
PLOT OF COVARIANCE MATRIX

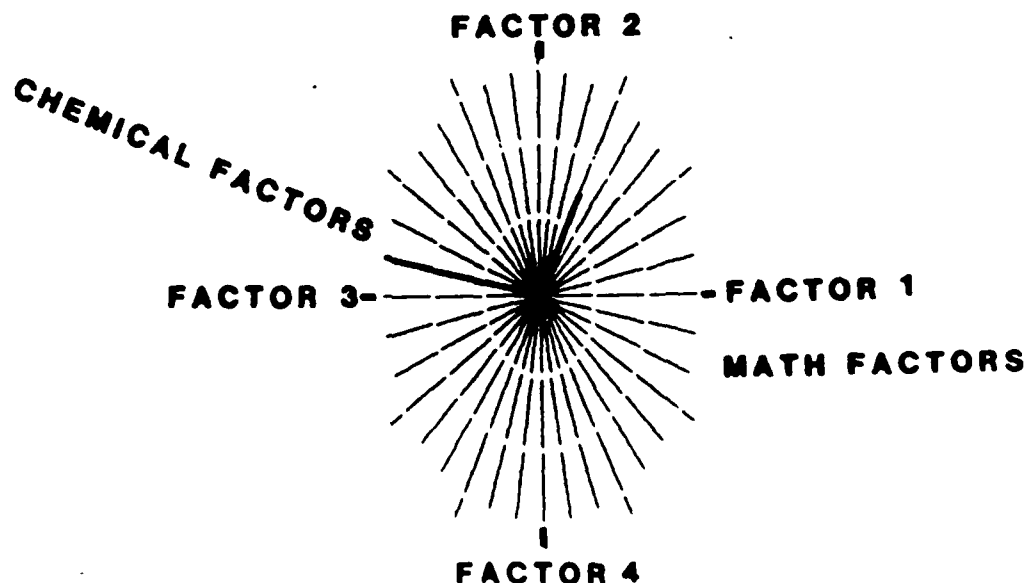

The allowable responses are shown below.

MAIN MENU		B Build	T Train	E Exercise
Type letter		I Inquire	G Generalize	C Compress
+ file name		M Modify	N Check consistency	P Check completeness
and return.		V Convert	R Recover	S Stop

Selection? e

First, some preliminaries ...

Do you want to record this in a print file? n

Tell me the names of all the files containing the expert systems you want to exercise. Enter one per line, and enter a "/" as the last line:

EXCONFIG
strategy
config
datgen1
datarea
datareal
interpt1
angoall
microcap ←

Shall I pick the test cases for you (Y or N) ? n

Shall I narrow down the possible choices as I go (Y or N) ? n

Do you want me to use terse (T) or verbose (V) mode? t

Ok, we're ready to start ...

SCOPE	IS screen
SAMPLE AMT	IS mg
SAMPLE FORM	IS powder
SAMPLING PROCESS	IS *
SAMPLE HISTORY	IS *
INSTR. AVAIL	IS all
ANALY STRATEGY	IS FTIR/SYS2(80) GC/SYS1(20) (Reliability = 100)

EXMAT - HELP "CHOICE SYS 8/CT"

FIGURE 17

TAGA 6000 Triple Quadrupole MS/MS Interfaced with Pyroprobe Model 122

CHOICE "SYS8/CT"

Operational settings for the Pyroprobe are under condition C; for the triple quadrupole Sciex TAGA 6000, under condition T. When using the CDS Pyroprobe Model 120 or 122, the probe is positioned in the glass tube support about 1.5 cm from the end with the metal outside rods just near the glass. The corona discharge needle in the APCI source should be ca. 2 mm from the inside quartz tube insert of the Pyroprobe which holds the samp. Carrier gas stream or air or N₂ of 80 ml/min sweeps the volatiles into the ionization chamber. Condition T: The TAGA may be operated in the MS-mode or MS/MS-mode, with each parameter set to its optimum in terms of peak resolution, shape, and transmission of the ions thru the mass filters of quad 1-quad 3. Recommended values for the conventional MS quad 1 scan parameters are given in (....) and for the daughter-ion mode (MS/MS) in [....]. MTL (sec) measurement time limit: (0.0100)/(0.0300); STEP (amu) mass increments: (1.00000)/(1.0000); PE (sec) delay period between scans: (0.00)/(0.00); C = constant time mode: (10)/(1); TH (ions/sec) Threshold (10)/(1); TT (min.) total time: (2.00)/(2.00); MO (amu) mass offset: (-0.1)/(-0.1); DI (uA) discharge current: (2.000)/(2.500); IN (volt) interface plate voltage: (650)/(650); L1 (volt) Lens 1: (63)/(70); CB (volt) Lens 2 (40)/(50); L3 (volt) Lens 3a and lens 3b: (40)/(50); L4 (volts) lens 4: (5.1)/(40); L5 (volt) lens 5: (-250)/(-250); L6 (volt) lens 6: (25)/(25); R1 (volt) rod offset of quad 1: (30)/(45); RE1 (volt) mass resolution of quad 1: (153)/(135); R2 (volt) rod offset of quad 2: (15)/(0); CG (volt) collision gas: (off)/(Ar); L10 (volt) lens 10: (60)/(60); MU (volt) multiplier voltage: (-3650)/(-4200); R3 (volt) rod offset of quad 3: (15)/(auto); RE 3 (volt) mass resolution of quad 3: (120)/(153). Parameters SY, SE, PE, G1, G2, G3, BAF, DC, MC1, DM1, Q2r, and DM3 are not important in the above compilation.

Box-Behnken Design for SPE

Block 1

Figure 18

	Initial Pressure (psi)	Time (min)	Temp	Rate/Final	Microcap
1	1000	10	R.T.	100/3500	1
2	1200	8	100	200/4000	2
3	1500	6	120	300/4500	3
4	1200	8	100	200/4000	2
5	1500	6	100	200/4000	2
6	1000	10	100	200/4000	2
7	1000	6	100	200/4000	2
8	1200	8	120	300/4500	2
9	1200	8	120	100/3500	2
10	1200	8	R.T.	300/4500	2
11	1200	8	R.T.	100/3500	2
12	1200	10	100	200/4000	3
13	1200	10	100	200/4000	1
14	1200	6	100	200/4000	3
15	1200	6	100	200/4000	1
16	1500	8	120	200/4000	2
17	1500	8	80 R.T.	200/4000	2
18	1000	8	120	200/4000	2
19	1000	8	80 R.T.	200/4000	2
20	1200	8	100	300/4500	3
21	1200	8	100	300/4500	1
22	1200	8	100	100/3500	3
23	1200	8	100	100/3500	1

Box-Behnken Design for SPE

Block 2

Figure 19

	Initial Pressure (psi)	Time (min)	Temp	Rate/Final	Microcap
24	1200	8	100	200/4000	2
25	1200	8	100	200/4000	2
26	1200	8	100	200/4000	2
27	1200	6	120	200/4000	2
28	1200	6	80 R.T.	200/4000	2
29	1200	10	120	200/4000	2
30	1200	10	80 R.T.	200/4000	2
31	1500	8	100	300/4500	2
32	1500	8	100	100/3500	2
33	1000	8	100	300/4500	2
34	1000	8	100	100/3500	2
35	1200	8	120	200/4000	3
36	1200	8	120	200/4000	1
37	1200	8	80 R.T.	200/4000	3
38	1200	8	80 R.T.	200/4000	1
39	1500	8	100	200/4000	3
40	1500	8	100	200/4000	1
41	1000	8	100	200/4000	3
42	1000	8	100	200/4000	1
43	1200	6	100	300/4500	2
44	1200	6	100	100/3500	2
45	1200	10	100	300/4500	2
46	1200	10	100	100/3500	2

Figure 20(a)

TIME-RESOLVED ANALYTICAL PYROLYSIS APCI MS/MS SCHEMATIC

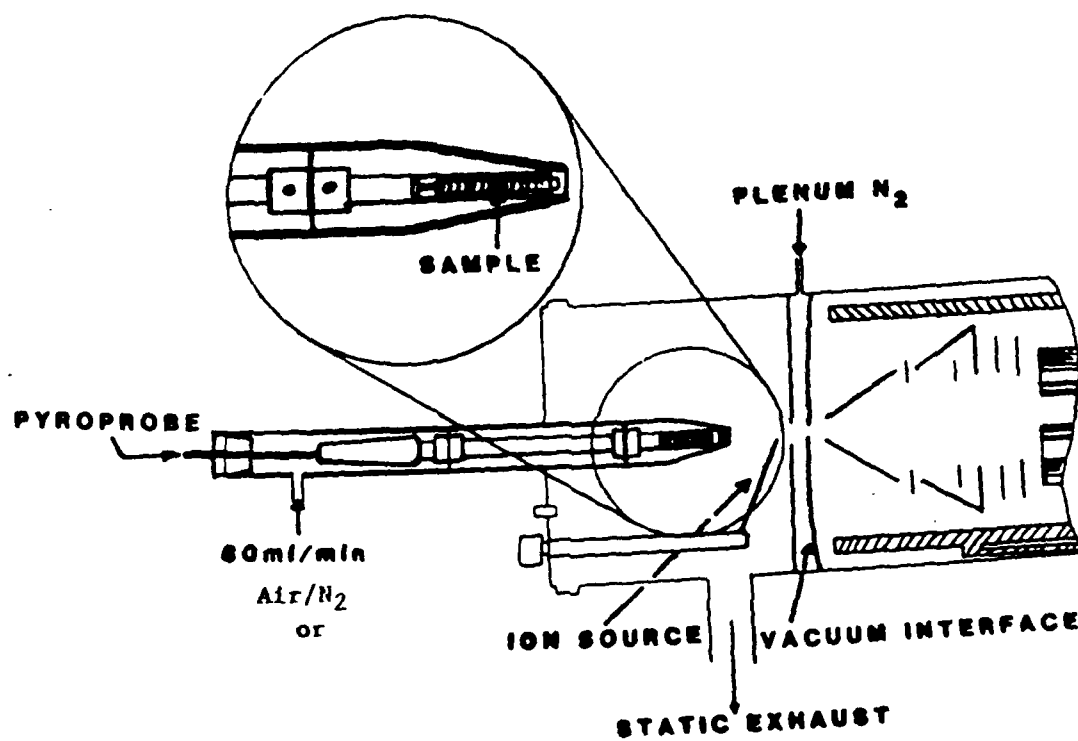


Figure 20(b)

Pyrogram Pulse 7000C/20 sec Polyethylene (Air, N₂), Oxidised PE (Air, N₂) and PEO (Air, N₂)

PE $\{CH_2\}_n$

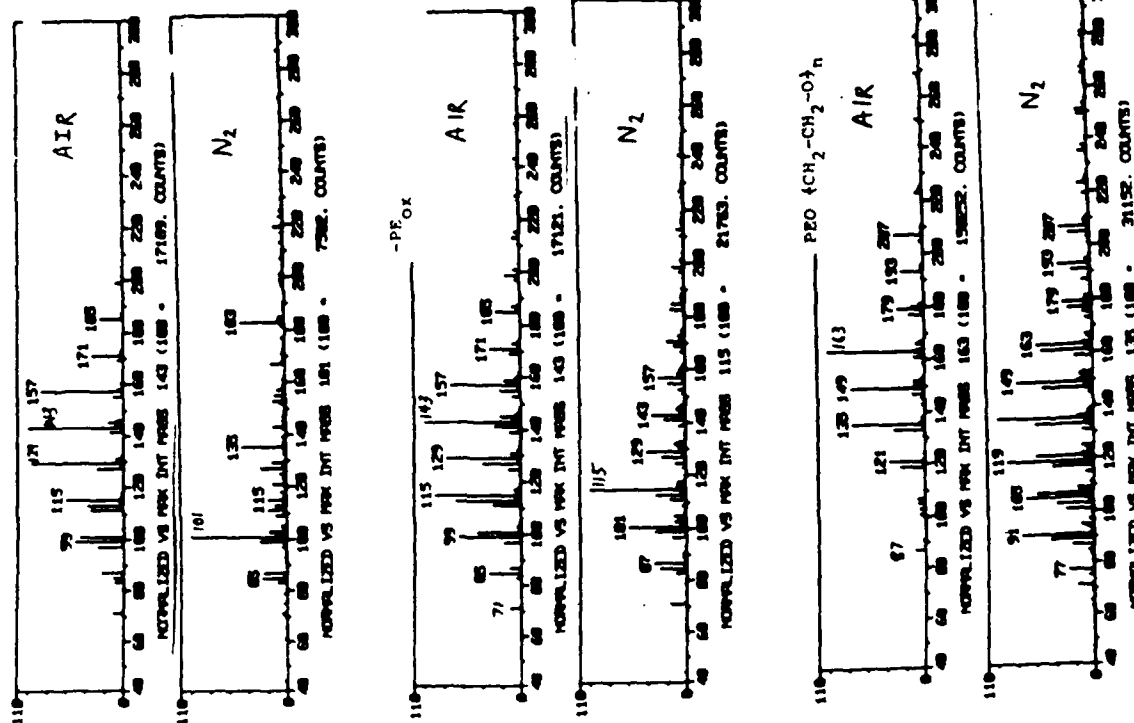


Figure 21(a)

Pyrogram - Microcep 91 (1800/min in 400°C) No Heat

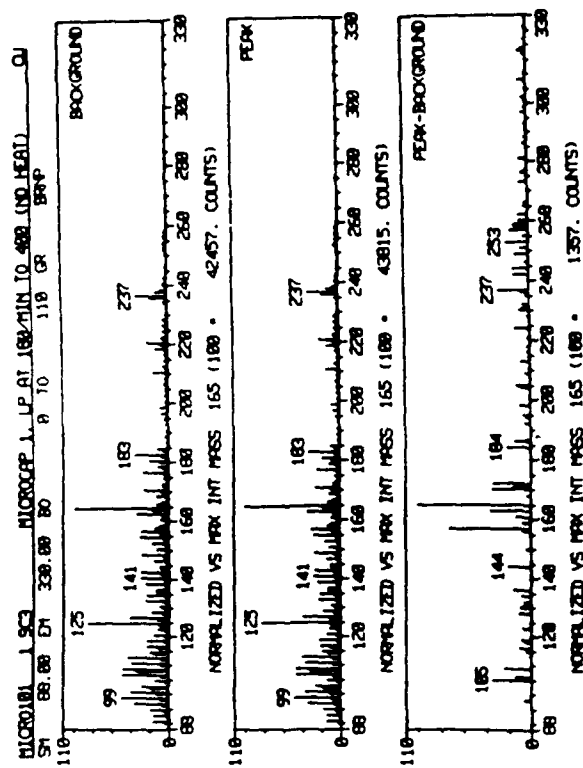


Figure 21(b)

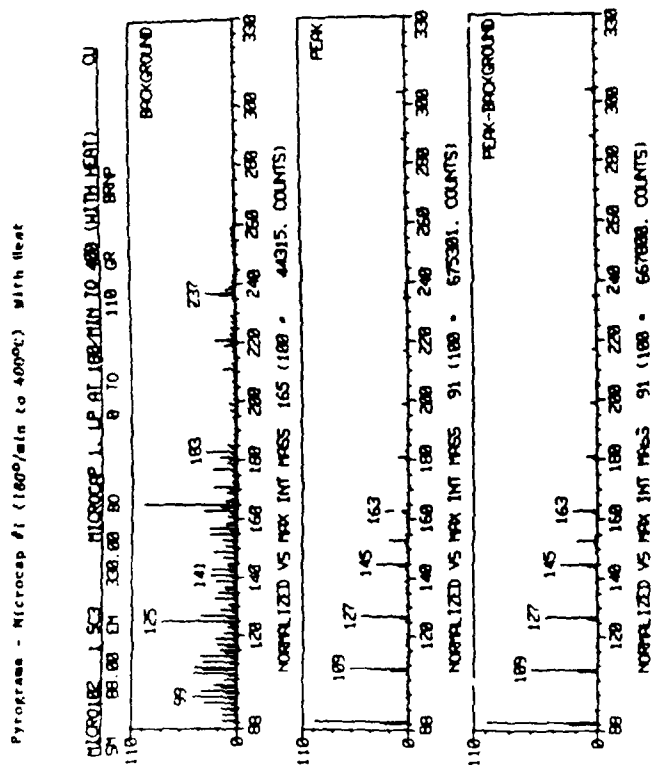


Figure 21(c)

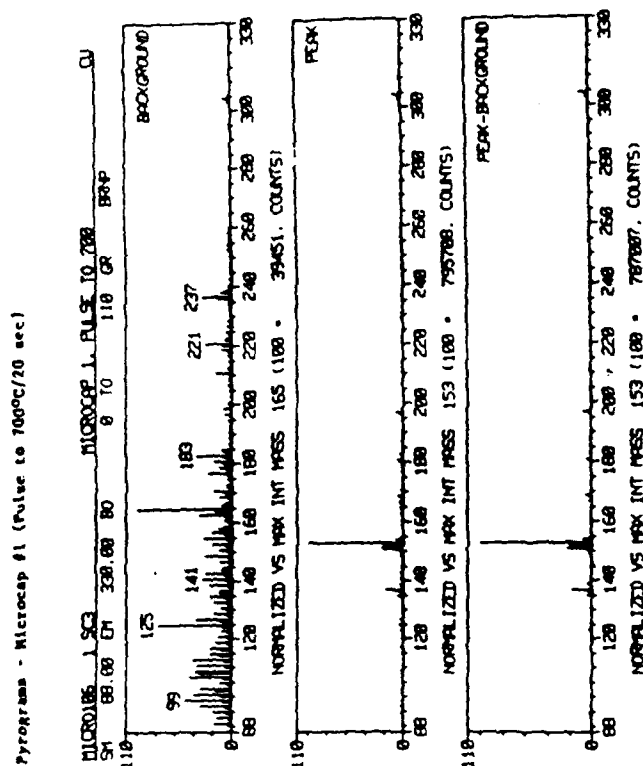


Figure 21(d)

Pyrogram - Microcap #1 (240°C/min to 500°C)

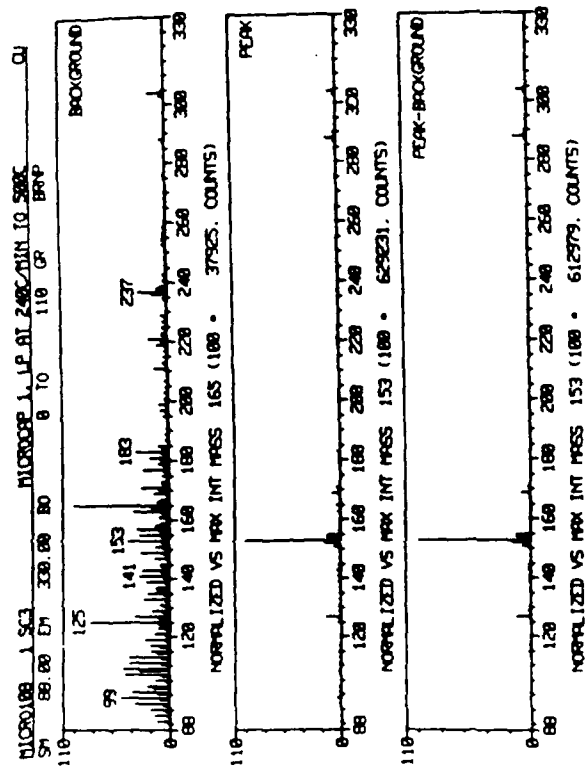


Figure 21(e)

Pyrogram - Microcap #1 (300°C/min to 200°C), 4 Min Interval

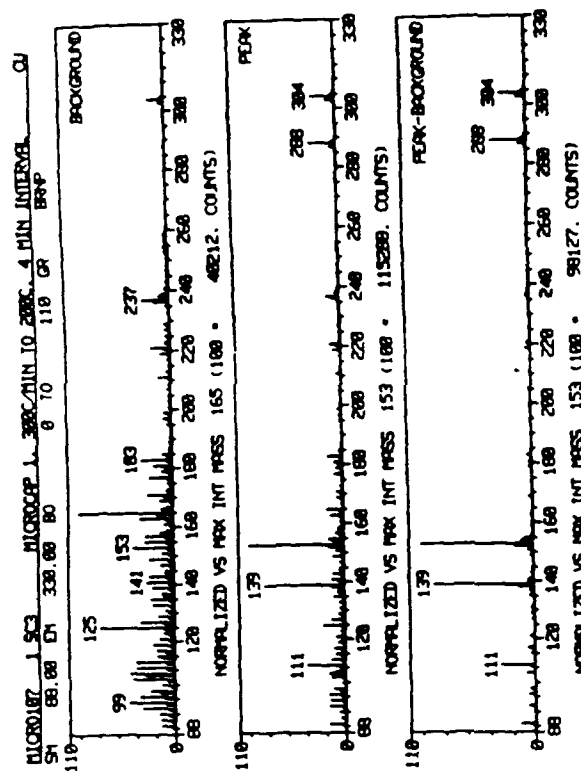


Figure 22

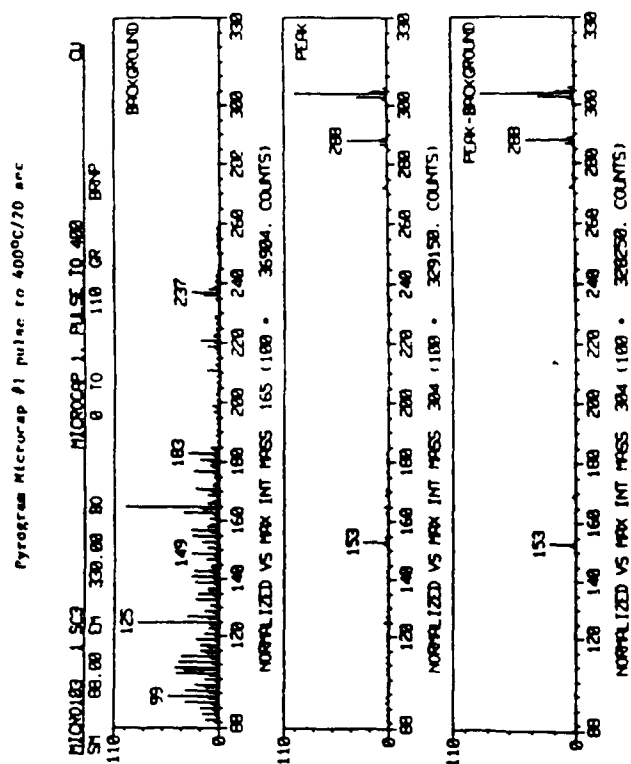


Figure 23(a)

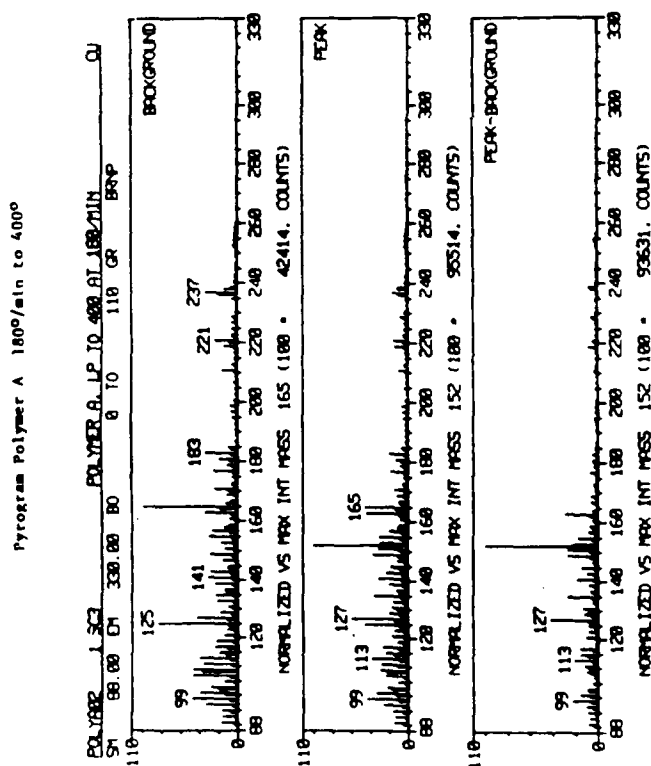


Figure 23(b)

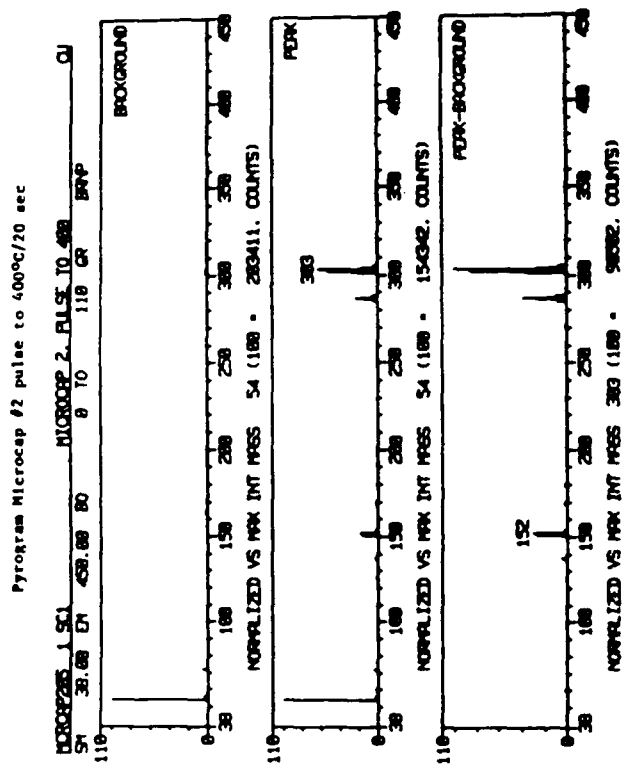


Figure 23(c)

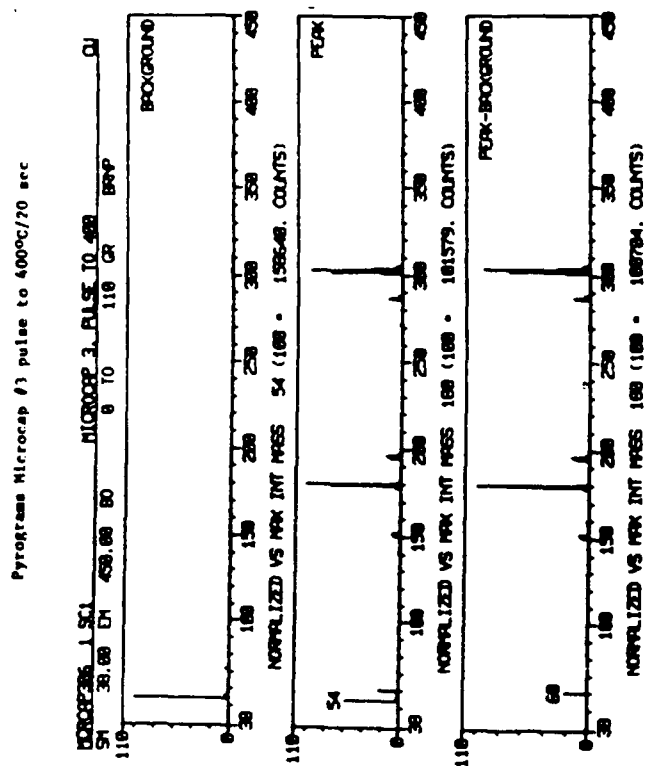


Figure 23(d)

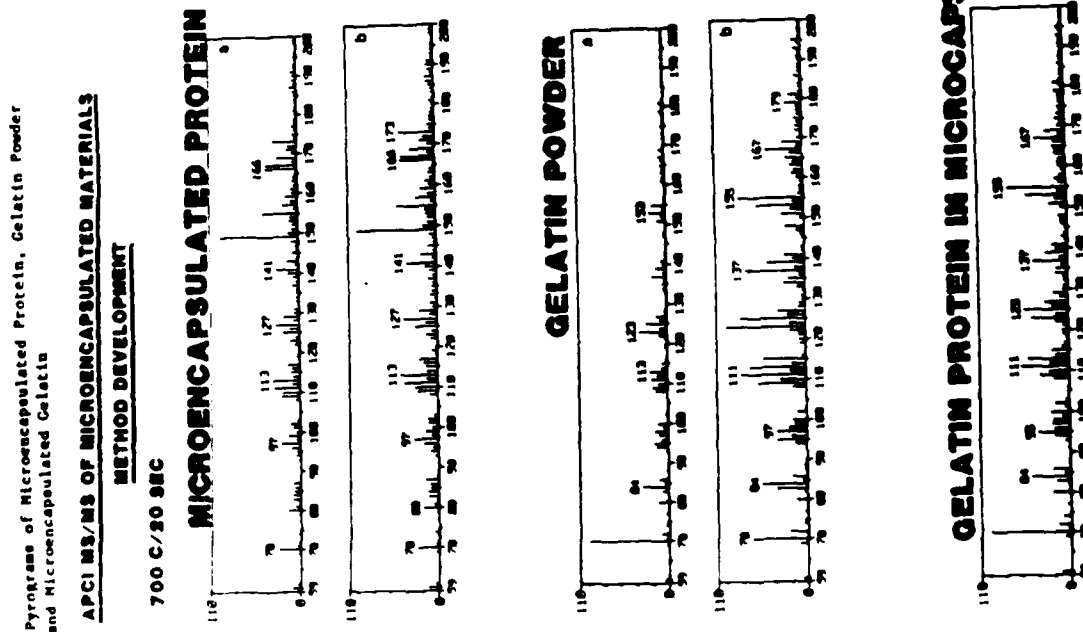


Figure 24(a)

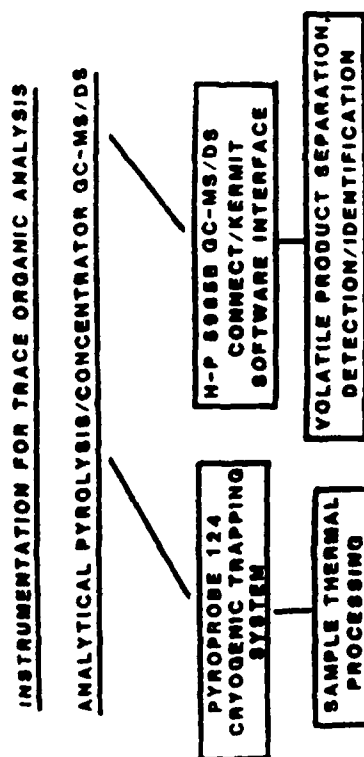


Figure 24(b)

Pyrograms - Concentrator/GC-MS (HP 5985B) Data Generation -

Groups I, II, III Organic Amine Salts

DATA GENERATION

PYROGRAMS (TIC)

GROUP I

GROUP II

GROUP III

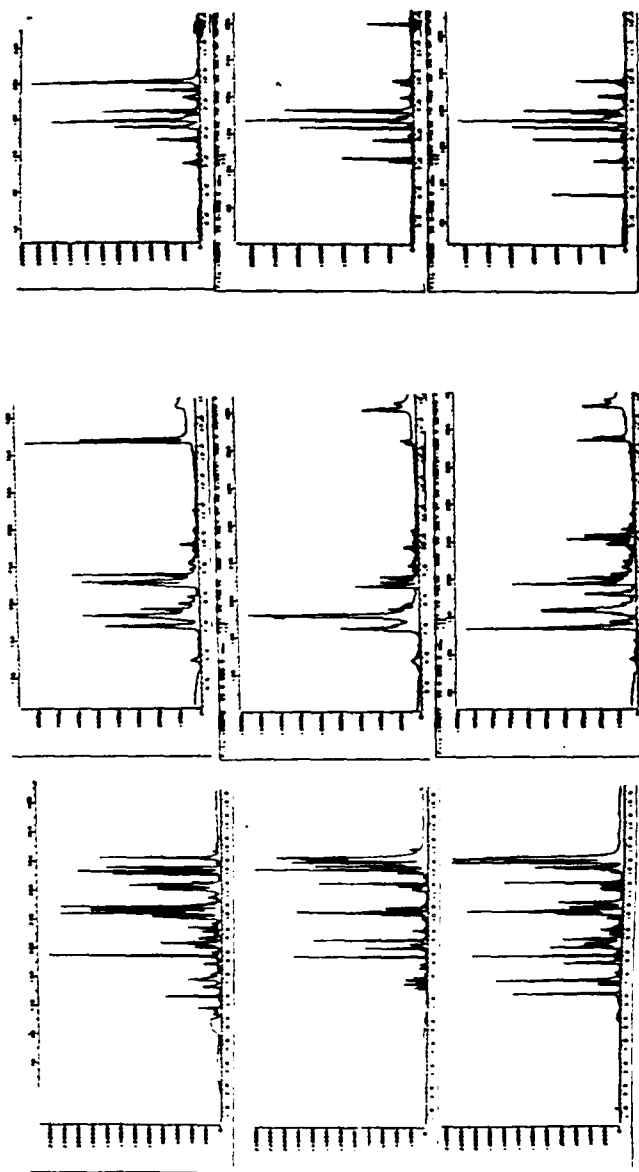


Figure 24(c)

MS Patterns of TIC Peaks from Pyrograms and Library Match

DATA GENERATION

MS PATTERNS OF TIC PEAKS

PYROGRAM SAMPLE 4

PEAK 1

GROUP 1

**IDENTIFICATION BY
LIBRARY MATCH**

MS SPECTRUM

PEAK 1


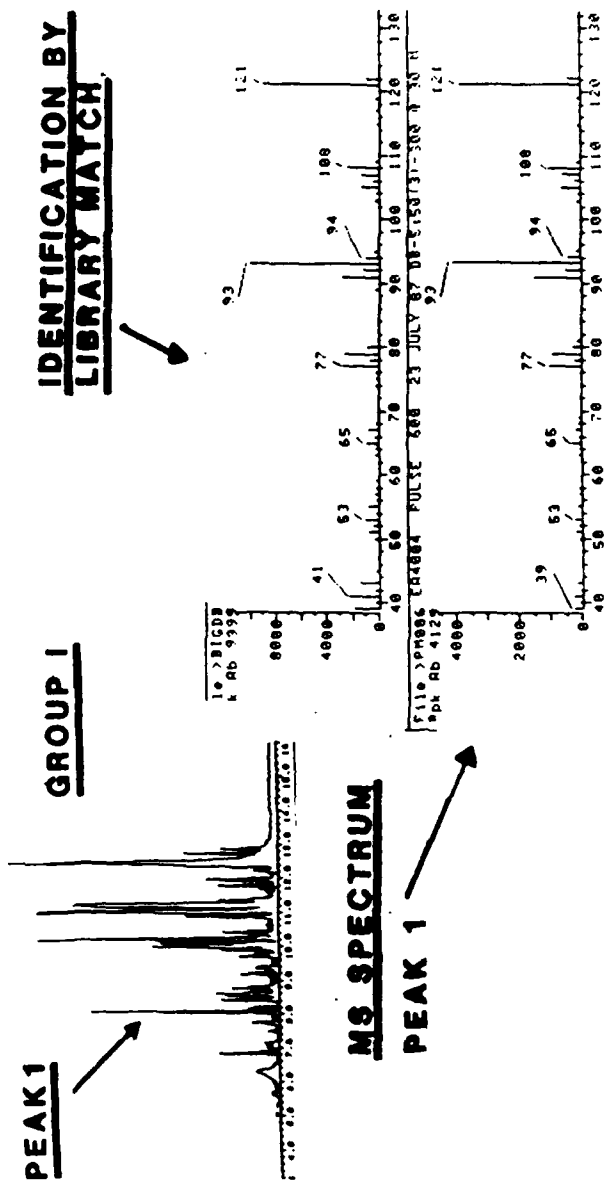


Figure 25


AMINE SALT STRUCTURES FOR TEST METHOD

(EXMAT SERIES)


GROUP I

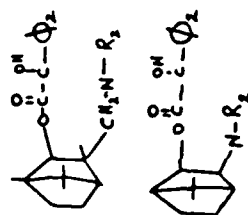

GROUP II

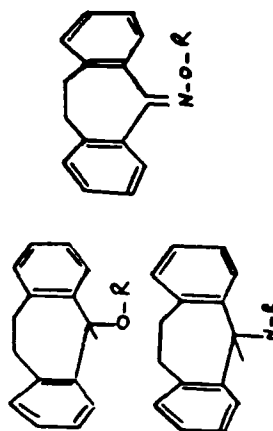

GROUP III

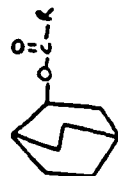


o HCl

o FUMARATE

o MALEATE


High Density PE, 800°C/20 sec Pyrogram - RIC

Figure 26(a)

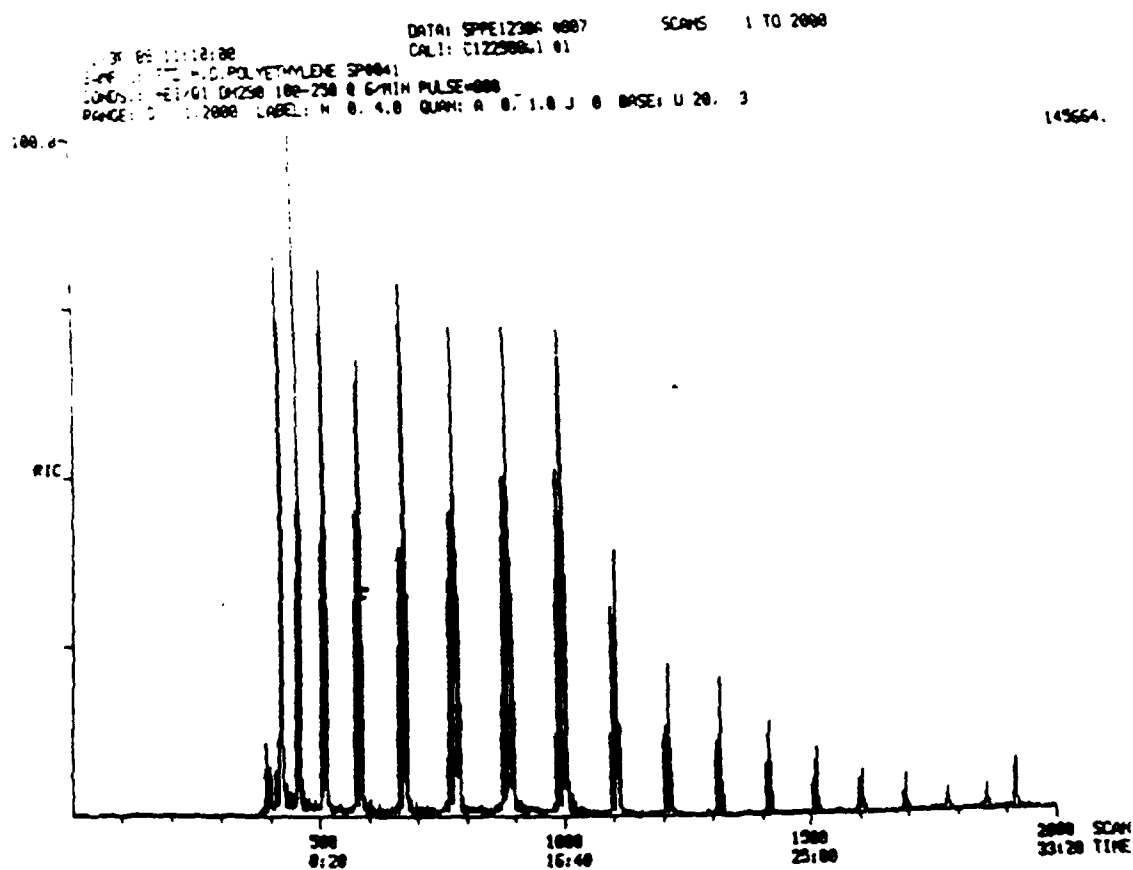

Scansets 650-950, Carbon Triplets ~C₁₃-C₁₅ of 26(a)

Figure 26(b)

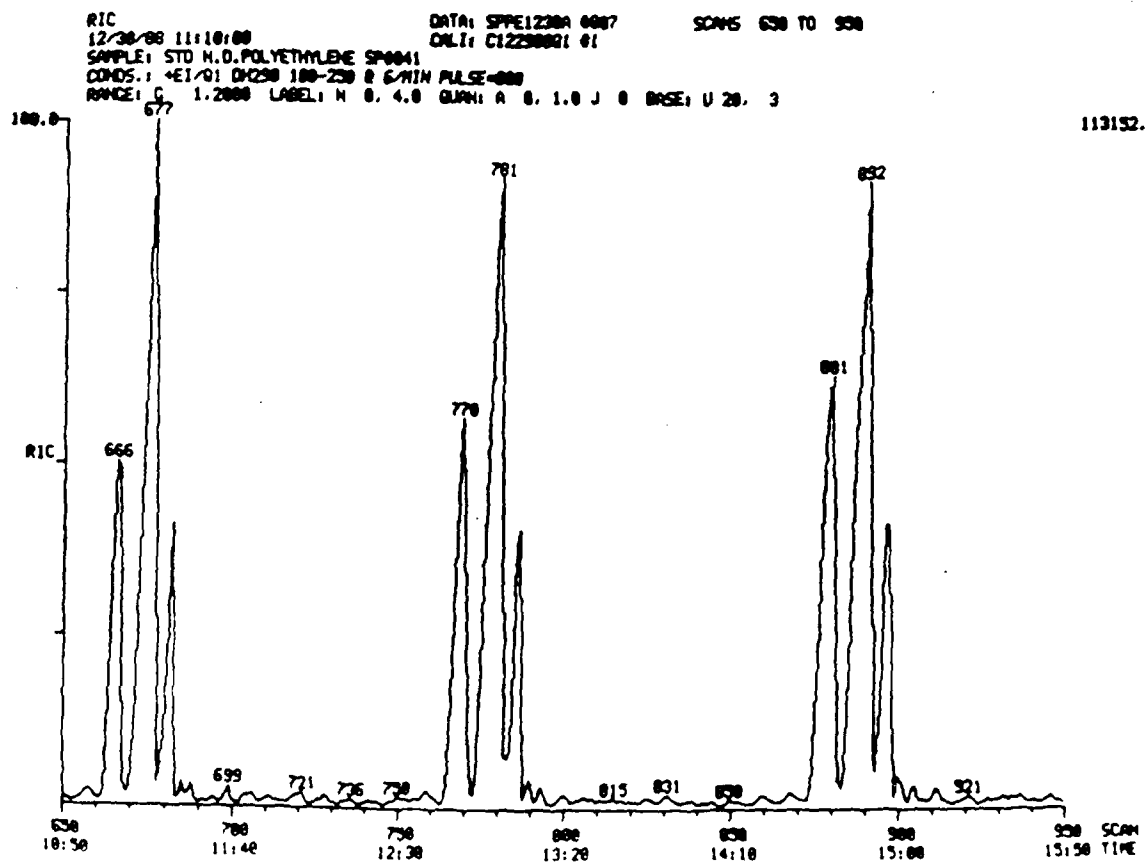



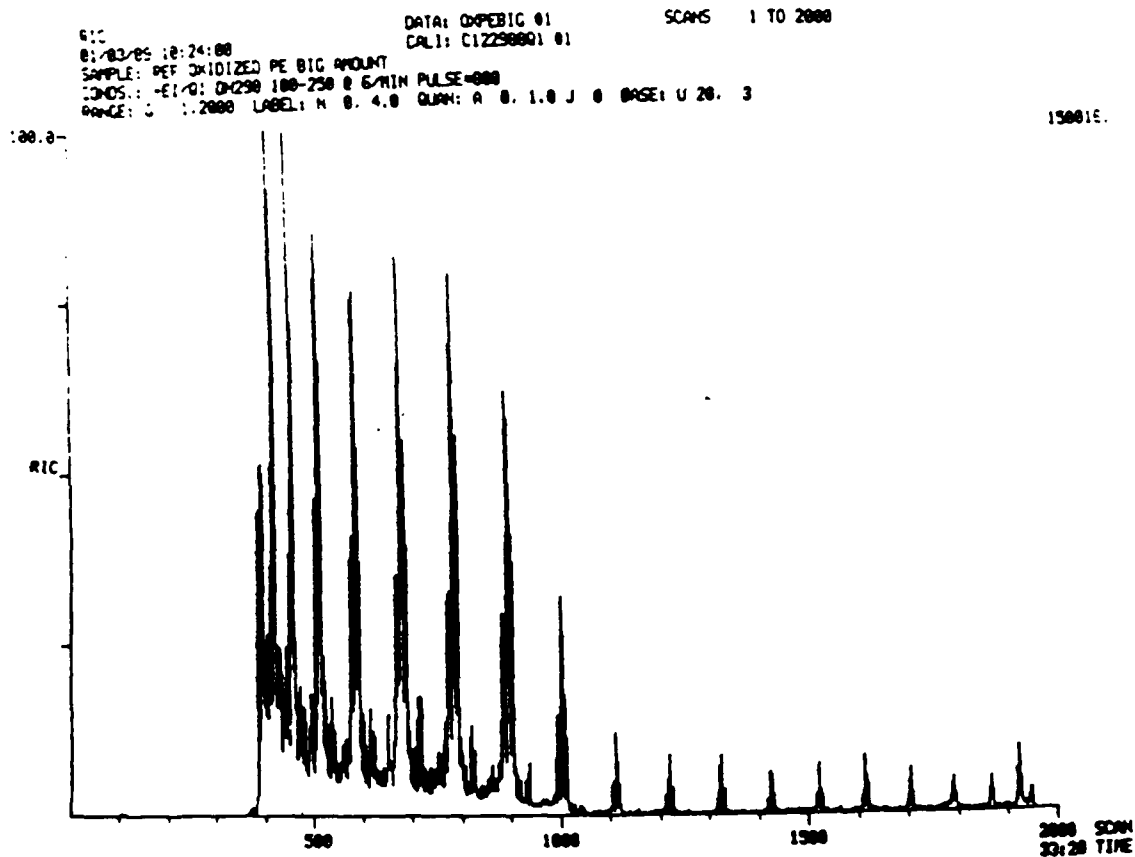
Scansets 650-950, ~C₁₃-C₁₅ of 26(c)

Figure 26(d)

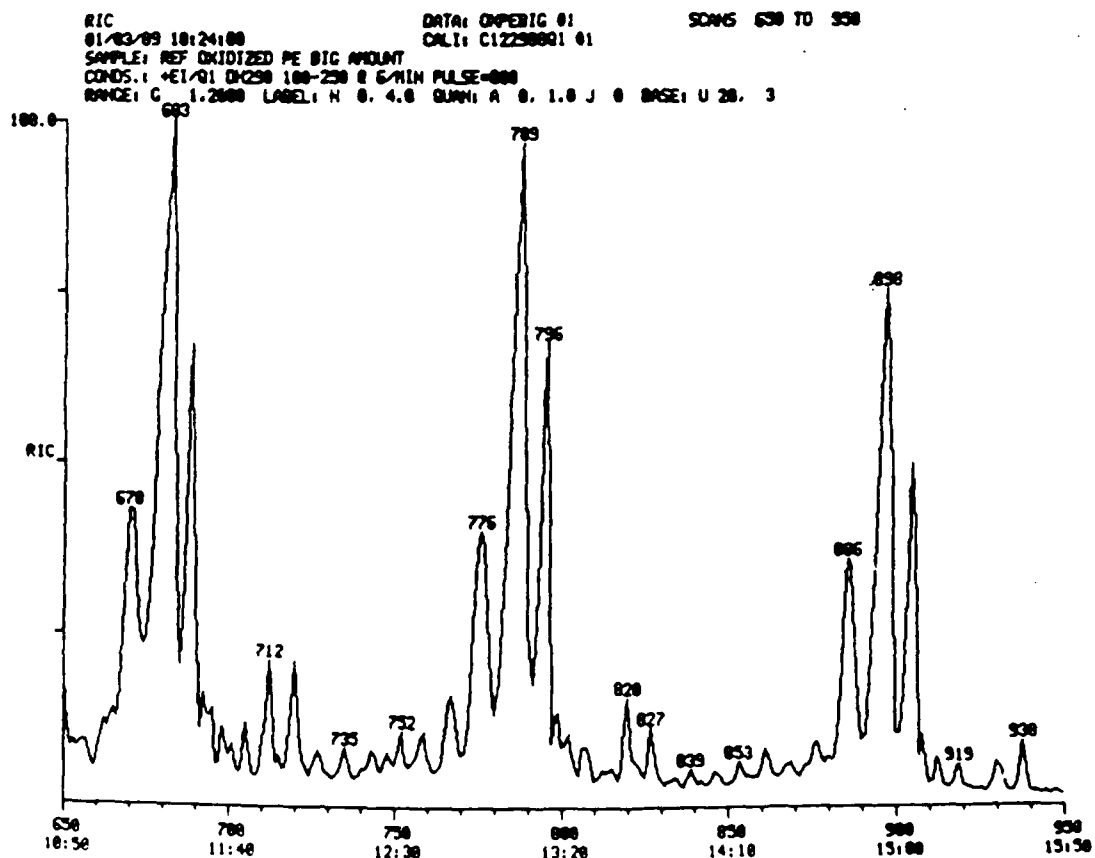

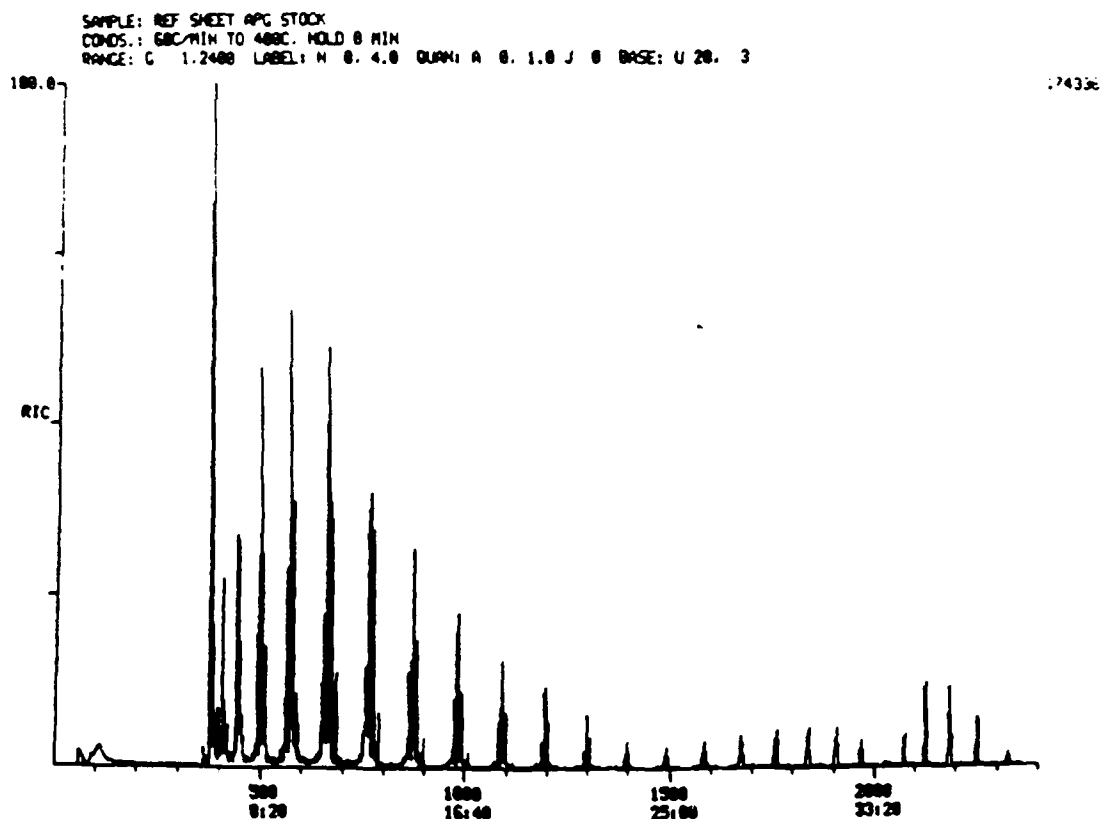

Scansets 650-950, ~ C₁₃-C₁₅ of 26(e)

Figure 26(f)

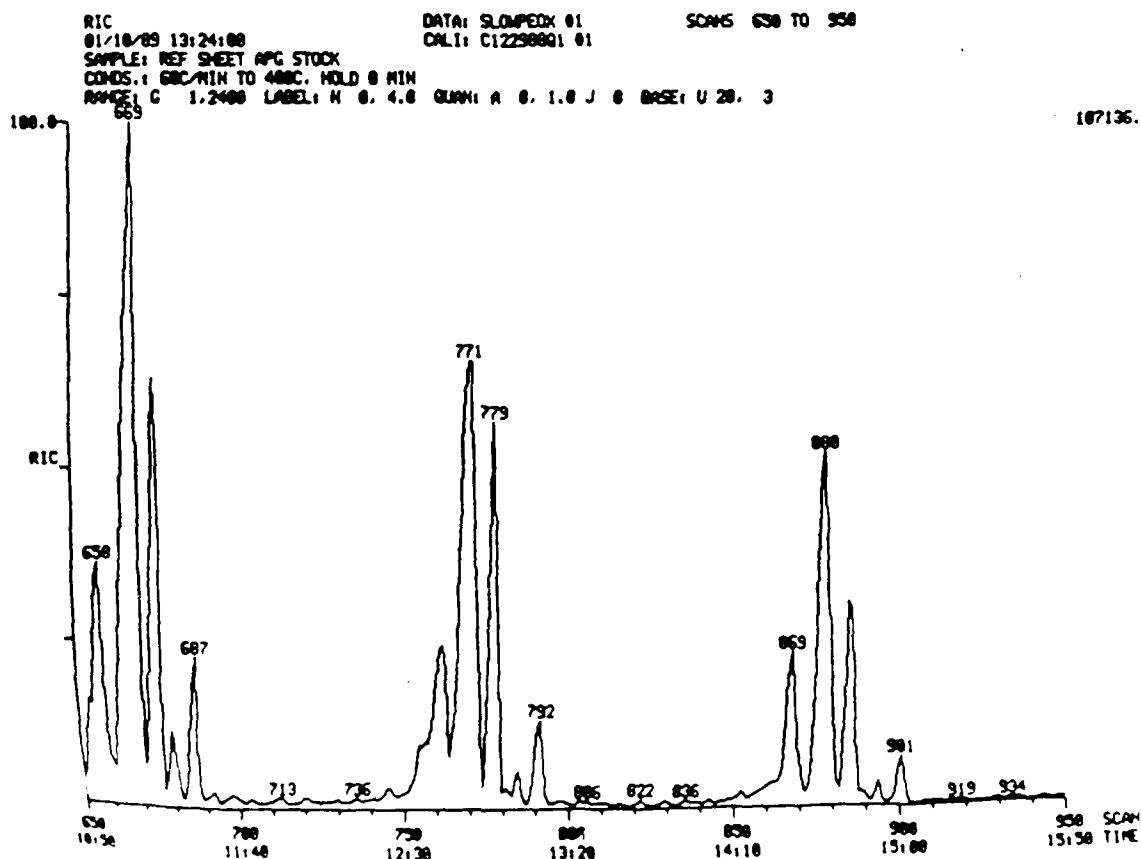


Figure 27

GC-MS Experimental Conditions (Finnigan TSQ GC/MS)

```

12/30/88 11:10:21 SCAN 1 OF 2000
ACQUISITION STARTED

ACQUIRE 12/30/88 11:10:00 + 0.03 RUN 0 SPPE1230A ACQUIRING
FREE SECTORS 20966 SCAN 3 OF 2000
SAMPLE STD M D POLYETHYLENE SP#041
CONDS +E1/Q1 DM290 100-250 @ 6/MIN PULSE=800
FORMULA - INSTRUMENT Q1 WEIGHT 0.000
SUBMITTED BY MW ANALYST MW ACCT NO P687

***** GC PARAMETERS *****
LOADED GC DESC PE CURRENT GC OVEN TMP 100 DEGC INJECTOR : 280 DEGC
CURRENT GC DESC PE GC ELAPSE TIME 0.7 MIN INT OVEN : 280 DEGC
SEQ # TEMP(C) RATE(C/M) TIME(MIN) TOTAL TIME(MIN) OPEN CLOSE
1 100 - 100 - 5.0 5.0 SWEEP/SPLIT 1.0 0.0
2 100 - 250 6.0 25.0 30.0 OVERT 17.0 0.0
3 250 - 300 25.0 2.0 32.0

***** SCAN PARAMETERS *****
LOW MASS 43 UP 0.93 L* TOP: 0.00
HIGH MASS 650 DOWN 0.00 L BOTTOM: 0.05
CENT S/P 10 ACTUAL 10 SAMP INT (MS) 0.130 PEAK WIDTH: 1000.
FRAG S/P 10 ACTUAL 10 SAMP INT (MS) 0.130 INTEN/ION: 2
MIN PEAK WIDTH 3 MIN FRAG WIDTH % 80 MIN AREA: 20
A/D THRESHOLD 1 BASELINE 0

***** NODE: CENTROID *****
INTERFACE NUMBER 0
SUB-INTERFACE NUMBER 0
# OF ACQ BUFFERS 16
INSTRUMENT TYPE 0
FULL SCALE MASS 1009 AMU
ZERO SCALE MASS 1 AMU
INTENSITY/ION 2
PEAK WIDTH 1000. MMU
OFFSET AT LOW MASS 0 MMU
OFFSET AT HIGH MASS 0 MMU
VOLTAGE SETTLE TIME(MS) 4

12/30/88 12:06:47
ACQUISITION COMPLETED
SCANS 1 TO 2000 CENTROID

NODE SCANS SECS OUT OF % PEAKS PER SCAN PER SEC
CENTROID 2000 638.0 2000.0 32.9 31365 26. 26.


PARAMETERS
12/30/88 11:10:00 + 33.20 DATA SPPE1230A T1 (SPPE1230A M1)
CALI C12298801 SCANS 1 TO 2000
SAMPLE STD M D POLYETHYLENE SP#041
CONDS +E1/Q1 DM290 100-250 @ 6/MIN PULSE=800
FORMULA - INSTRUMENT Q1 WEIGHT 0.000
SUBMITTED BY MW ANALYST MW ACCT NO P687

ACC VOL 8000 THRESHOLD 1 INTEN/ION 2
A/D S I 0.025 CENT S I 0.130 FRAG S I 0.130
PEAK WIDTH 1000 CENT SAMP/PK 10 FRAG SAMP/PK 10
MIN WIDTH 3 MIN FRAG WIDTH (%) 80 MIN AREA 20

2000 SCANS (433 SECTORS) OF LINEAR UP CENTROID DATA
LOW MASS 43 SCAN TIMES (SECS) UP 0.93 TOP 0.00
HIGH MASS 650 DOWN 0.00 BOTTOM 0.05


```


Diazinon-Pyrogram, 240°/min to 700°C, 4 min (Scanset #672) Figure 28(a)

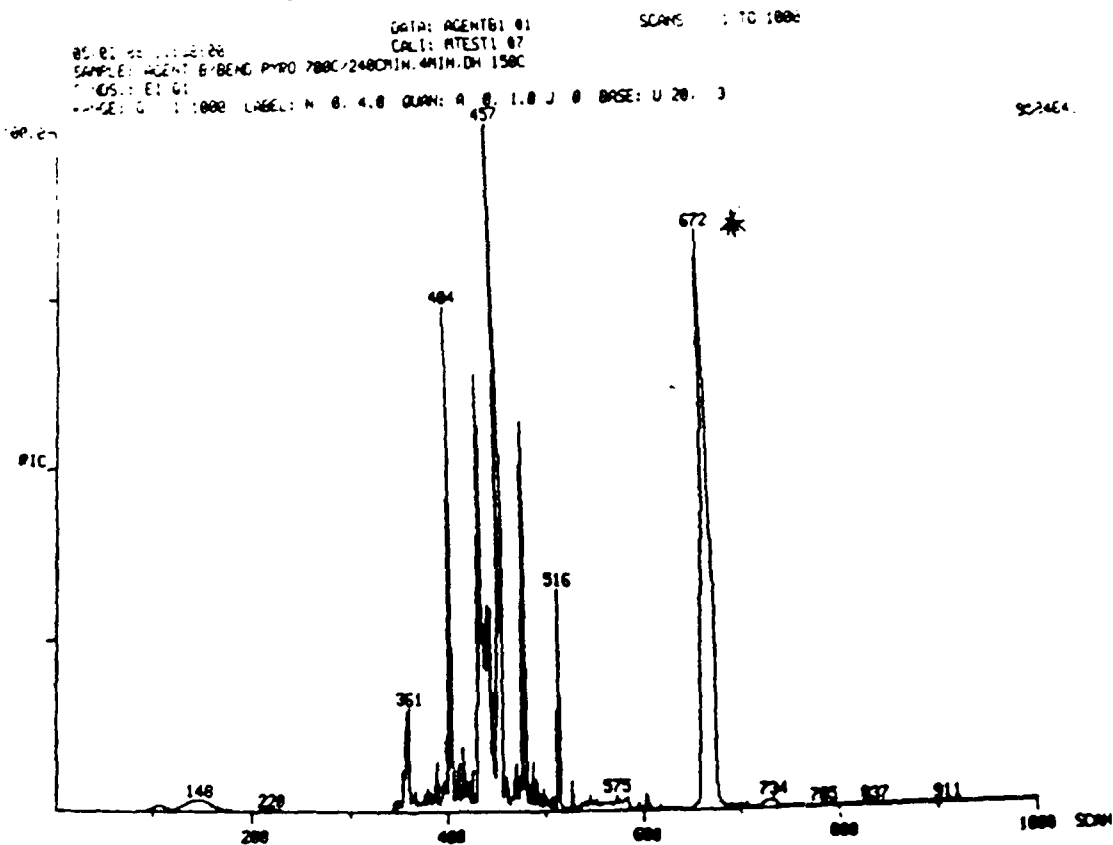

EI-MS of Scanset #672 - Diazinon


Figure 28(b)

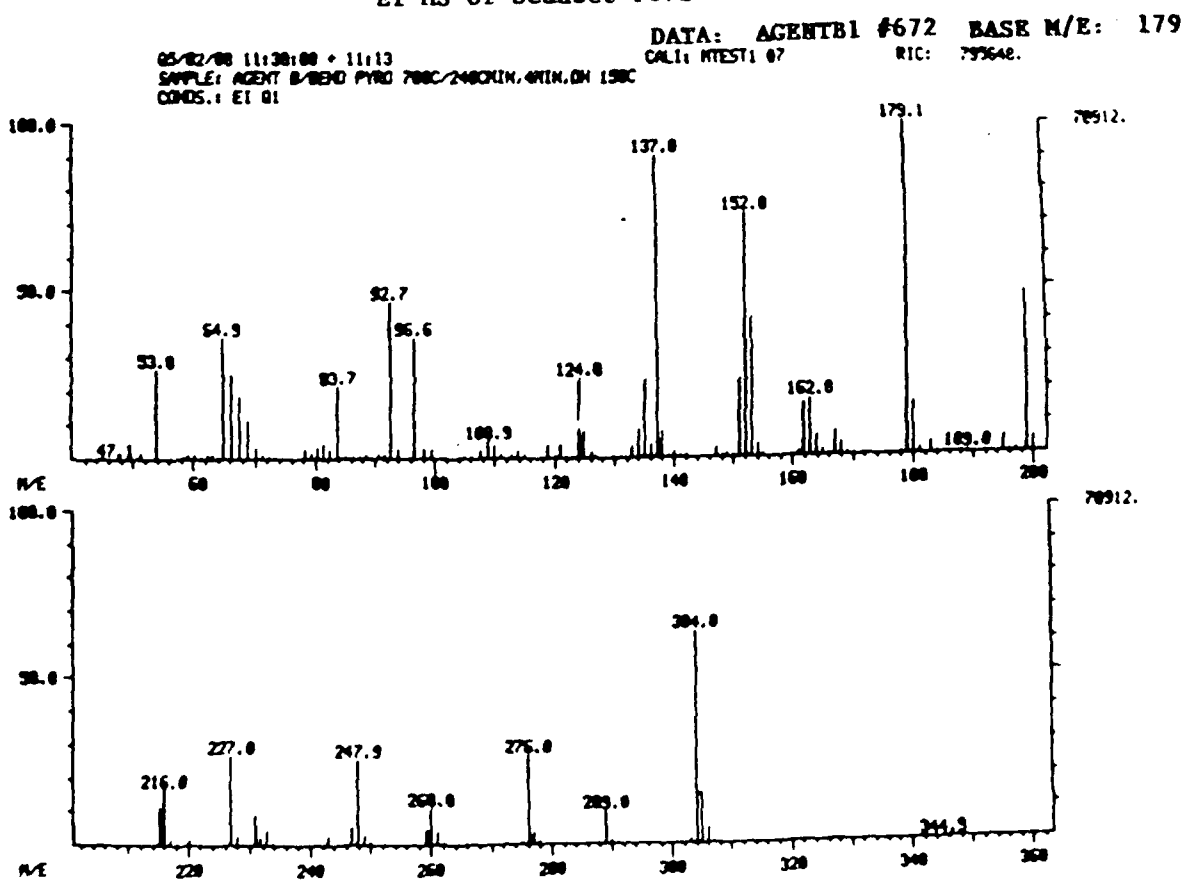

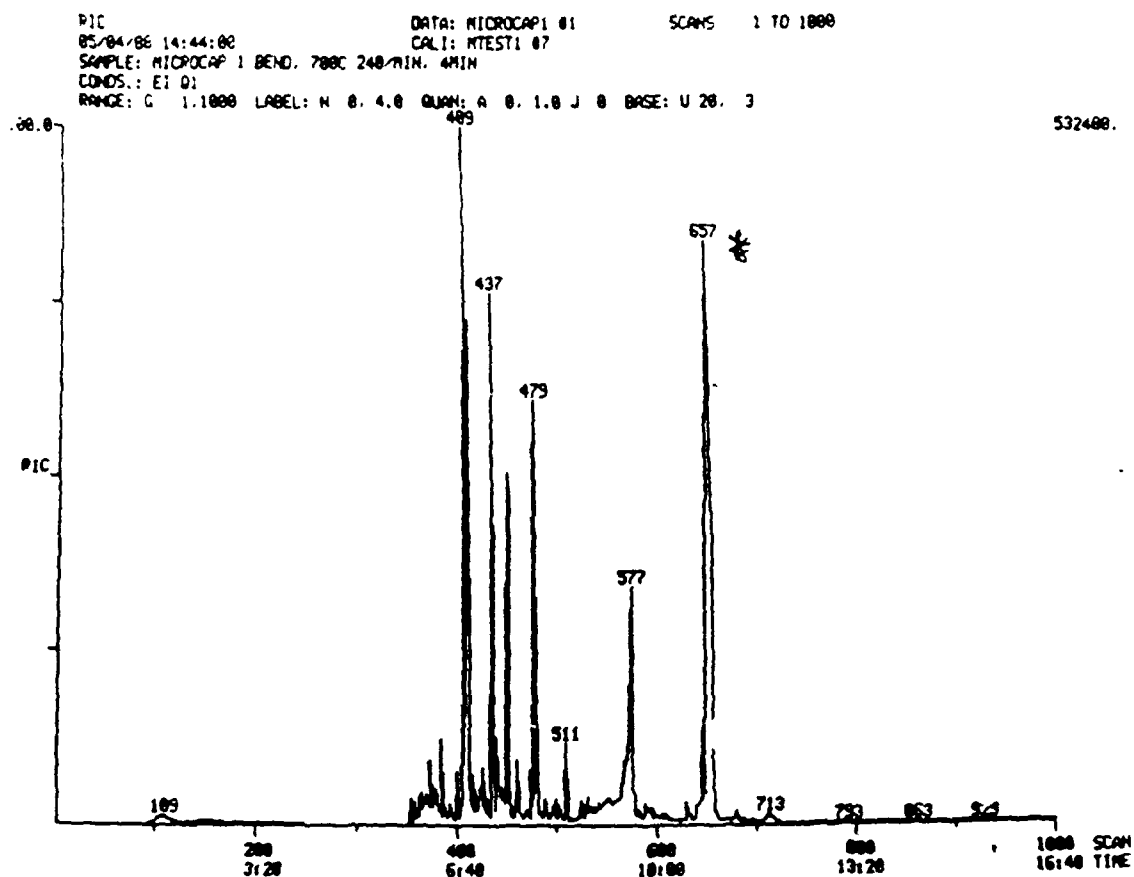

EI-MS of Scanset #656 - Diazinon core


Figure 29(b)

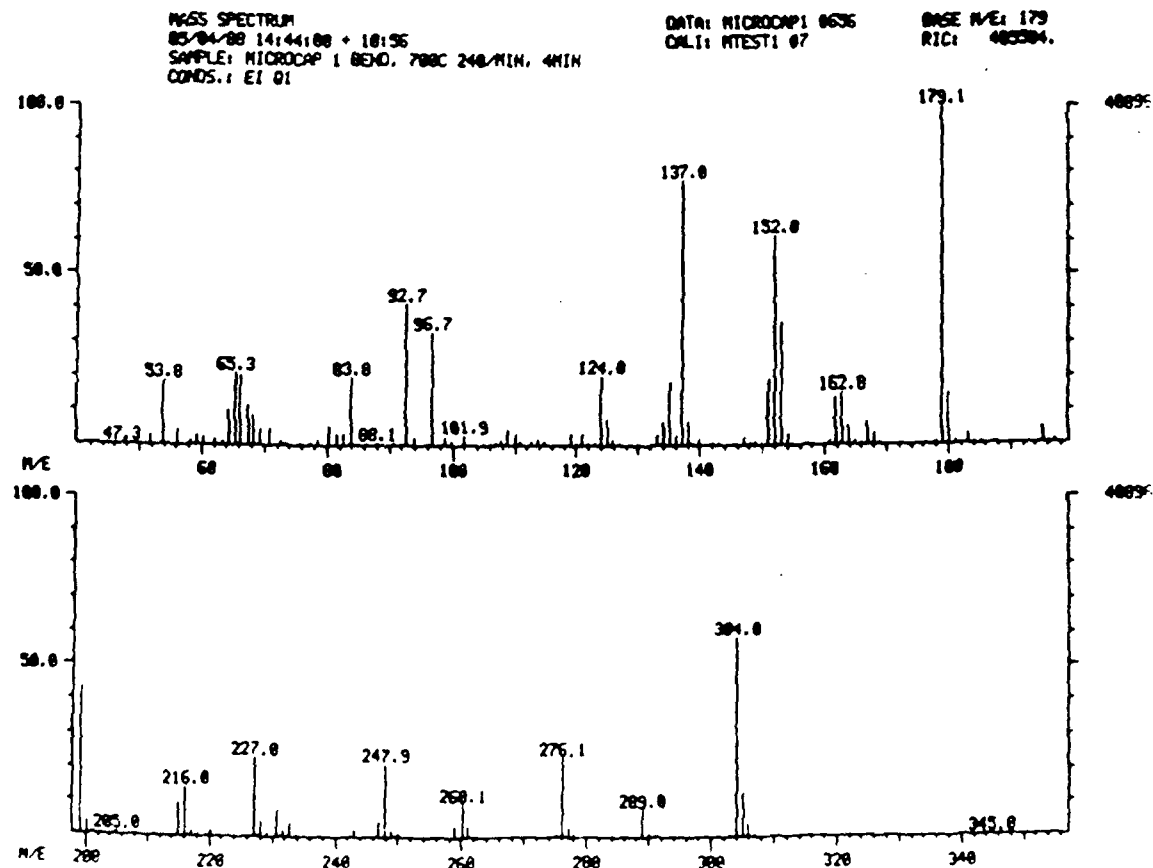

Microcap #1 RIC 180°C/min to 600°C (Scanset #662) Figure 30(a)

RIC
 05/04/88 15:13:02
 SAMPLE: MICROCAP
 COND.: EI 01
 RANGE: C 1.1000 LABEL: N 0. 4.0 QUAN: A 0. 1.0 J 0 BASE: U 20. 3
 DATA: MICROCAPIR 01
 CALI: MTEST1 07
 100/MIN. 0MIN
 SCANS 1 TO 1000

783368.

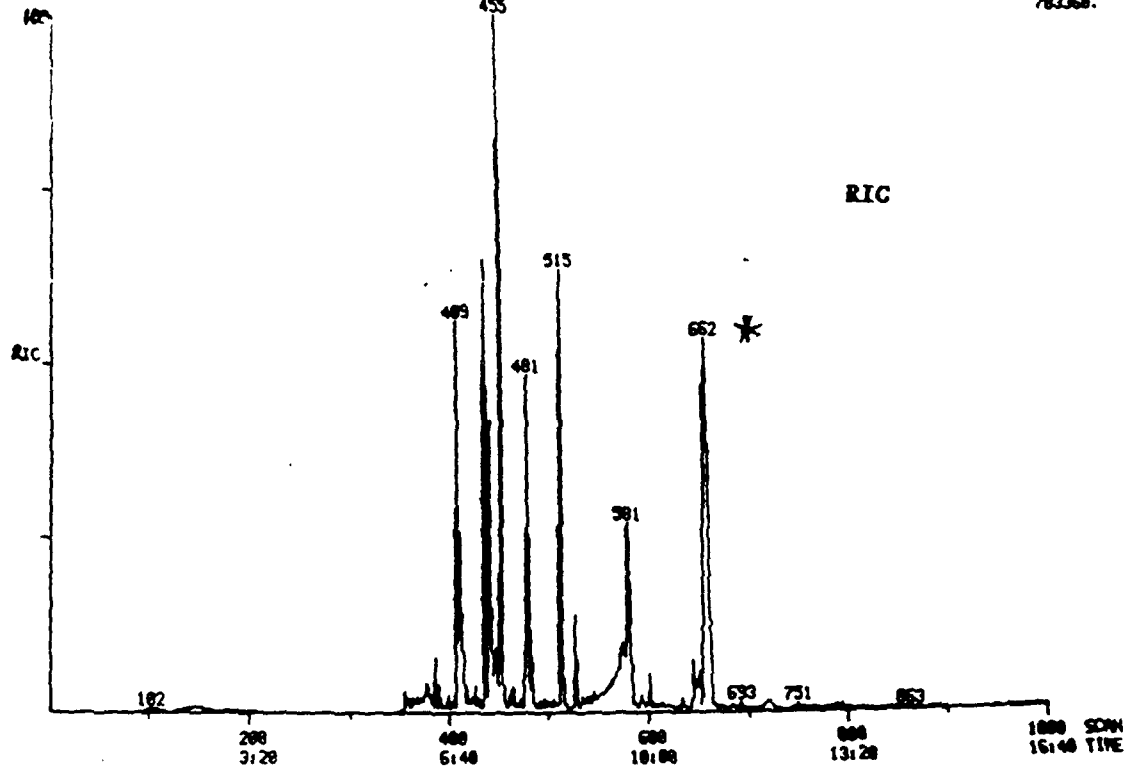


EI-MS of Scanset #662 - Diazinon core

Figure 30(b)


MASS SPECTRUM
 05/04/88 15:13:00 + 11:02
 SAMPLE: MICROCAP
 COND.: EI 01

DATA: MICROCAPIR 0662
 CALI: MTEST1 07
 100/MIN. 0MIN
 BASE M/E: 179
 RIC: 418384.

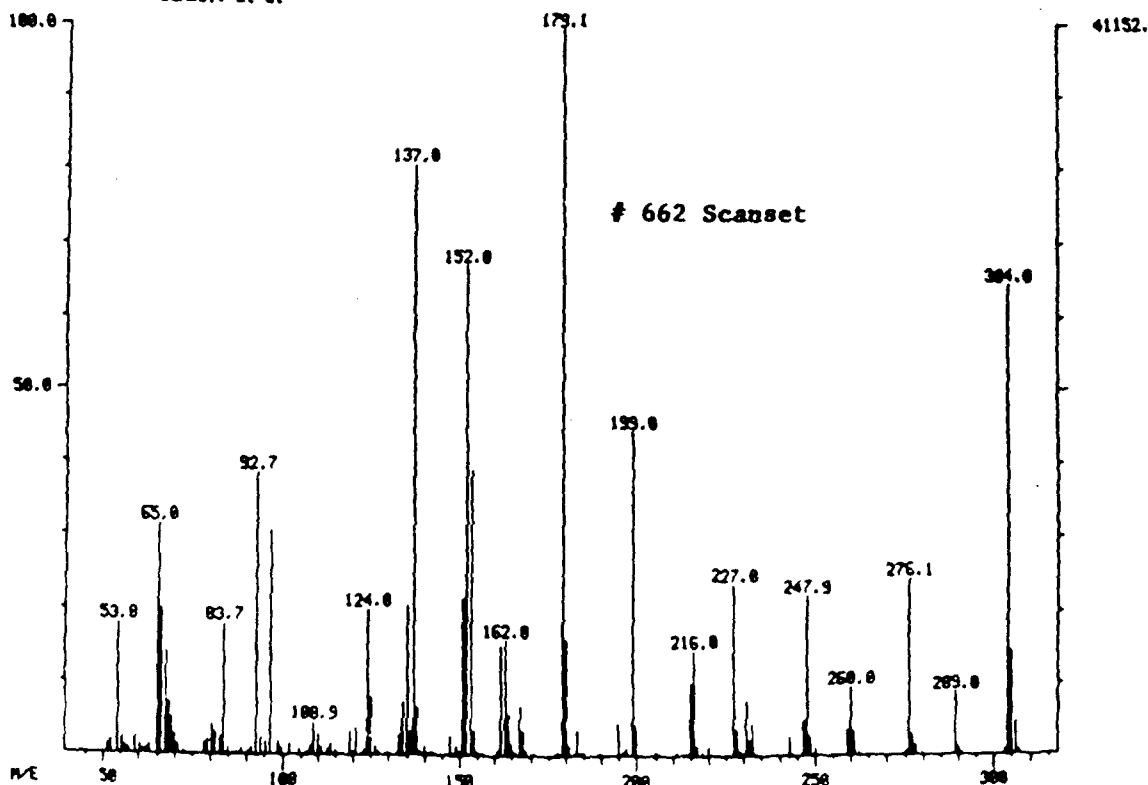


Microcap #1 RIC 240°/min to 700°C (Scanset #668)

Figure 30(c)

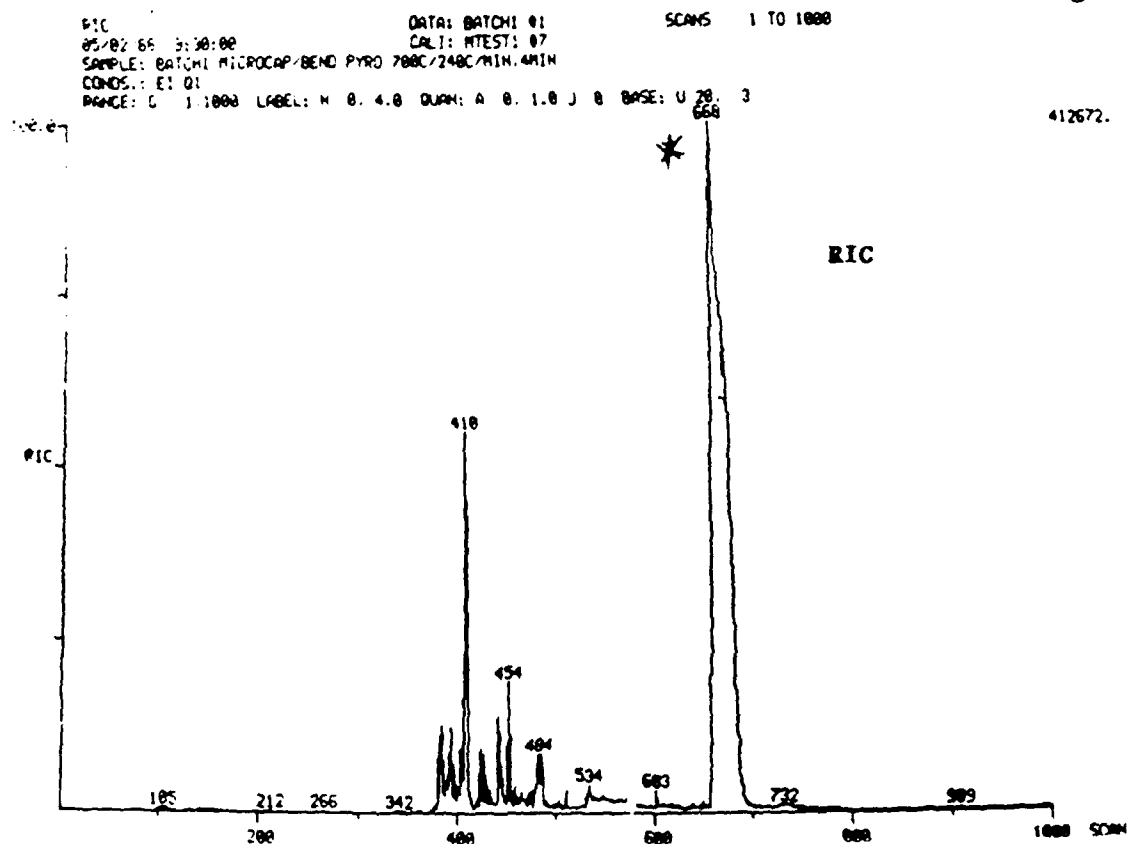

Microcap #1 RIC Pulse 500°/20 sec (Scanset #753)

Figure 31(a)

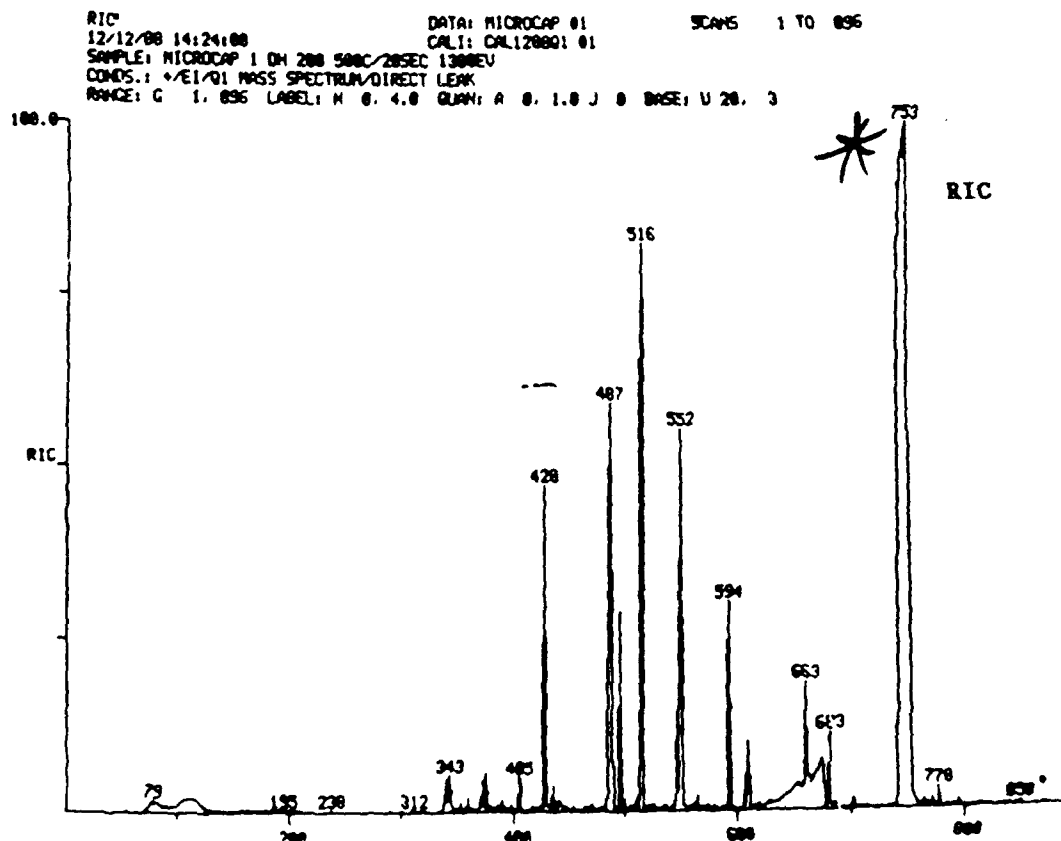


EI-MS of Scanset #753 - Diazinon Core

Figure 31(b)

MASS SPECTRUM
12/12/88 14:24:00 + 13:11
SAMPLE: MICROCAP 1 DN 200 900C/20SEC 1300EV
CONDOS.: +E1/01 MASS SPECTRUM/DIRECT LEAK

DATA: MICROCAP #753
CALI: CAL120801 01


BASE M/E: 179
RIC: 408384.

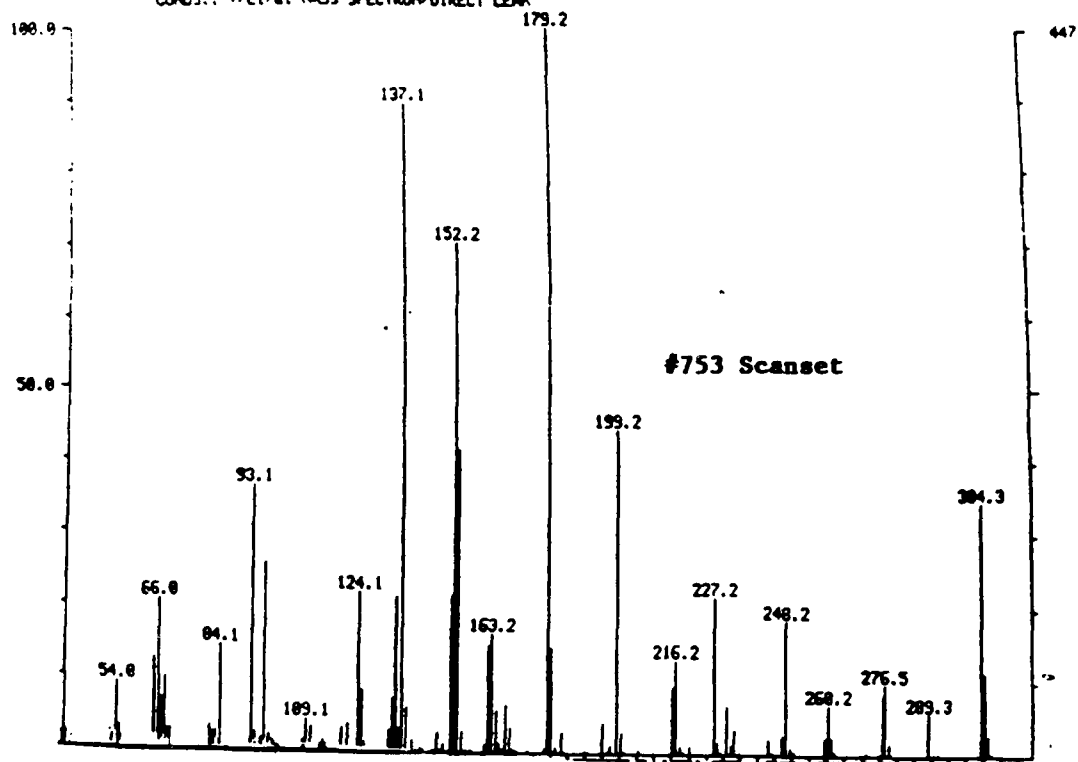

MS-MS Daughter Ions from Parent Ion of Microcap #1 Sample Figure 31(c)


SAMPLE: MICROCAP 1 DN200C, 900C/20SEC, 1000EV
CONDOS.: +E1/CAO — DAUGHTERS OF 304 /CE1-18.6/DIRECT LEAK

RIC: 7328.

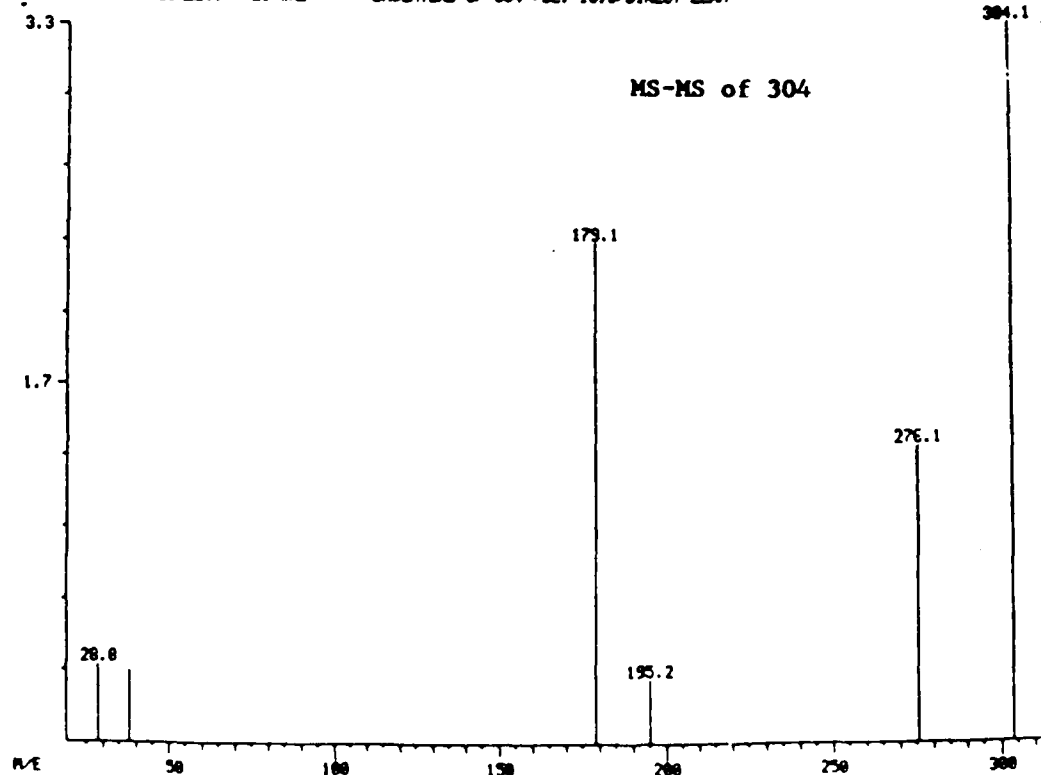


Microcap #4 - Cypermethrin Core
Pulse 500°/20 sec (Scanset #608)

Figure 32(a)

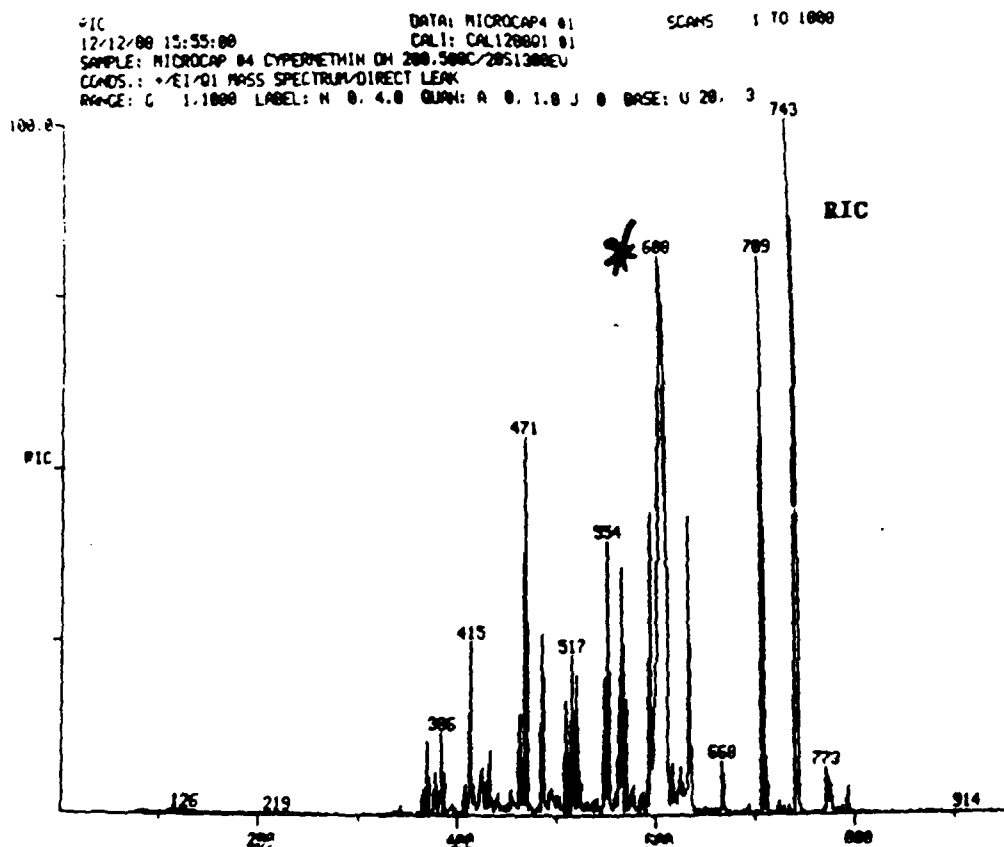


EI-MS of Scanset #608 - Cypermethrin Core

Figure 32(b)

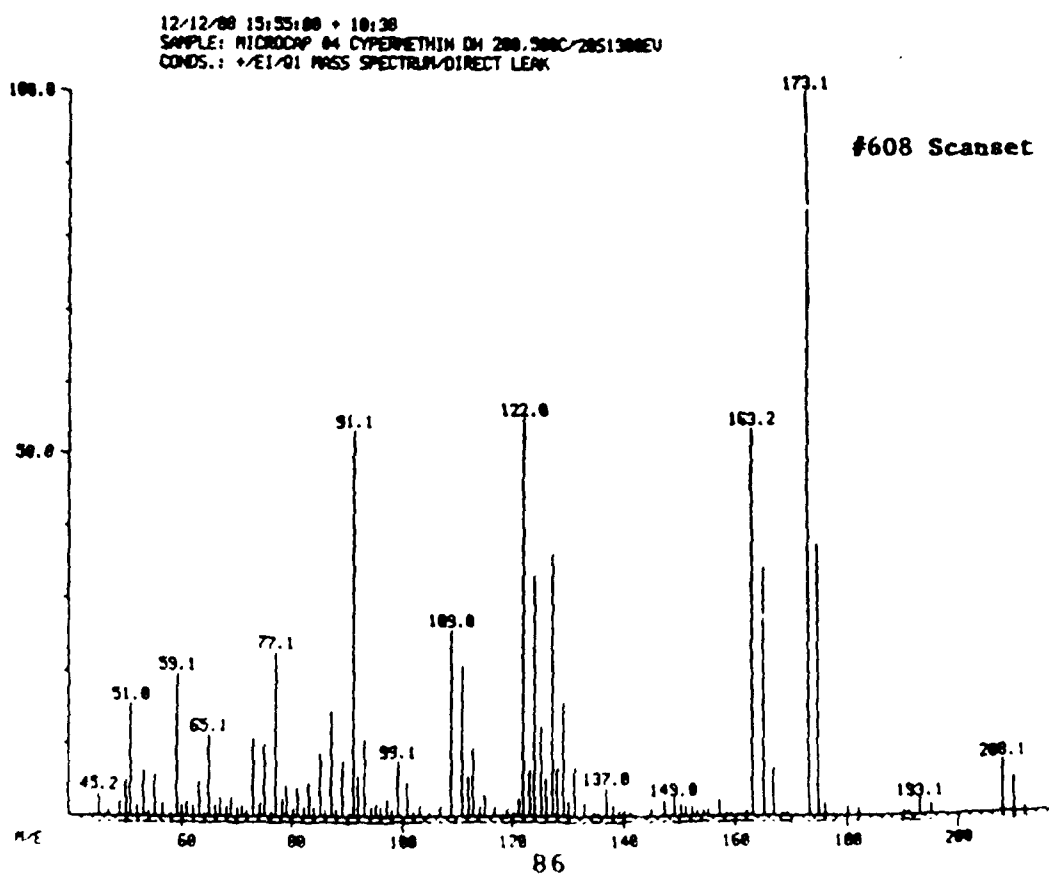


Microcap #4 (RIC) Cypermethrin Core (Scanset #395) Figure 33(a)

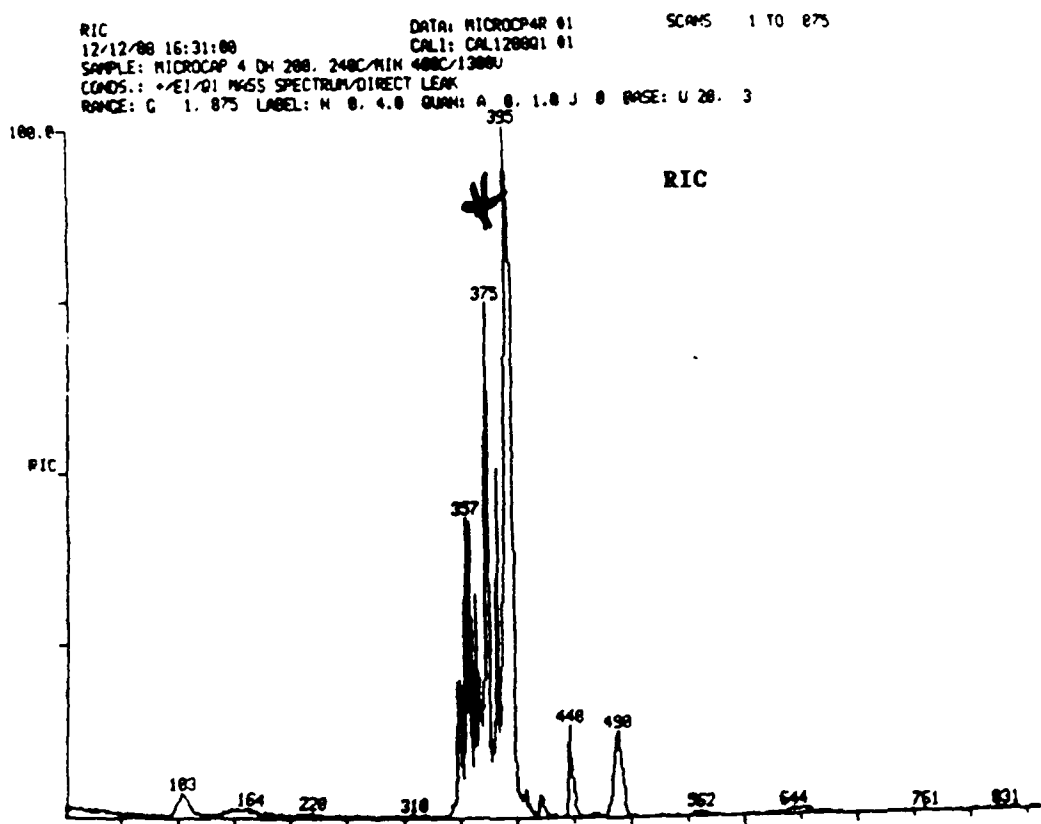


EI-MS of Scanset #395 from Microcap #4 RIC

Figure 33(b)

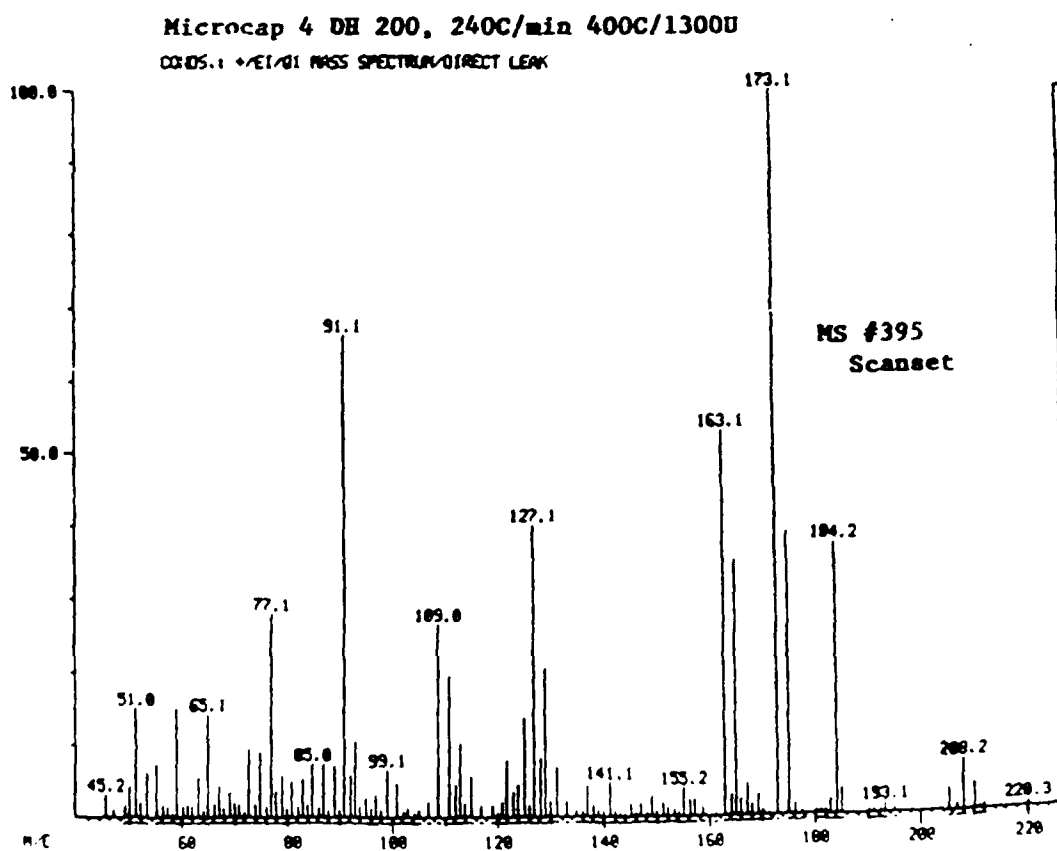


Microcap #5 (RIC) Dursban Core (Scanset #831)

Figure 34(a)

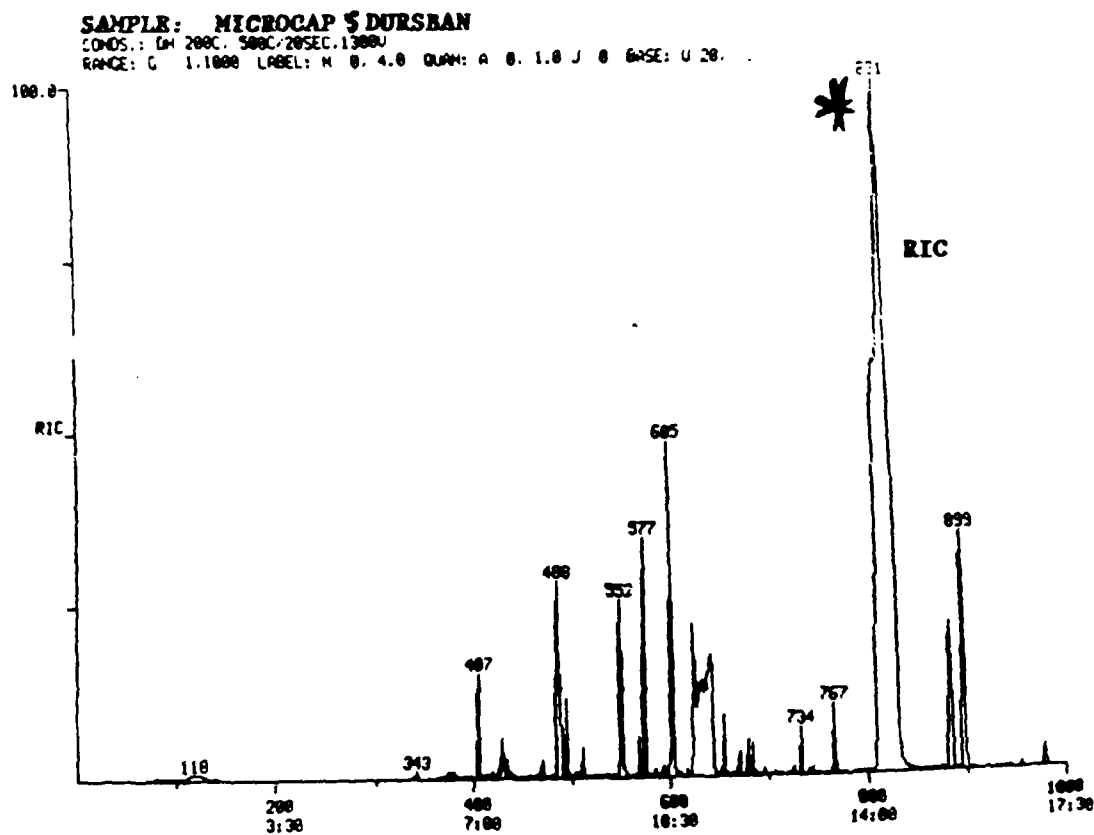


EI-MS of Scanset #831 Dursban Core

Figure 34(b)

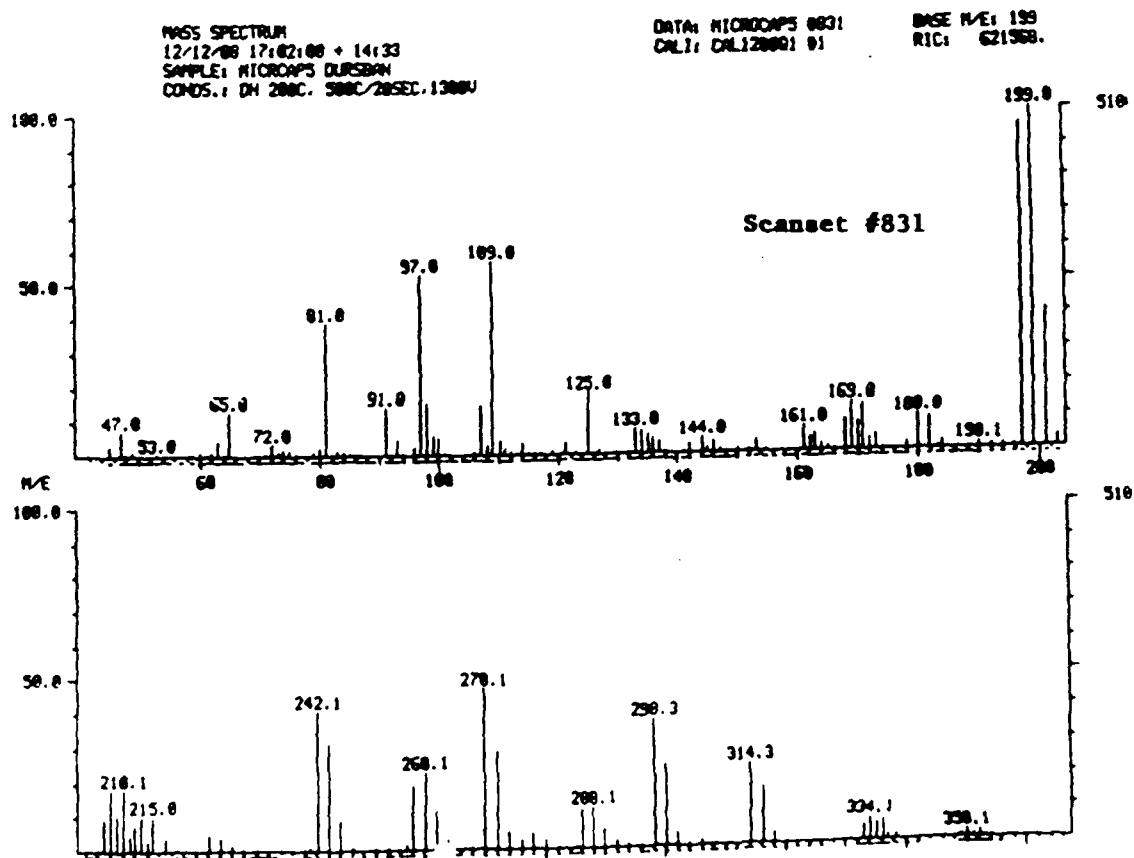


MS-MS Daughter Ions from Parent 199 Ion
of Microcap #5 Sample

Figure 34(c)

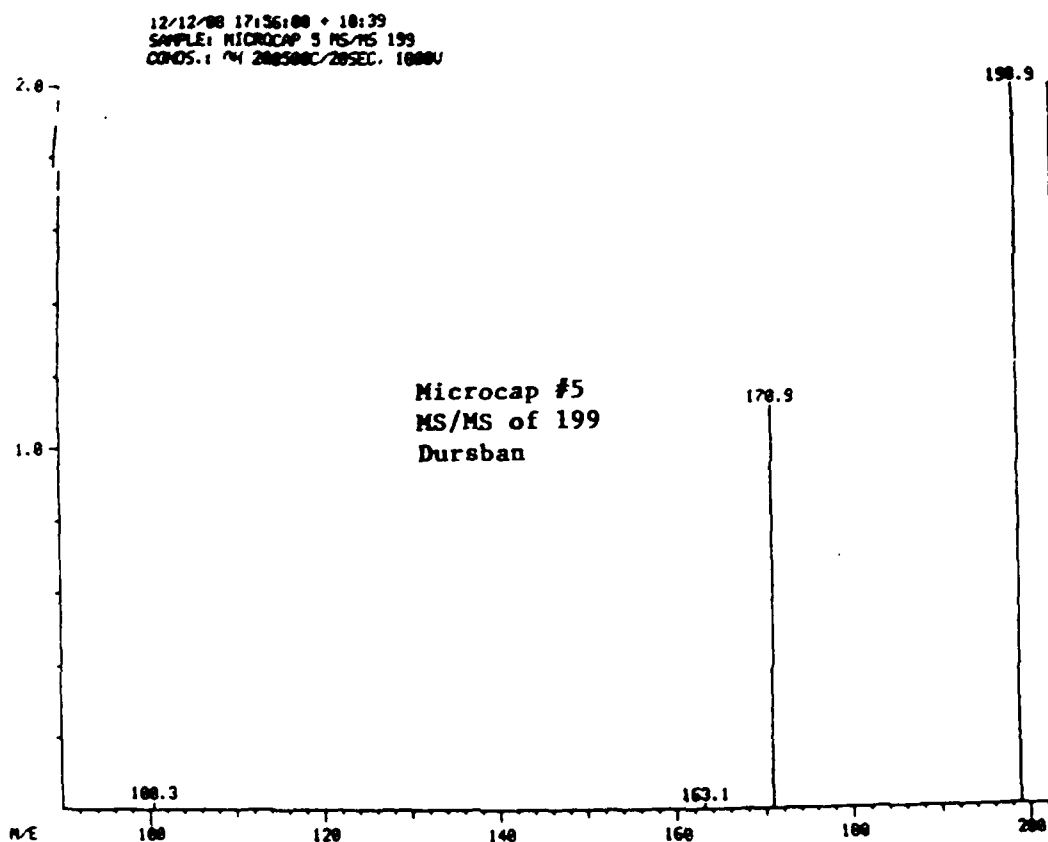


Pyrogram - RIC of Dodge Gelatin in PS-PES Shell

Figure 35(a)

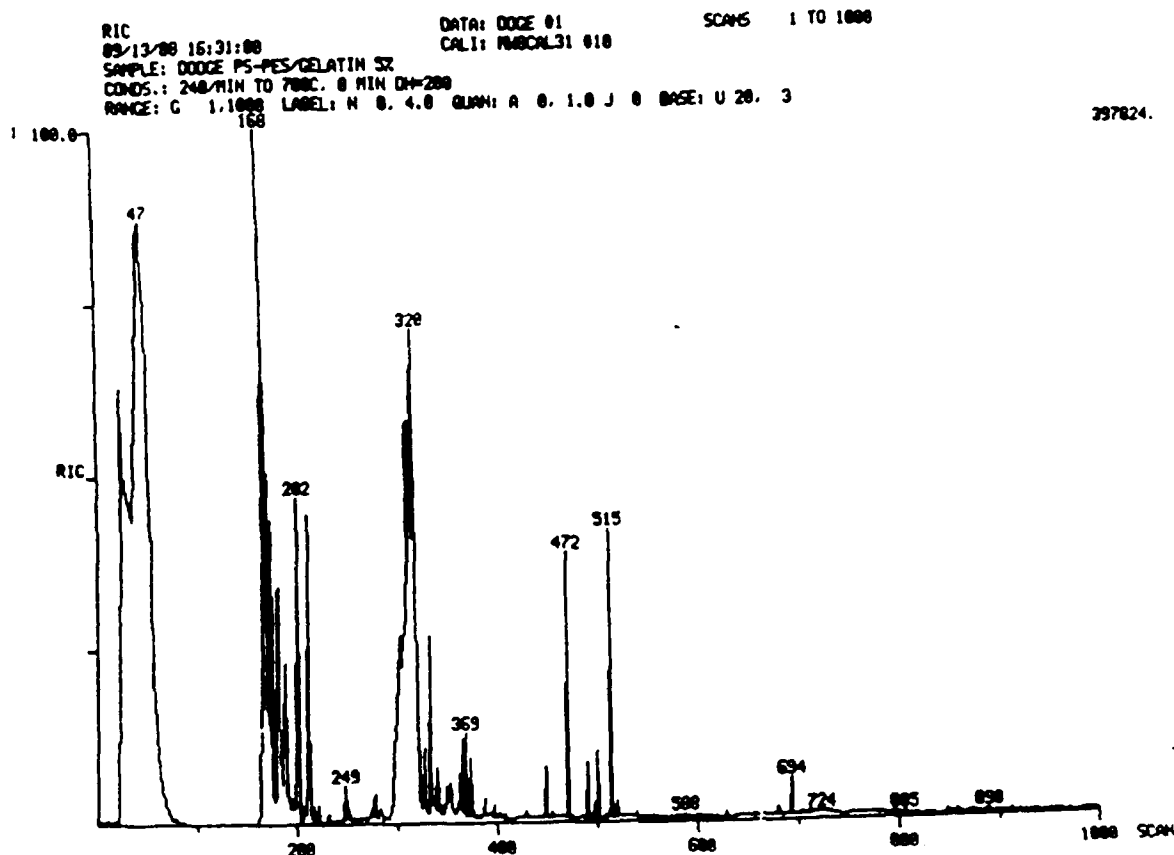


Pyrogram - RIC of Emulsifier Gelatin (Bend)

Figure 35(b)

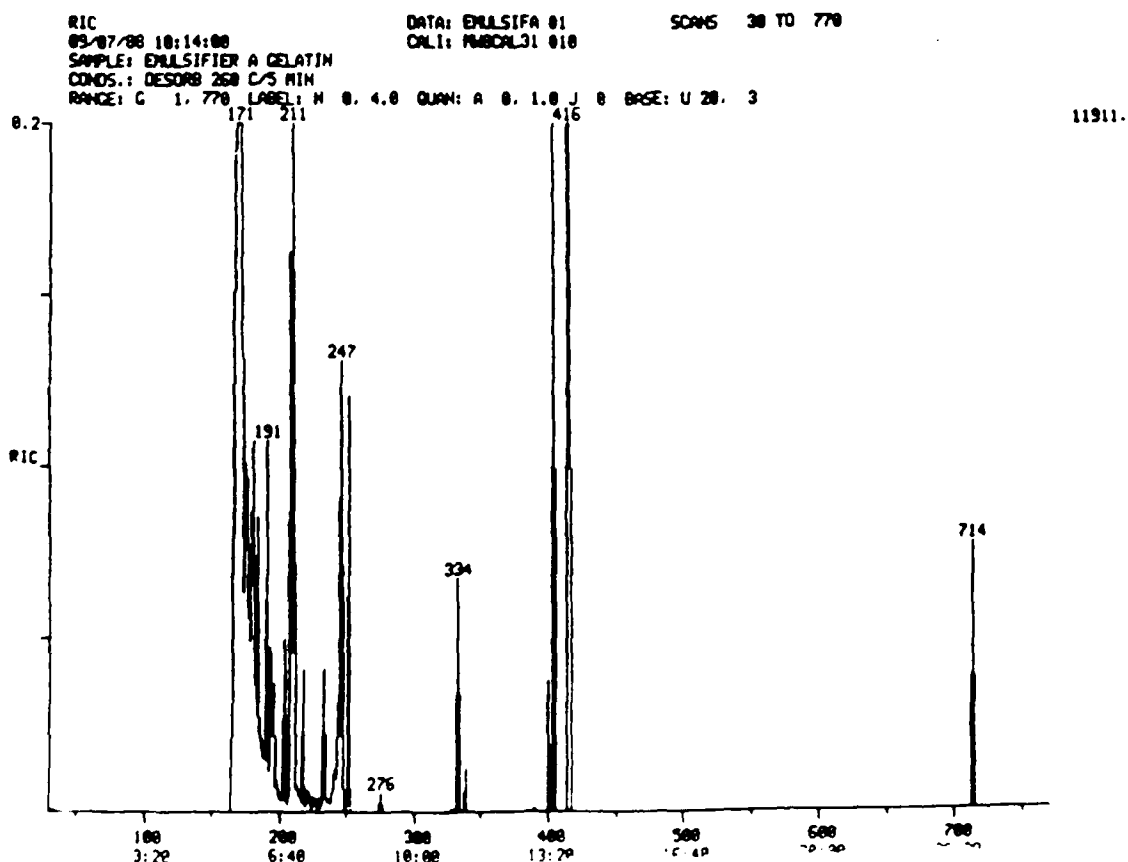


Pyrogram - RIC of Gelatin (P. Snyder)


Figure 35(c)

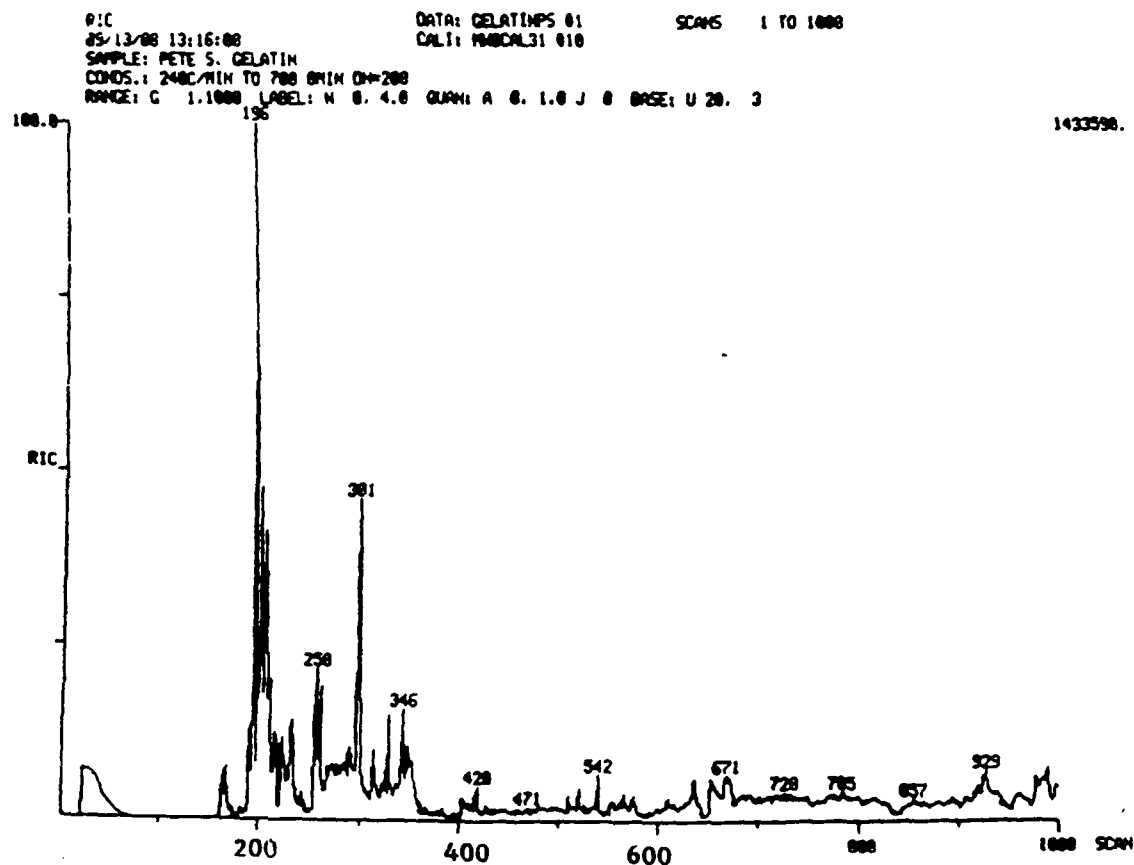


Pyrogram - RIC of Gelatin (~500 μ Particles - Dodge) Figure 35(d)

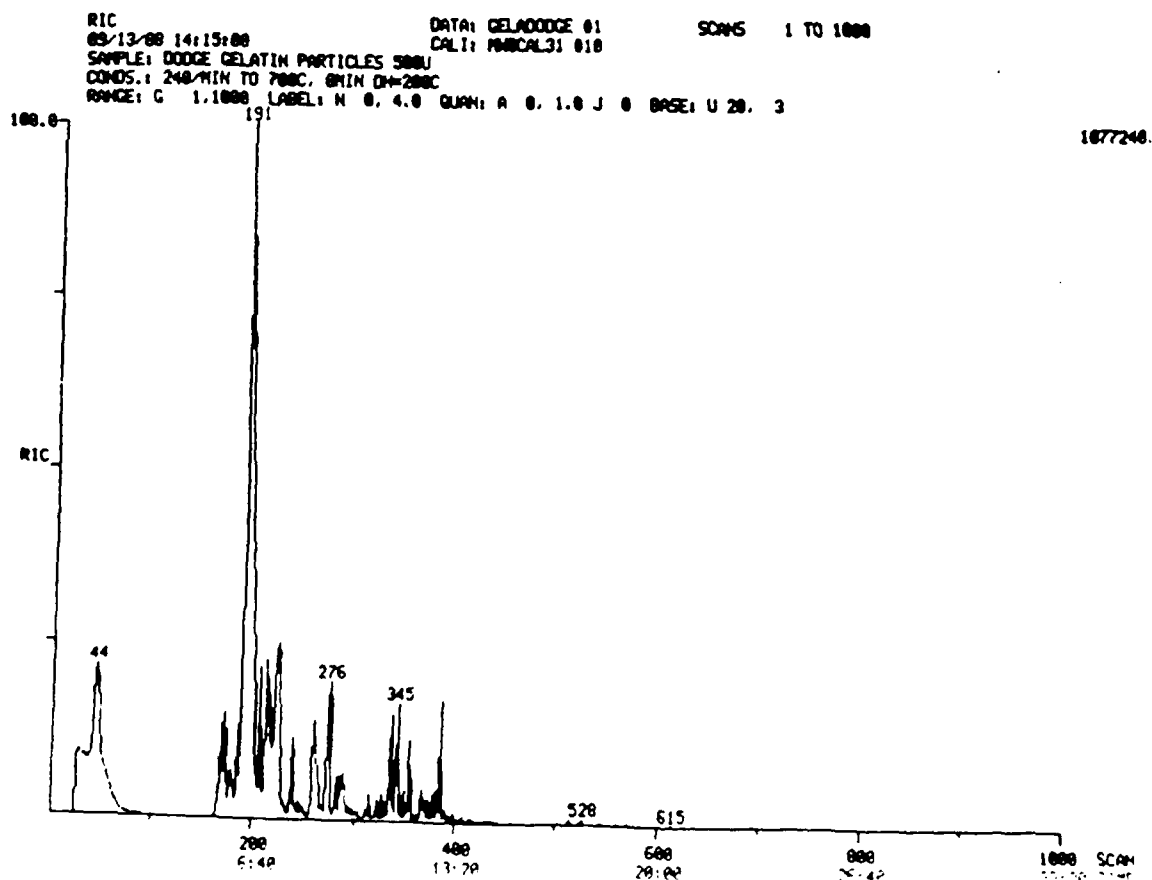

Pyrogram - RIC of Gelatin (~500 μ Particles - Dodge) Figure 35(e)
Desorb 260°C/5 min (Scanset #206)

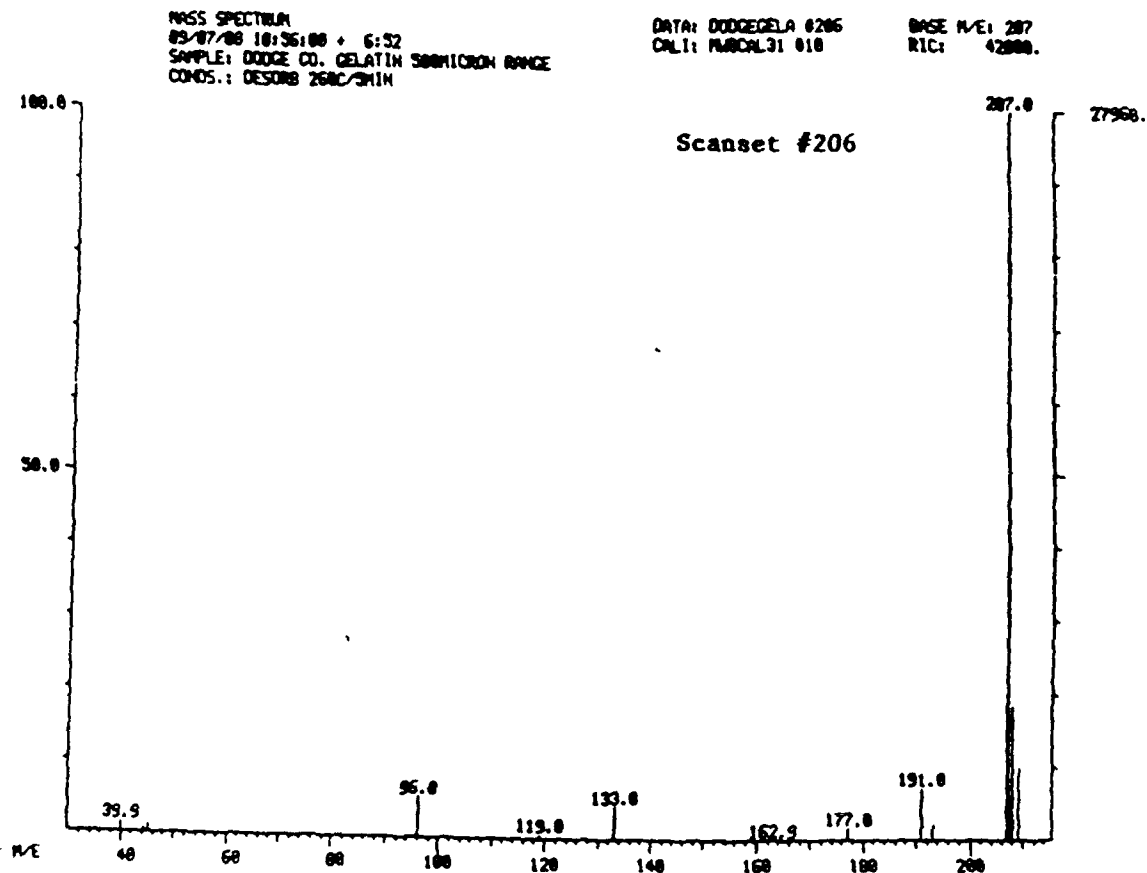


EI-MS of Scanset #206 - Gelatin Particles Figure 35(f)

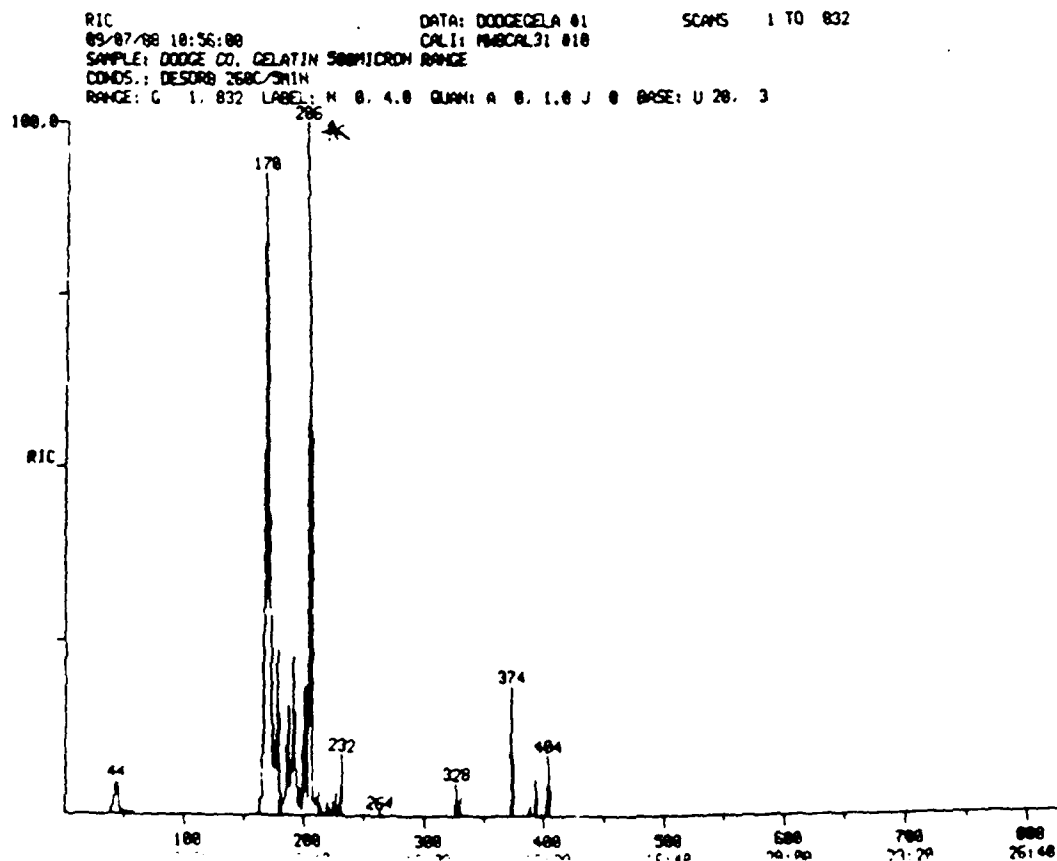

PYROLYSIS GC-FTIR


Total Response Chromatogram (TRC)
Microcap #1 (Diazinon Core)

Figure 36(a)

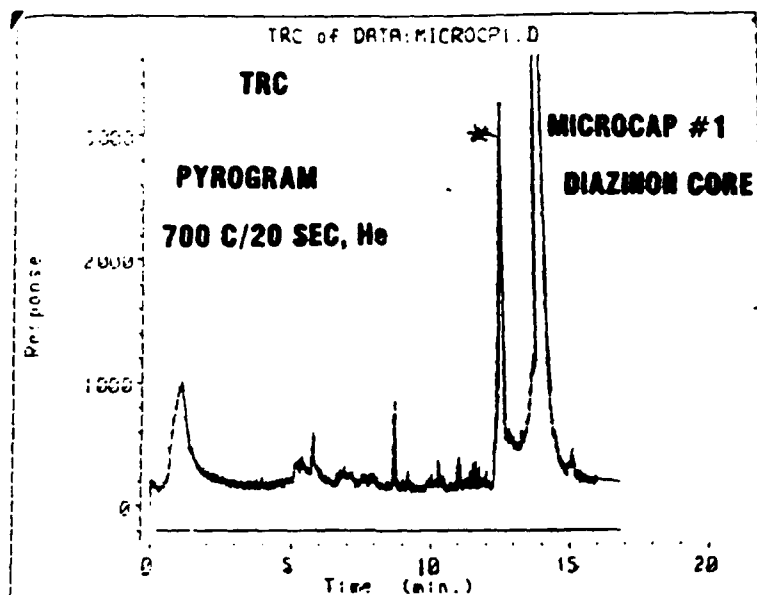

FTIR of ~12.5 min Peak with Library Match
to Phosphorothioic Derivatives

Figure 36(b)

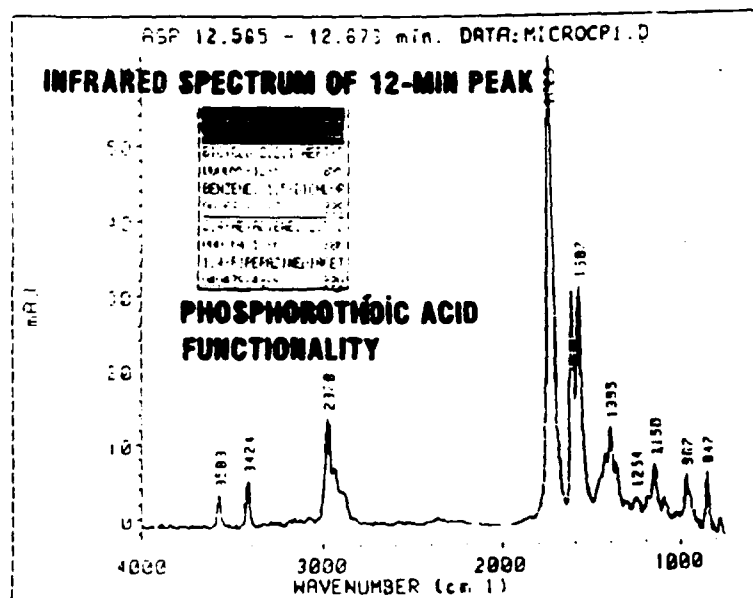

SPE-SFC of Propellant Formulations Figure 37

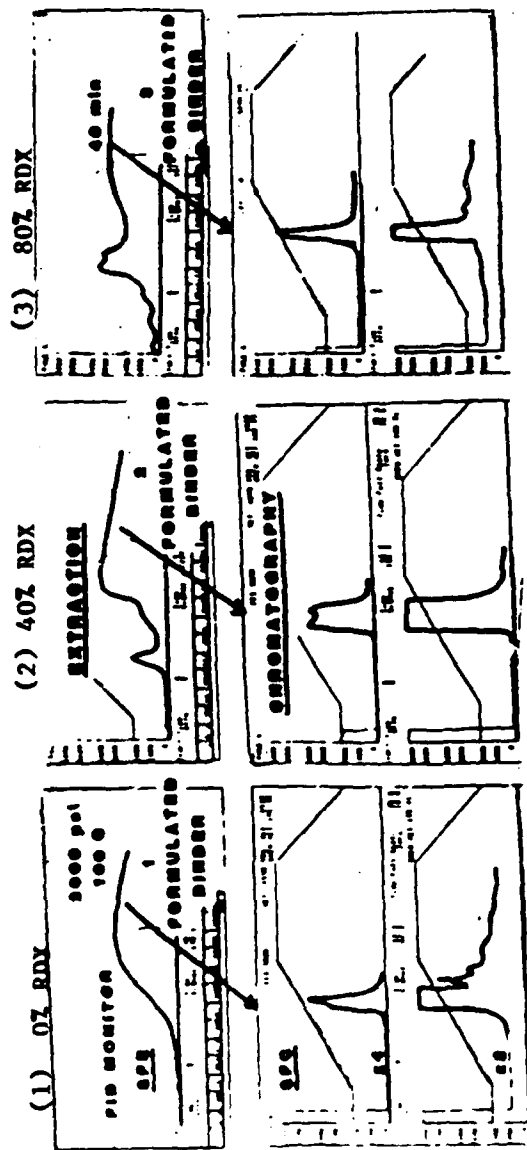


Figure 39(a)

SPE METHOD DEVELOPMENT
BOX-BENKEN EXPERIMENTAL DESIGN

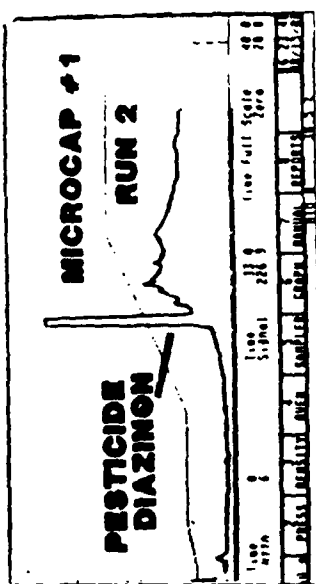


Figure 38
Two Samples from Same Batch

SPE-SFC of Microcap #1,

SPE-SFC

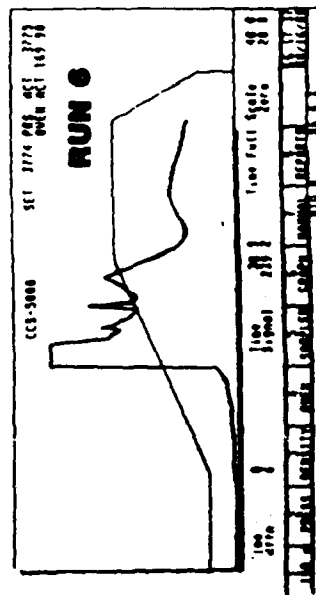
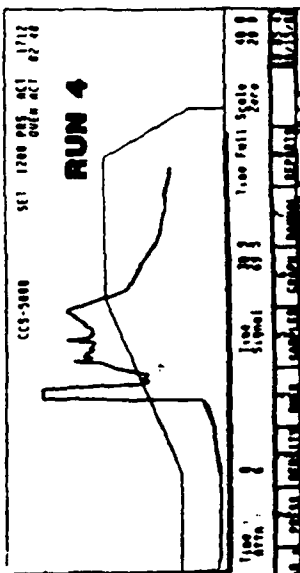
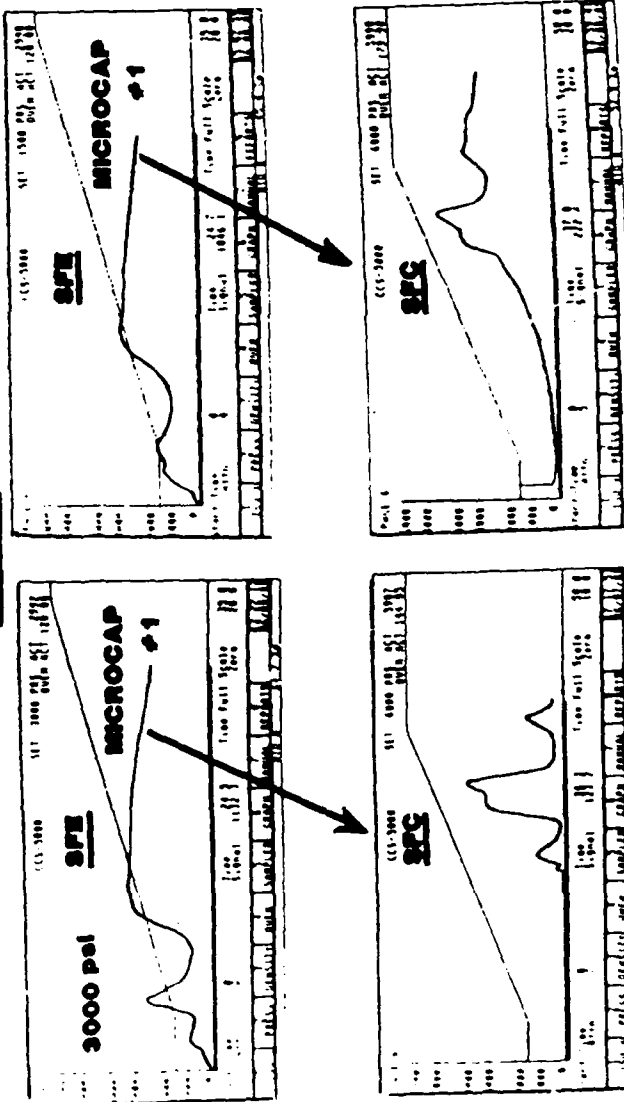


Figure 39(b)

SFE PROFILES FOR MICROENCAPSULATED PESTICIDES

BEND RESEARCH, INC.
BEND, OREGON

SFE

1200 PSI, 8 MIN

100°C, 200 PSI/MIN

TO 4000 PSI



(1) POLYMER SHELL A

DIAZINON CORE

FID

MONITOR




(2) POLYMER SHELL B

DIAZINON CORE




(3) POLYMER SHELL C

DIAZINON CORE

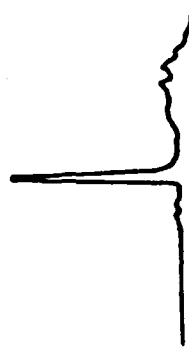

(4) POLYMER SHELL B

CYPERMETHRIN

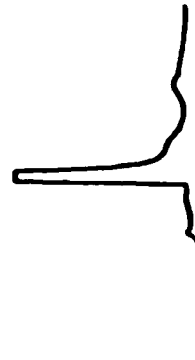

(5) POLYMER SHELL A

DURSBAN

Figure 39(c)

SFE PROFILES FOR MICROCAP #1

BOX-BEHNKEN EXPERIMENTAL DESIGN

FID

MONITOR



Figure 39(d)

SFE PROFILES MICROCAP #2


BOX-BEHNKEN EXPERIMENTAL DESIGN


FID
MONITOR

SFE PROFILES FOR MICROCAP #3

BOX-BEHNKEN EXPERIMENTAL DESIGN

FID
MONITOR

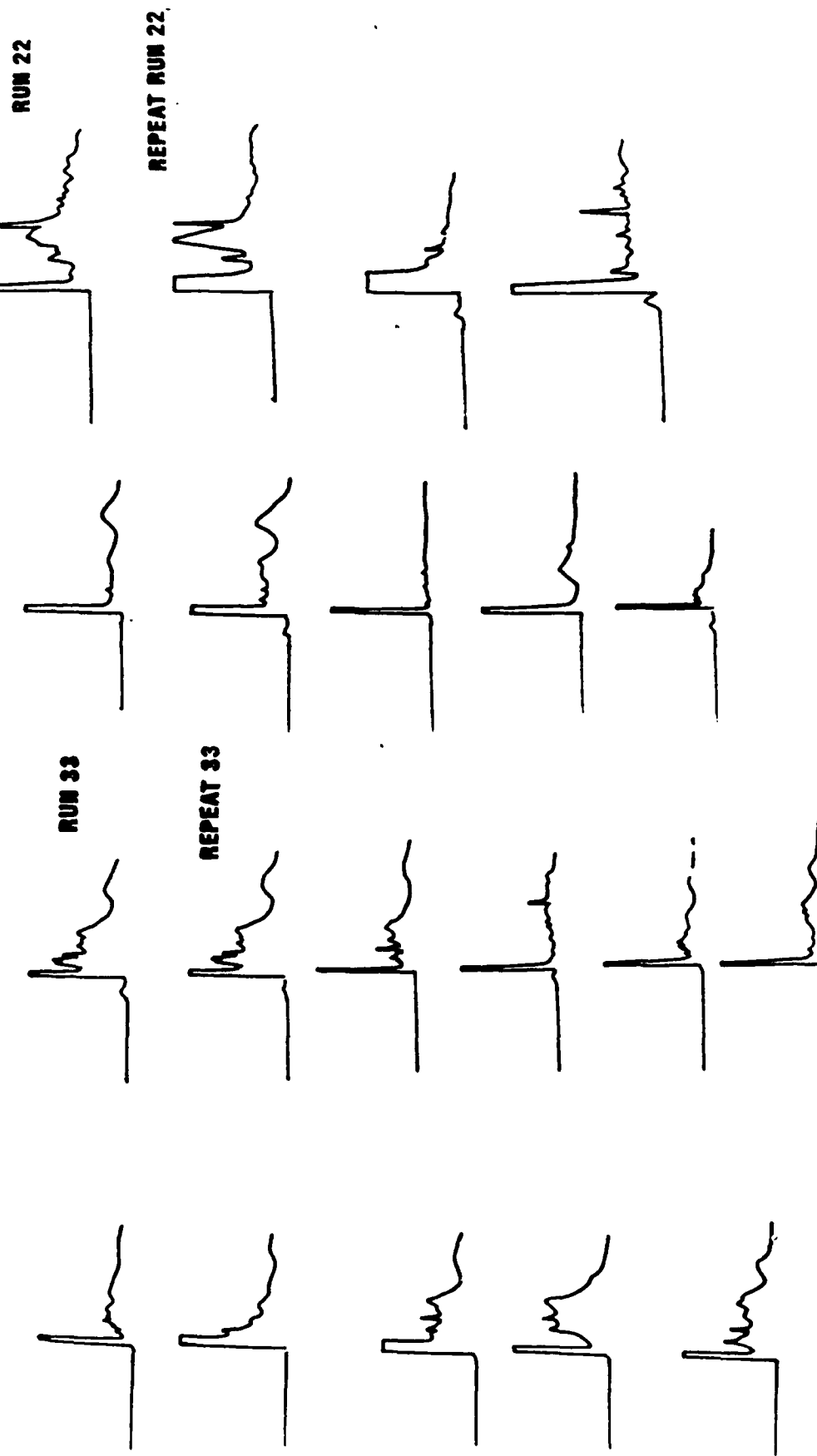


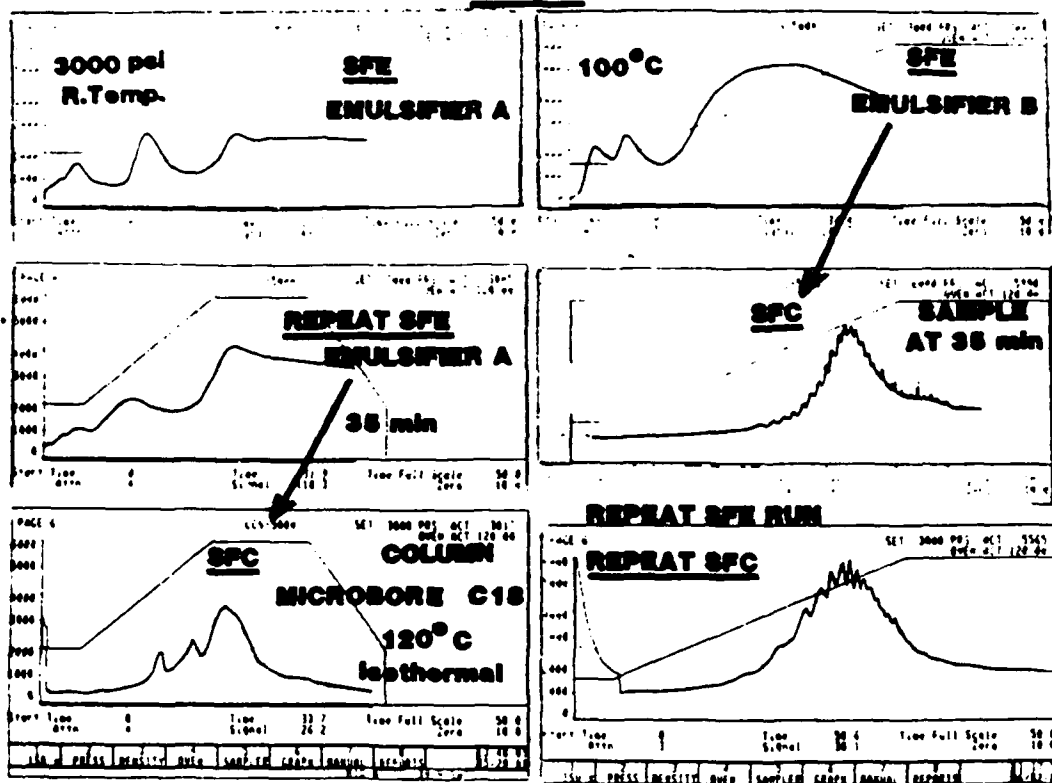SUPERCritical FLUID EXTRACTION-CHROMATOGRAPHY**SFE-SFC**

SFR-SFC of RDX and RDX and Catalyst

Figure 41


Revised Chromatographic Sciences Vol 27 March 1999

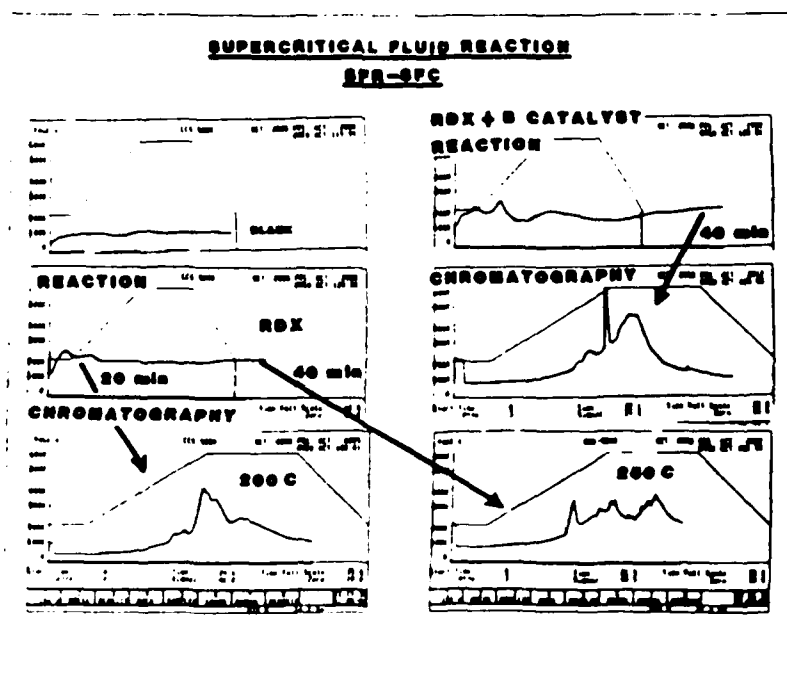


SFD (SFE) - SFC of Simulant in Soil

Figure 42(a)

SUPERCritical FLUID EXTRACTION-CHROMATOGRAPHY

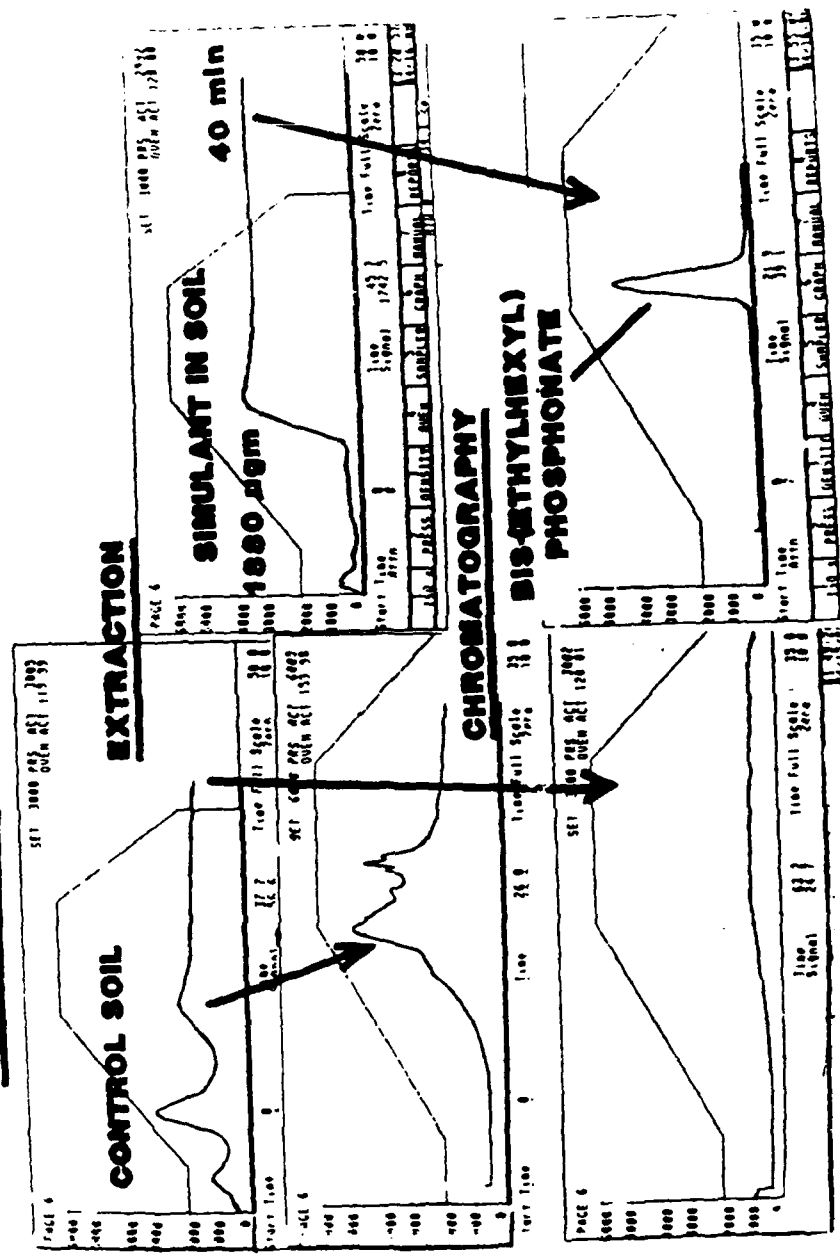




SFD-SFC and Capillary GC of Charcoals from Air Sampling Figure 42(b)

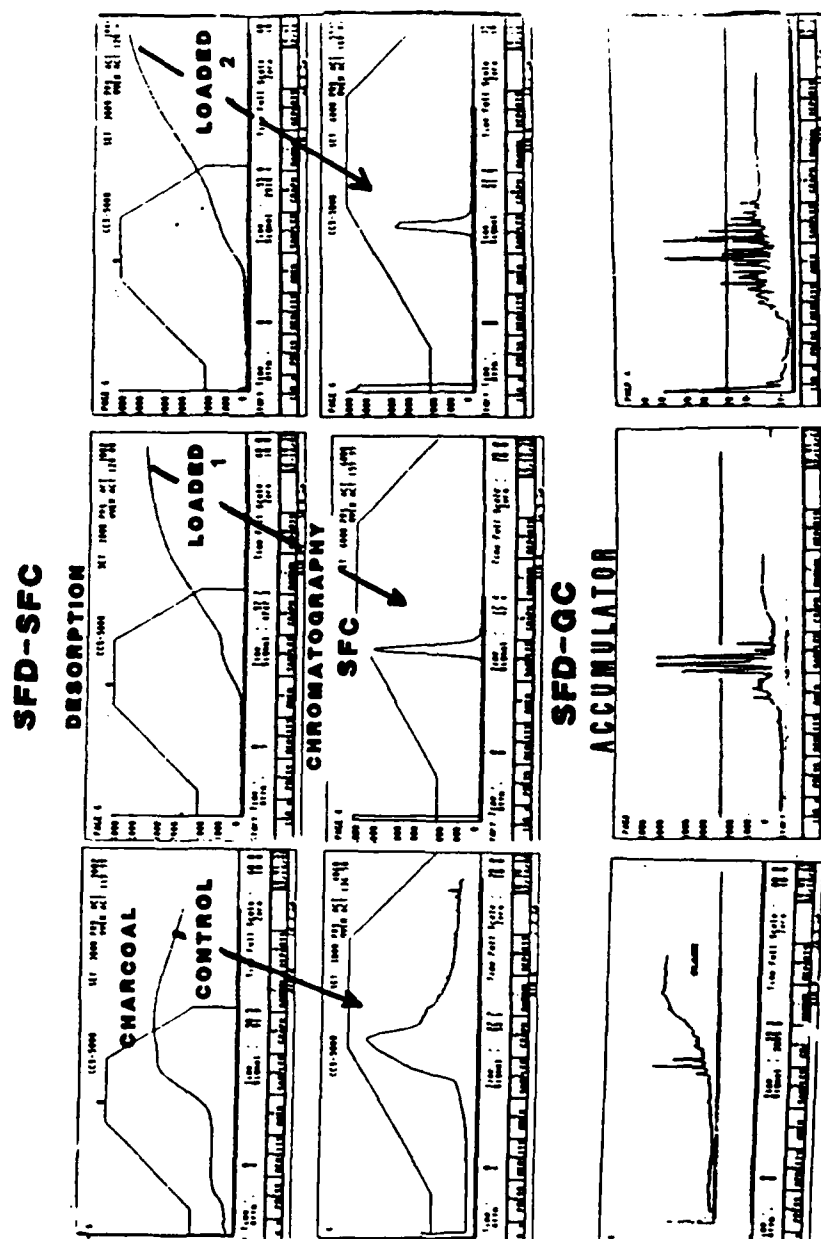


CONCLUSIONS

- PESTICIDE CORES (DIAZINON, CYPERMETHRIN, DURSABAN) ARE READILY EXTRACTED FROM FORMULATED MICROCAPSULES HAVING THREE DIFFERENT POLYMER SHELLS OVER A SELECTED SFE RANGE OF PRESSURE (DENSITY), TEMPERATURE, AND TIME.
- SURFACTANT/EMULSIFIER ADDITIVES ARE MORE SENSITIVE TO SFE CONDITIONS. QUALITATIVE AND QUANTITATIVE ANALYSES FOR FULLY FORMULATED PRODUCTS (SHELL, CORE, ADDITIVES) REQUIRE FURTHER METHOD DEVELOPMENT USING SFE/ACCUMULATOR-SFC AND CAPILLARY GC AUTOMATED/INTERFACED INSTRUMENTATION.
- BOX-BEHNKEN EXPERIMENTAL DESIGN HAS BEEN EXTENDED TO INCLUDE A MULTIVARIATE RESPONSE (FACTOR ANALYSIS) FOR COMPLEX CHROMATOGRAPHIC/SPECTRAL PATTERN EVALUATIONS.
- SFE PROFILES PROVIDE IMPORTANT COMPOSITIONAL INFORMATION FOR PROCESS AND QC/QA OPERATIONS USED IN THE AGRICULTURAL, FOOD, BIOMEDICAL, AND CONSUMER PRODUCT INDUSTRIES IN DEVELOPMENT OF CONTROLLED-RELEASE MATERIALS.

REMOTE SPECTROSCOPY Figure 43
FIBER OPTIC MONITOR (FOM) PROTOTYPES
FOR CCS MODEL 5000 SFC/GC
SUPERCRITICAL FLUID EXTRACTOR/ACCUMULATOR MODULE

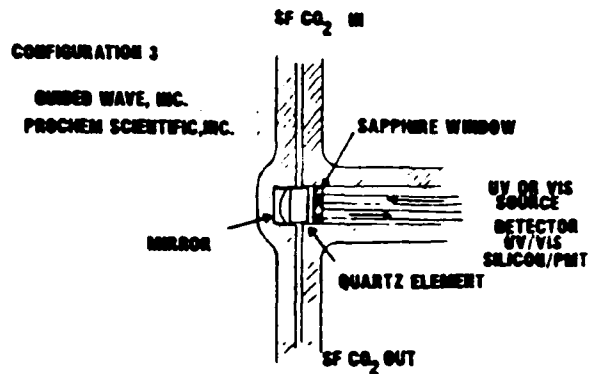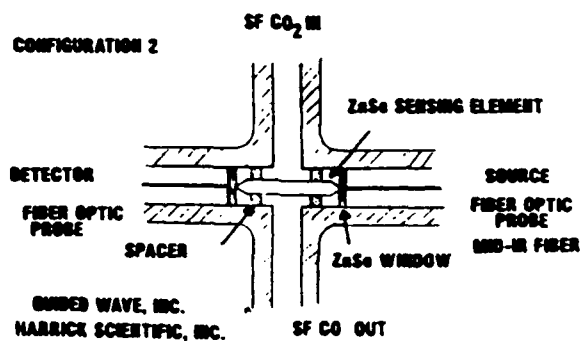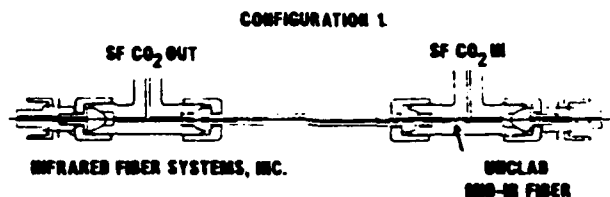


TRACE AROMATICS - CONFIGURATION 2 Figure 44(a)

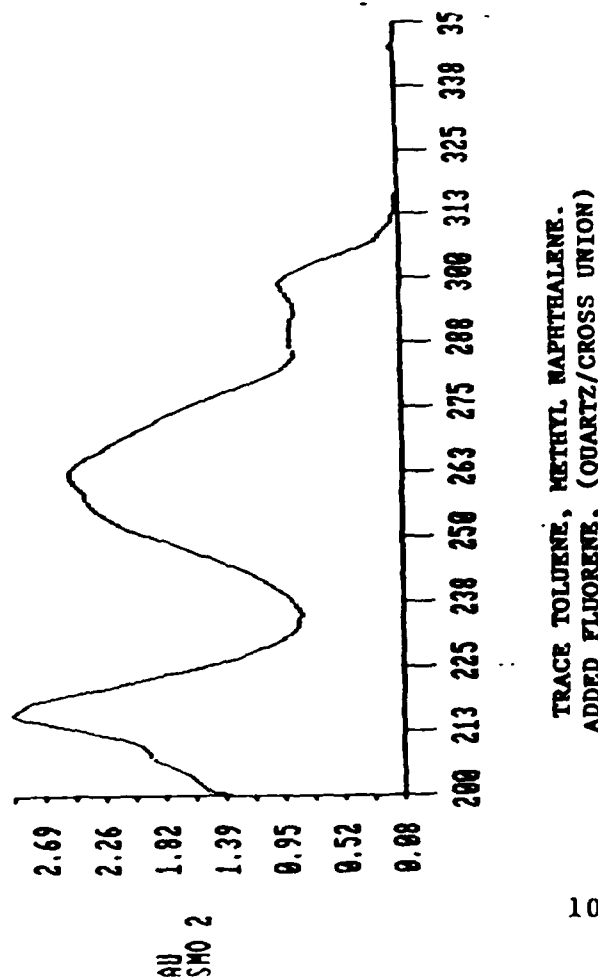


102


TRACE AROMATICS - CONFIGURATION 2 Figure 44(b)

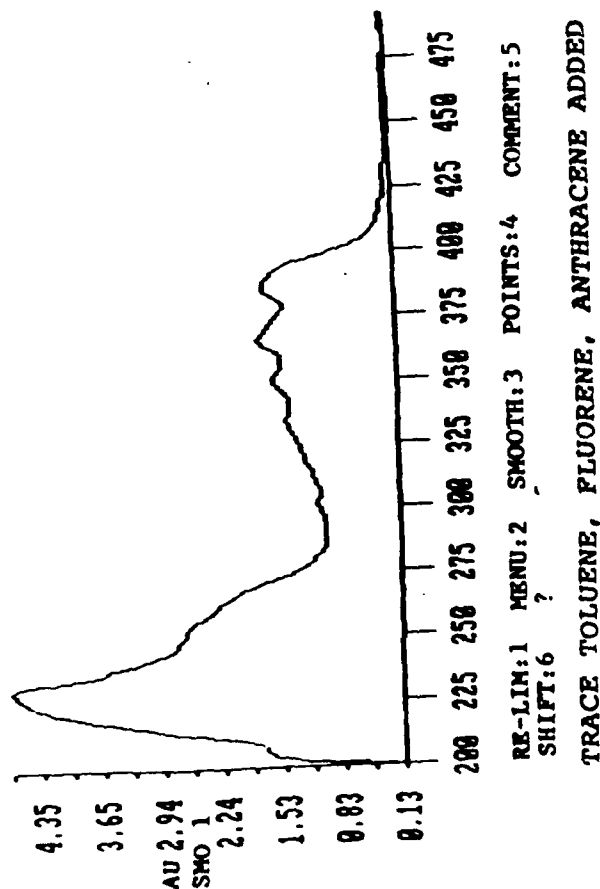

Mid-Infrared FOM Prototype in Configuration 2 (Modified) for Fabric Softener on ZnSe Rod and Fiber Optic Cable Figure 45

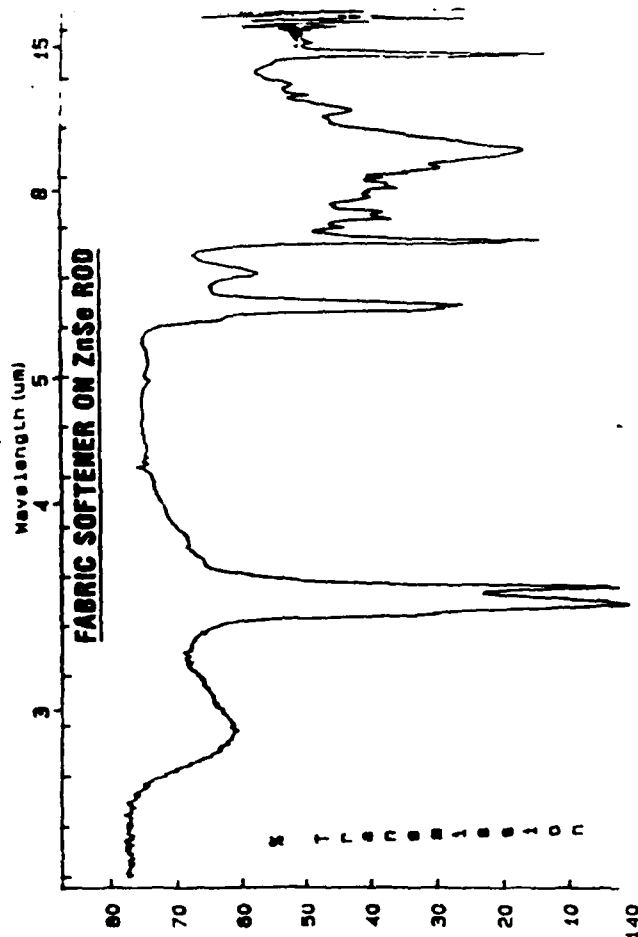

FABRIC SOFTENER ON ZnSe ROD WITH FIBER OPTIC CABLE (1.1 METER)

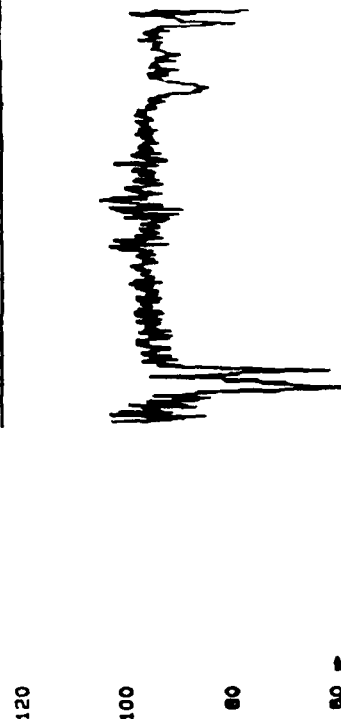

SUMMARY

Figure 46

Remote Spectroscopy and FOM Applications


- o ADVANTAGES OF REMOTE SPECTROSCOPY ARE EXTENDED TO SUPERCRITICAL FLUID EXTRACTION WITH THREE FOM PROTOTYPES IN UV-VIS-NIR AND MID-IR REGIONS
- o VIS-FOM DEMONSTRATED WITH CCS MODEL 5000 SFC/GC WITH SPE MODULE FOR ELUENTS AT HIGH CONCENTRATIONS
- o UV-FOM DEMONSTRATED WITH SPE/ACCUMULATOR MODULE FOR AROMATICS AT LOW CONCENTRATIONS WITH QUARTZ ELEMENTS
- o MID-IR FOM PROTOTYPES DEMONSTRATED WITH SPE/ACCUMULATOR MODULE FOR POLYOL SURFACTANTS WITH ZnSe ELEMENTS AND CHALCOGENIDE FIBERS/CABLES/PROBES

CONCLUSIONS

- o FOM APPLICATIONS TO SPE TECHNOLOGY PERMIT ON-LINE NON-DESTRUCTIVE ANALYSES TO SUPPLEMENT FID PROFILES WITH SELECTIVE DETECTION/IDENTIFICATION UNITS IN AUTOMATED SUPERCRITICAL FLUID INSTRUMENTATION
- o INTEGRATED SPE/ACCUMULATOR-FOM ELEMENTS WITH SOLID STATE SPECTRAL SYSTEMS WILL PROVIDE TRACE ON-LINE REAL-TIME ANALYSES FOR PROCESS R&D, QC/QA, AND ENVIRONMENTAL APPLICATIONS. POTENTIAL APPLICATION ACOUSTO-OPTIC TUNABLE FILTER (AOTF) SOLID STATE SPECTROMETER WITH FOM UNITS IN UV TO MID-IR RANGE

Applied AI - MicroEXMAT

Decision Structure for Expert System - Configuration "EXCONFIG" Figure 47(a)

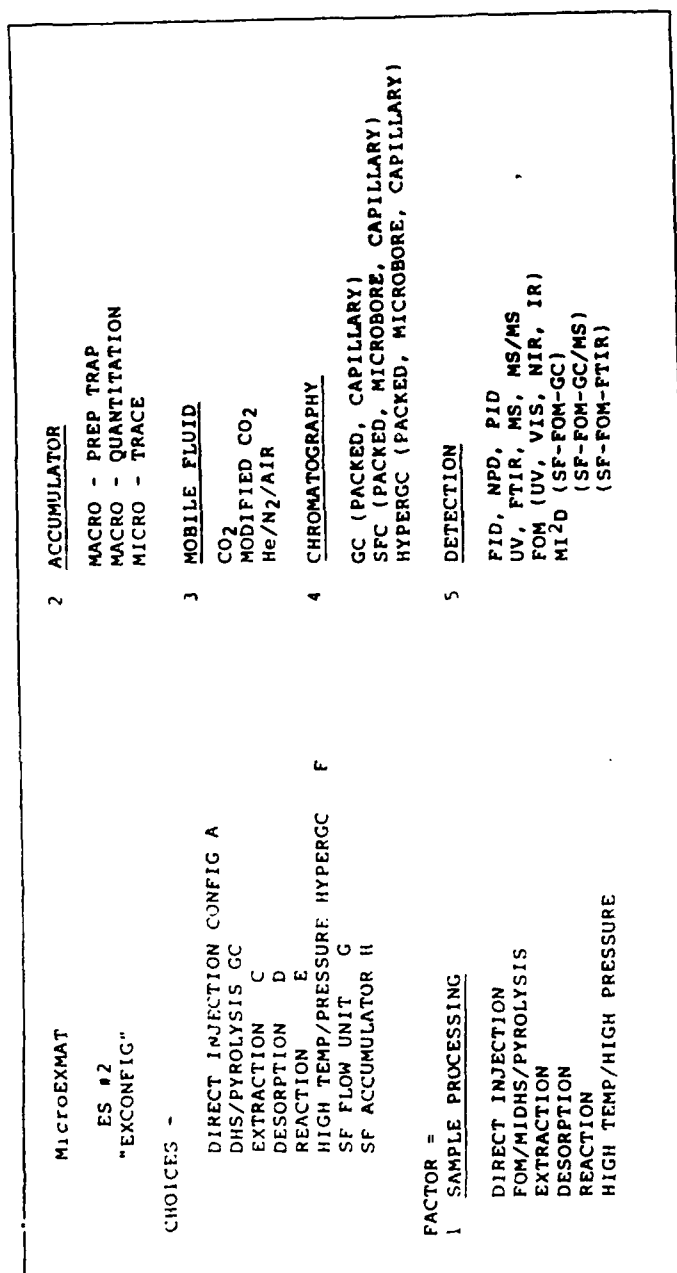

Partial Decision Structure for Expert System - Figure 47(b) Method Development "EXMETH"

Journal of Chromatographic Science, Vol. 27, March 1989

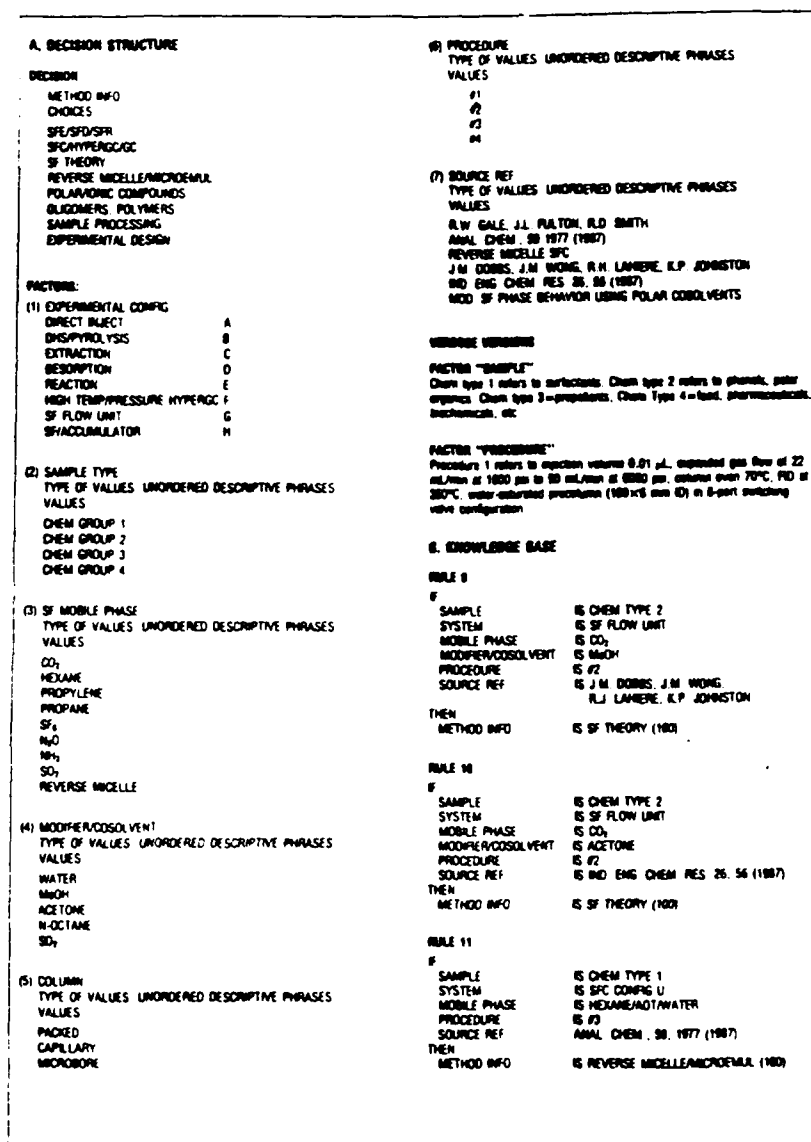


Figure 47(c)

EMBEDDED INTELLIGENCE AND APPLIED AI NETWORKS

GRC

